

3M™ Tegaderm™ High Performance Foam Non-Adhesive Dressing

Tegaderm™ High Performance Foam Non-Adhesive Dressing

Absorbent
Breathable
& Conformable



Absorbs & Evaporates


Innovative layer technology absorbs and evaporates moisture to maintain an optimal wound healing environment.

- Innovative fluid handling technology
- Unique multi-layer design provides high absorbency with breathability to reduce the risk of maceration
- Adapts to changing levels of exudate to maintain moisture balance for optimal wound healing
- Effective under compression
- Easy to use



3M™ Tegaderm™ High Performance Foam Non-Adhesive Dressing

3M™ Tegaderm™ High Performance Foam Non-Adhesive Dressing is a highly absorbent, breathable, non-adherent wound dressing. It is constructed from a conformable polyurethane foam pad, absorbent nonwoven layers, and a top layer of printed adhesive film. This waterproof film is moisture vapor permeable which prevents wound exudate strike-through and acts as a barrier to outside contamination, including bacteria and viruses.* The dressing maintains a moist wound environment, which has been shown to enhance wound healing.

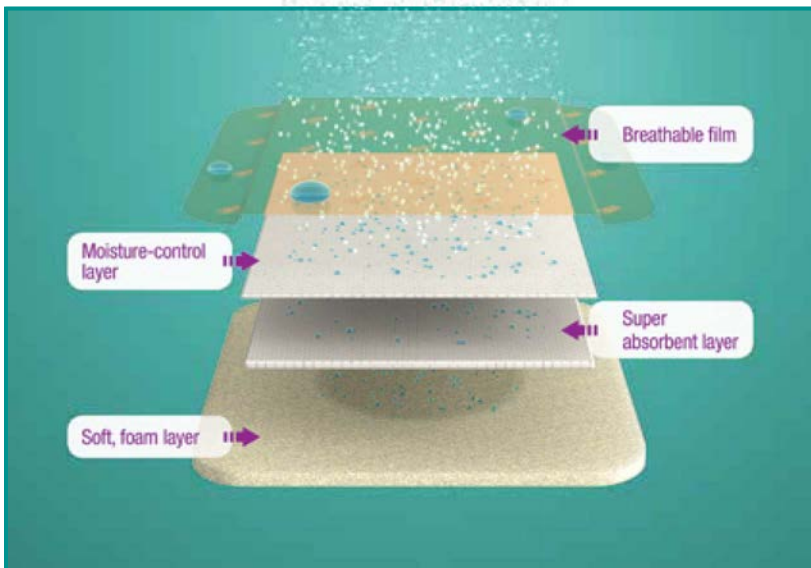

The Science Behind 3M's Excellent Fluid Handling Capacity

First, the soft foam layer rapidly wicks away exudate, which is then retained in the super absorbent layer. Next, the exudate is distributed across the moisture control layer before evaporating through the breathable film.

Indications for Use

Tegaderm™ Foam Dressing is indicated for use as a primary or secondary dressing for low- to highly-exuding partial and full thickness wounds such as:

- Pressure ulcers
- Ulcers of the lower extremity including; venous leg ulcers, arterial ulcers and neuropathic (diabetic) ulcers
- Skin tears
- Donor sites
- Surgical wounds
- Abrasions
- Superficial and partial thickness burns
- Suitable for use around tube exit sites
- Suitable for use with compression therapy


*In vitro testing shows that the transparent film provides a viral barrier from viruses 27nm in diameter or larger while the dressing remains intact without leakage.

Follow Universal Precautions and facility guidelines for infection control.


Application Instructions

- 1 Cleanse wound and surrounding skin according to facility policy.
- 2 Select appropriate size dressing to ensure foam pad is larger than the wound area. The dressing may be cut to accommodate wound contours.

- 3 Remove the Tegaderm™ High Performance Foam Non-Adhesive Dressing from package.
- 4 Position the dressing, printed side up, over the wound with edges overlapping onto intact skin. Secure with elastic or cohesive wrap, tape or other methods (see Dressing Secural Options).

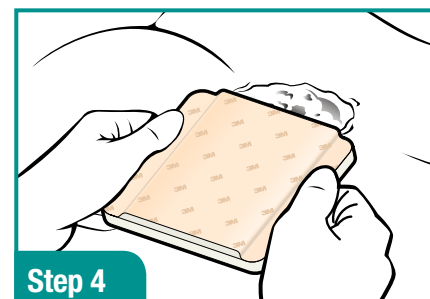


Dressing Secural Options

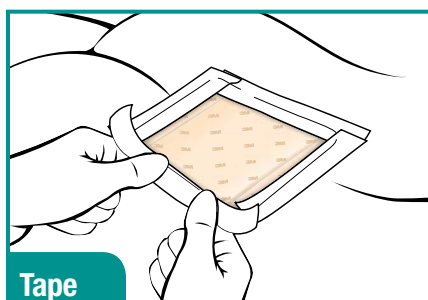

Tape


Tape: Picture frame the dressing with tape (3M™ Kind Removal Silicone Tape, 3M™ Medipore™ Soft Cloth Tape or 3M™ Micropore™ Tape).

Note: To maintain moisture vapor evaporation, tape only edges of the dressing in picture frame manner.

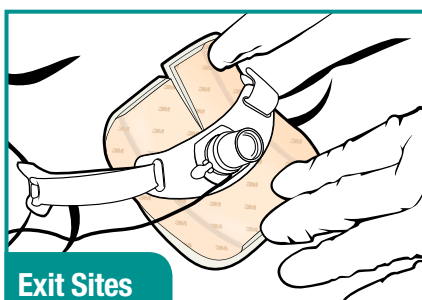

Exit Sites

Tube Exit Sites: When using the dressing for absorption and cushioning of tube sites (such as tracheostomy or gastric tubes), the dressing may be held in place by the tube itself or secured by taping the edges of the dressing in a picture frame manner.

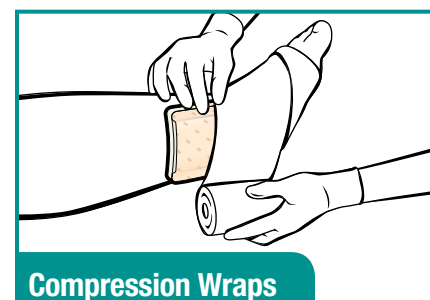

Compression Wraps

Lower Extremity Compression Wraps:


The dressing may be used under compression wrap systems such as 3M™ Coban™ 2 Layer Compression System. Additional taping is typically not required.



Wraps

Wraps: When dressing an extremity, secure the dressing with a gauze wrap, a self-adherent wrap (3M™ Coban™ Self-Adherent Wrap), tubular stockinette or an elastic wrap.

Removal Instructions

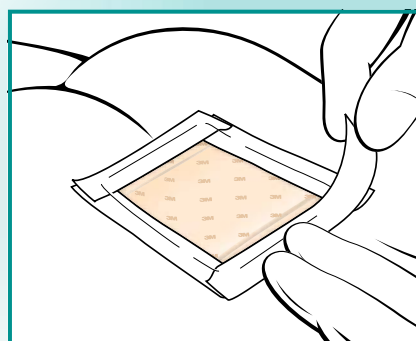

- 1 Gently remove adhesive tape border or overwrap.
- 2 Lift foam dressing from the wound bed. If dressing is adhered to the wound surface, saturate with normal saline and gently loosen.
- 3 Discard dressing according to your institution's policy.







3M™ Tegaderm™ High Performance Foam Non-Adhesive Dressing

Considerations

Observe the dressing frequently. As the dressing absorbs, the exudate will wick to the top of the dressing, and discoloration may be noticeable. When the exudate spreads to the edges of the dressing or the dressing begins to leak, a dressing change is indicated.

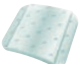- Tegaderm™ Foam Dressing may remain in place for up to seven days. Dressing change frequency will depend on the type of wound, volume of exudate, and clinical situation. Change at least every seven days, or as indicated per treatment protocol.
- The dressing may be used as an absorbent cover dressing over wound fillers such as alginates or gauze packing.
- The dressing may be cut as necessary to accommodate wound size and contours. Use of non-sterile product may increase the risk of contamination and infection.
- Covering the entire foam dressing with an occlusive tape may reduce the breathability of the film backing. To maintain moisture vapor exchange, tape only edges of the dressing in picture frame manner.

Ordering Information - [Click here to see options and pricing](#)

	Catalog No.	Size	Dressings/Box	Boxes/Case	HCPCS Code
	90600	2 in x 2 in (5cm x 5cm)	10	4	A6209
	90601	4 in x 4 in (10cm x 10cm)	10	4	A6209
	90602	4 in x 8 in (10cm x 20cm)	5	6	A6210
	90603	8 in x 8 in (20cm x 20cm)	5	6	A6211
	90604	3.5 in x 3.5 in (8,8cm x 8,8cm) <i>Fenestrated</i>	10	4	A6209
	90605	4 in x 24 in (10cm x 60cm), Roll	1 roll/pkg.	6 rolls	A6211

Disclaimer:

HCPCS and CPT codes have been provided to assist you in the preparation of insurance claims. Please note, however, that the reimbursement information provided by 3M Health Care and its representatives is intended to provide general information relevant to coverage and coding for 3M products. Insurers' reimbursement policies can vary and the use of the codes discussed here does not guarantee that an insurer will cover or pay at any particular level. Health care providers should exercise independent clinical judgment in choosing the codes which most accurately describe the products provided.

For more information visit our website at www.3M.com/tegaderm, contact your 3M Critical & Chronic Care Solutions representative, or call the 3M Health Care Customer Helpline at 1-800-228-3957.

These products can be ordered from your local distributor. Outside the United States contact the local 3M subsidiary.



Critical & Chronic Care Solutions Division
3M Health Care
2510 Conway Avenue
St. Paul, MN 55144
USA
1-800-228-3957
www.3M.com/HealthCare

3M™ Tegaderm™ SIMPLE. DEPENDABLE. TRUSTED.
Dressings

3M, Coban, Medipore, Micropore and Tegaderm are trademarks of 3M. Used under license in Canada. Please recycle. Printed in U.S.A.
© 3M 2014. All rights reserved. 70-2011-5477-3