

Philips Healthcare is part of
Royal Philips Electronics

How to reach us

www.philips.com/healthcare
healthcare@philips.com

Asia

+49 7031 463 2254

Europe, Middle East, Africa
+49 7031 463 2254

Latin America

+55 11 2125 0744

North America

+1 425 487 7000

800 285 5585 (toll free, US only)

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service

+1 724 387 4000

800 345 6443 (toll free, US only)

Philips Respironics International
Headquarters

+33 1 47 28 30 82

CoughAssist, Cough-Trak and Respironics are trademarks of Koninklijke Philips Electronics N.V. and its affiliates. All rights reserved.
DisCide is a registered trademark of AllPied, Inc.

Please visit www.philips.us/CoughAssistT70



©2012 Koninklijke Philips Electronics N.V.

All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.

Gezys JH 03/26/12 MCI 4100796 PN 1 099607



My CoughAssist

A patient guide to CoughAssist T70

PHILIPS
RESPIRONICS
sense and simplicity


A natural part of life

Coughing. We often think of a cough as a symptom of a cold or a means to clear our throat. Truth is, the ability to cough is essential to life. Normally, the lining of the airways produces small amounts of mucus which trap dirt and bacteria. If allowed to accumulate, these secretions can result in infection. Coughing is the body's way of clearing mucus from the lungs.

Muscles located in our rib cage, neck and abdomen play an important role during coughing. With many diseases, a weakened ability to remove secretions increases the chance of a respiratory infection.

If coughing becomes difficult, your physician may prescribe the use of a manual and/or mechanical-assisted coughing technique. A manual cough assist involves applying firm and rapid pressure to the upper abdomen, forcing air out of the lungs.

CoughAssist T70 is a mechanical insufflation-exsufflation device that does just that. It simulates a natural cough by gradually delivering a large volume of air to your lungs, similar to a normal deep breath, when you breathe in (positive pressure), and then quickly reversing the flow of secretions out (negative pressure). CoughAssist T70 makes your cough more effective to help keep your airways clear, reducing the chance of respiratory infections. It also incorporates a triggering feature called Cough-Trak that allows you to synchronize treatment with your own breathing patterns.

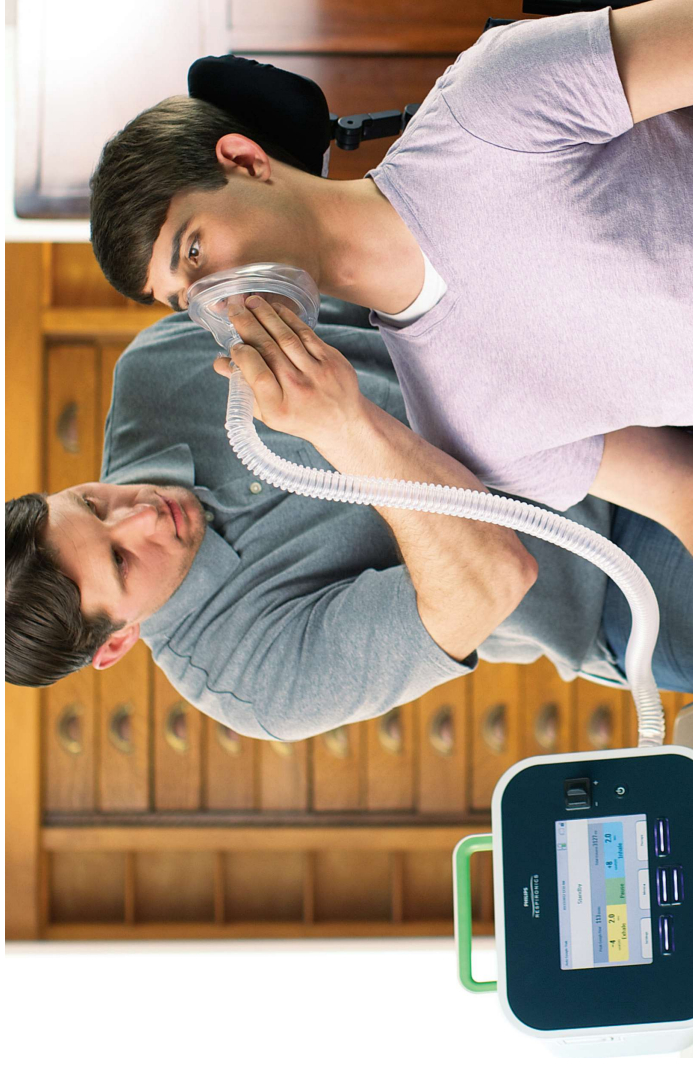

Product features

- 1 — Left button. Allows you to select display options or perform certain actions specified on the screen.
- 2 — Up/down button. Allows you to navigate the display menu.
- 3 — Right button. Allows you to select display options or perform certain actions specified on the screen.
- 4 — Start/stop button. Turns the device on or off.
- 5 — Manual switch. Activates the insufflation and exsufflation phases. Pressing switch to the right (+) activates the insufflation phase. Pressing switch to the left (-) activates the exsufflation phase. The center position activates the pause phase.
- 6 — Breathing circuit connection. Connect your circuit to this connector on the device.
- 7 — Display screen. Allows you to view settings and system status information.




- 1 — Fan exhaust. Air from inside the device is expelled here.
- 2 — AC power inlet. Connect the AC power cord here.
- 3 — Remote control connector. Connect the optional remote control accessory (foot pedal) here.
- 4 — Connect a USB cable to this connector for service only.
- 5 — SpO₂ connector. If you are using the optional oximetry accessory connect it here.
- 6 — DC power inlet. Connect an external battery here using the Philips Respironics DC power cord.
- 7 — Airpath outlet where air exits the device.
- 8 — Airpath inlet where outside air enters the device. Insert the filter supplied with the device here.
- 9 — Detachable battery slot. Insert optional Lithium Ion detachable battery here.
- 10 — Tubing retainer. Route the tubing and mask through this bracket for proper tubing management when the device is not in use.



How to set up and use my device

- Plug the device into an AC power source, DC power source, or use the detachable battery
- Connect the filter and tubing to the device and the mask, mouthpiece or trach adaptor to the tube
- Press the start/stop button to turn the device on and enter Standby
- Position the mask, mouthpiece or trach adaptor on the patient
- Press the right (Therapy) button to begin therapy – in the Auto mode, the device will cycle from inhale (positive) to exhale (negative); if in Manual mode, slide the control lever from the right (inhale) to the left (exhale)
- Press right (Standby) button to end therapy

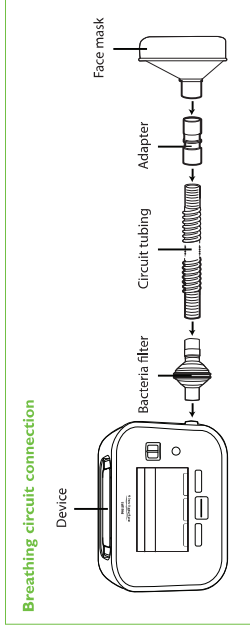

Getting started


When starting treatment with CoughAssist T70, you may want to begin with lower pressures and low inhale flow to become comfortable with the therapy. Be aware that lower pressures may not be effective. You may also want to agree on a signal with your caregivers in the event you wish to alert them that you want to stop treatment.

A cough cycle is composed of an inhale time, an exhale time, and a pause time. Four to six cycles represent a sequence. Between each sequence you may rest for 20 to 30 seconds or more. This resting time allows you or your caregiver to remove secretions from the mask, mouthpiece, or trach tube. Four to six sequences represent a typical treatment. The number of treatments you take each day may vary.

During a treatment you should sit upright as tolerated.

Treatment delivery can be manual or automatic. This means that cycles can be started either by the manual switch, foot pedal, or automatically, based on your inspiratory effort or a set time. Cough-Trak is available in Auto mode, allowing you to control when you want to start the inspiratory phase. CoughAssist T70's Cough-Trak will detect when you are ready to inhale and will deliver the set pressure.

Your device allows up to three prescription settings. Your clinician will advise you on which prescription preset to use.

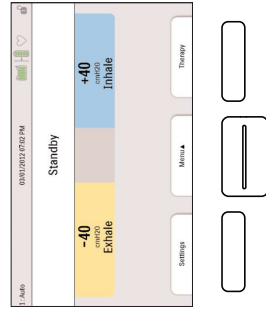

Patient interface

- When using a mouthpiece, nose clips may be required.
- A well-fitting, soft-seal mask is important to eliminate possible leakage.
- Some patients may need to have their head supported from behind to prevent leaks. In order to be effective, firm pressure may be required.
- If the patient has a cuffed tracheostomy tube, inflate the tube cuff during the CoughAssist treatment.

Always follow the treatment plan that your physician has prescribed.

Information log

If the Information log picture **1** is displayed on the screen, it means information is available in the information log screen but it does not necessarily affect therapy. To access these messages:

- From the Standby screen, press up (Menu) key, press down (Navigate) key to "Information Log."
- Press right (Select) key. One of the following messages will appear:

Message	User action
Battery not charging – Temp	The detachable battery is too hot and cannot be charged. Allow device and/or battery to cool to resume charging. Remove battery and charge with battery charger accessory. If battery still cannot be charged, contact your healthcare provider.
Check external battery	Power is being drawn from the detachable battery even though the external battery is usable. Replace external battery cable or external battery. If problem persists, contact your healthcare provider.
Detachable battery not charging	The detachable battery cannot be charged. Replace battery. If problem persists with a different battery, contact your healthcare provider.
Replace detachable battery	Detachable battery has failed or has reached end of life. Replace battery. If problem persists with a different battery, contact your healthcare provider.
Internal fan failure – see User Manual	The internal fan is not working. Device should not be used with oxygen added to the patient circuit. Before using oxygen with this device, contact your healthcare provider.
Card error	The device cannot write to or read from the SD card. Remove the SD card and use another card, if available. If problem persists, contact your healthcare provider.

Cleaning and maintenance

The device

Unplug the device and clean the external surface using one of the following cleaning agents:

- A clean cloth dampened with water and a mild detergent
- 70% Isopropyl alcohol
- DisCide towelettes
- 10% chlorine bleach solution

Allow the device to dry completely before plugging in the power cord. The device does not require routine servicing.

The breathing circuit, interface and bacteria filter

- After each use, the breathing hose and patient interface should be washed thoroughly with mild liquid detergent and water. The parts need to be completely dried before each use.
- The bacteria filter cannot be washed. It should be replaced when blocked with mucus or trapped moisture.
- Replace the mask if it becomes torn or if you can no longer maintain a good seal.

The air filter

Under normal usage, you should clean the air filter at least once every two weeks and replace it with a new one every six months.

- If the device is operating, stop the airflow. Disconnect the device from the power source.

- Remove the filter from the enclosure.
- Examine the filter for cleanliness and integrity.
- Wash the filter in warm water with a mild detergent. Rinse thoroughly to remove all residue.
- Allow the filter to air dry completely. Replace the filter if it is torn or damaged. Only Philips Respironics filters should be used as replacement filters.
- Reinstall the filter.

SD card

The CoughAssist T70 comes with an SD card inserted in the SD card slot on the side of the device. Your clinician may use this card to record information and may ask you to remove the card and forward for evaluation. To remove the SD card:

- Select the "Safely Remove SD Card" option from the Main menu.
- After the "Safely remove SD Card" message appears, remove the card.



Optional accessories

Description	Item no.
Detachable battery	1043570
Detachable battery charger	1072182
DC auto adapter	1067424
Foot pedal	1059017
Roll stand	1098655

My doctor:

My homecare provider:

My mask and tubing size:

My prescription:

Preset 1
Mode: Manual Automatic
 On Off
Cough-Trak (inspiratory trigger):
Inhale pressure: _____ cm H₂O High NA
Inhale flow: Low Medium High
Inspiratory time: _____ sec. NA
Exhale pressure: _____ cm H₂O NA
Expiratory time: _____ sec. NA
Pause time: _____ sec. NA

Preset 3
Mode: Manual Automatic
 On Off
Cough-Trak (inspiratory trigger):
Inhale pressure: _____ cm H₂O High NA
Inhale flow: Low Medium High
Inspiratory time: _____ sec. NA
Exhale pressure: _____ cm H₂O NA
Expiratory time: _____ sec. NA
Pause time: _____ sec. NA

Preset 2
Mode: Manual Automatic
 On Off
Cough-Trak (inspiratory trigger):
Inhale pressure: _____ cm H₂O High NA
Inhale flow: Low Medium High
Inspiratory time: _____ sec. NA
Exhale pressure: _____ cm H₂O NA
Expiratory time: _____ sec. NA
Pause time: _____ sec. NA

Recommendations:
