INSTRUCTIONS FOR USE

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JOBST® Bella™ Lite & Strong

- | | | |
|--|--|--|
| (EN) Medical compression garment | (NO) Medisinsk kompresjonsplagg | (HU) Gyógyászati kompressziós ruhadarab |
| (DE) Medizinische Kompressionsversorgung | (DA) Medicinsk kompressionsprodukt | (RO) Articole de îmbrăcămințe compresivă medicinală |
| (FR) Vêtement de compression médicale | (FI) Lääkinnällinen kompressiovaate | (CS) Zdravotní kompresní výrobky |
| (ES) Prenda de compresión médica | (IS) Þrýstifatnaður í lækningaskyni | (SL) Medicinsko kompresijsko oblacilo |
| (PT) Peça de vestuário de compressão médica | (ET) Meditsiinilised kompressioonrõivad | (SK) Lekárske kompresné odevy |
| (IT) Indumenti compressivi medicali | (LT) Medicininė kompresinė apranga | (RU) Медицинская компрессионная одежда |
| (NL) Medische compressiekleding | (LV) Medicīniskais kompresijas apģērbs | (TR) Tibbi kompresyon giysisi |
| (SV) Medicinskt kompressionsplagg | (PL) Medyczna odzież uciskowa | (AR) رداء ضاغط طبي |
| | | (HE) כפפות ושרוול ל לחץ רפואיים |

FIBER CONTENT

- EN

FIBER CONTENT
- DE

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COMPOSITION
- ES

CONTENIDO DE FIBRA
- PT

TEOR DE FIBRA
- IT

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VEZELGEHALTE
- SV

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NO

FIBERINNHOLD

DA

FIBERINDHOLD

FI

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IS

TREFAJAINNIHALD

ET

KIULINE KOOSTIS

LT

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LV

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PL

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HU

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RO

CONȚIN FIBRE

CS

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SL

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SK

OBSAH VLÁKIEŇ

RU

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Brand ¹ /CCL Marque/CCL	Elastane (Spandex) ² Élasthanne (Spandex)	Polyamide (Nylon) ³ Polyamide (Nylon)	Silicone ⁴ Silicone
JOBST® Bella Strong Ready-to-Wear Armsleeve with Silicone Band, 15-20 mmHg	18 %	74 %	8 %
JOBST® Bella Strong Ready-to-Wear Armsleeve with Silicone Band, 20-30 mmHg	20 %	72 %	8 %
JOBST® Bella Strong Ready-to-Wear Armsleeve with Silicone Band, 30-40 mmHg	23 %	69 %	8 %
JOBST® Bella Strong Ready-to-Wear Armsleeve, 15-20 mmHg	20 %	80 %	–
JOBST® Bella Strong Ready-to-Wear Armsleeve, 20-30 mmHg	24 %	76 %	–
JOBST® Bella Strong Ready-to-Wear Armsleeve, 30-40 mmHg	26 %	74 %	–
JOBST® Bella Strong Ready-to-Wear Gauntlet, 15-20 mmHg	9 %	91 %	–
JOBST® Bella Strong Ready-to-Wear Gauntlet, 20-30 mmHg	9 %	91 %	–
JOBST® Bella Strong Ready-to-Wear Glove, 15-20 mmHg	30 %	70 %	–

COMPOSITION

Brand ¹ /CCL Marque/CCL	Elastane (Spandex) ² Élasthanne (Spandex)	Polyamide (Nylon) ³ Polyamide (Nylon)	Silicone ⁴ Silicone
JOBST® Bella Strong Ready-to-Wear Glove, 20-30 mmHg	35 %	65 %	–
JOBST® Bella Lite Ready-to-Wear Armsleeve, 15-20 mmHg	16 %	84 %	–
JOBST® Bella Lite Ready-to-Wear Armsleeve with Silicone Band, 15-20 mmHg	17 %	80 %	3 %
JOBST® Bella Lite Ready-to-Wear Armsleeve, 20-30 mmHg	20 %	80 %	–
JOBST® Bella Lite Ready-to-Wear Armsleeve with Silicone Band, 20-30 mmHg	22 %	75 %	3 %
JOBST® Bella Lite Ready-to-Wear Gauntlet, 15-20 mmHg	16 %	84 %	–
JOBST® Bella Lite Glove, 15-20 mmHg	41 %	59 %	–
JOBST® Bella Lite Ready-to-Wear Gauntlet, 20-30 mmHg	14 %	86 %	–
JOBST® Bella Lite Ready-to-Wear Combined Armsleeve, 15-20 mmHg	19 %	81 %	–
JOBST® Bella Lite Ready-to-Wear Combined Armsleeve with Silicone Band, 15-20 mmHg	15 %	77 %	8 %
JOBST® Bella Lite Ready-to-Wear Combined Armsleeve, 20-30 mmHg	21 %	79 %	–
JOBST® Bella Lite Ready-to-Wear Combined Armsleeve with Silicone Band, 20-30 mmHg	17 %	75 %	8 %
JOBST® Bella Lite Glove, 20-30 mmHg	41 %	59 %	–

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INSTRUCTIONS FOR USE

EN

JOBST® Bella Lite & Strong Ready-to-Wear Armsleeve

JOBST® Bella Lite: 15-20 mmHg and 20-30 mmHg

JOBST® Bella Strong: 15-20 mmHg, 20-30 mmHg and 30-40 mmHg

UK: Class 1 (CI 1) = 15-20 mmHg, Class 2 (CI 2) = 20-30 mmHg, Class 3 (CI 3) = 30-40 mmHg

Your healthcare professional has recommended this JOBST® Bella Armsleeve for the management of your lymphedema or other edema.

Indications

Every JOBST® medical compression product has a specific compression level designed to help specific venous and lymphatic conditions.


The compression level is printed on the package label. The required compression class must be determined by your physician.

For the management of lymphedema and other edema. Phlebitis, post thrombotic syndrome, vascular malformation. JOBST® Bella Lite and JOBST® Bella Strong 20-30 mmHg for post liposuction management.

Post surgery management and management of oedema (20-30 mmHg JOBST® Bella Lite only).

Instructions for Donning JOBST® Bella Lite Armsleeve


1. Fold armsleeve on itself, from the top down. Note positioning of the wide band. Insert hand.
2. Gently pull complete armsleeve on to the elbow. The armsleeve should now be positioned over the wrist and elbow portion of the arm.
3. Gently pull the top of the armsleeve up to its correct end location (top of the arm).
4. Smooth out any wrinkles or folds that may appear.

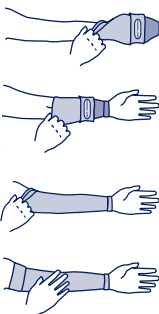
Optional: A donning aid may be used to help slide the sleeve over your hand.

Instructions for Donning JOBST® Bella Strong Armsleeve

1. Pull Donning Aid onto arm.

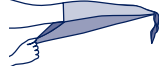2. Fold back band on armsleeve so silicone is exposed.
3. Pull armsleeve onto arm until wrist cuff is on or above wrist bone.
4. Place loop over door knob or other secure, stationary object. Holding silicone band in place, slowly step away from door - allowing Donning Aid to be removed. Adjust armsleeve as needed.
5. Flip silicone band into place using a donning glove to adjust the armsleeve as needed.

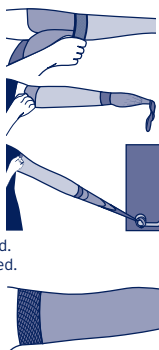JOBST® Bella Lite & Strong Ready-to-Wear Gauntlet

15-20 mmHg and 20-30 mmHg

UK: Class 1 (CI 1) = 15-20 mmHg, Class 2 (CI 2) = 20-30 mmHg

Your healthcare professional has recommended this JOBST® Bella Gauntlet for the management of your lymphedema or other edema.

Indications

Every JOBST® medical compression product has a specific compression level designed to help specific venous and lymphatic conditions.

The compression level is printed on the package label. The required compression class must be determined by your physician.

For the management of lymphedema and other edema.

Instructions for Donning

1. Slip fingers into gauntlet starting at the end that has the widest band.
2. Position gauntlet so that it sits comfortably around the palm and thumb.
3. The gauntlet should be positioned so that it ends around the knuckles.

JOBST® Bella Lite & Strong Ready-to-Wear Glove

15-20 mmHg and 20-30 mmHg

Your healthcare professional has recommended this JOBST® Bella Glove for the management of your lymphedema or other edema.

Indications

Every JOBST® medical compression product has a specific compression level designed to help specific venous and lymphatic conditions.

INSTRUCTIONS FOR USE

EN

The compression level is printed on the package label. The required compression class must be determined by your physician.

Bella Lite Glove 15-20, 20-30: For the control of mild to moderate hand swelling from acute and chronic edema, lymphatic and other swelling disorders.

Bella Strong Glove 15-20: For the control of mild hand swelling from acute and chronic edema, lymphatic and other swelling disorders.

Bella Strong Glove 20-30: For the control of moderate to severe hand swelling from acute and chronic edema, lymphatic, and other swelling disorders.

Follow your healthcare provider's instructions.

Instructions for Donning JOBST® Bella Lite Glove

Open wounds should be covered before applying glove. Remove jewelry before applying. The glove can be used on either hand:

- A. Right hand: label is visible on the wrist.
- B. Left hand: turn the glove "inside out", inverting the fingers. The label is on the inside.

Follow these steps when putting on your glove:

1. Slip fingers into glove.
2. Position glove so that it sits comfortably around the palm and thumb.
3. The glove should be positioned so that the fingers end near the beds of the fingernails.
4. If the finger lengths are too long, mark the finger at the appropriate length, remove the glove, then trim the length of the glove finger directly above the marked area and remove the loose thread to lock down the seam. Be careful not to remove too much of the finger length.

Instructions for Donning JOBST® Bella Strong Glove

Follow these steps when putting on your glove:

1. Slip your hand through the wrist of the glove and into the fingers.
2. Adjust the fingers by carefully pulling on the fabric until the glove fits snugly into the web spaces between your fingers.
3. Be sure the fabric is smoothed out over the fingers and hand.

JOBST® Bella Lite Ready-to-Wear Combined Armsleeve

15-20 mmHg and 20-30 mmHg

Your healthcare professional has recommended this JOBST® Bella compression garment for the management of your lymphedema or other edema.

Indications


Every JOBST® medical compression product has a specific compression level designed to help specific venous and lymphatic conditions.

The compression level is printed on the package label. The required compression class must be determined by your physician.

For the management of lymphedema and other edema. Phlebitis, post thrombotic syndrome, vascular malformations, post liposuction (20-30 mmHg only).

Instructions for Donning

1. Fold armsleeve part on itself, from the top down. Slip fingers into gauntlet and position it so that it sits comfortably around the palm and thumb.
2. Gently pull complete armsleeve part on to the elbow. The armsleeve part should now be positioned over the wrist and elbow portion of the arm.
3. Gently pull the top of the Armsleeve part up to its correct end location (top of the arm).
4. Smooth out any wrinkles or folds that may appear.

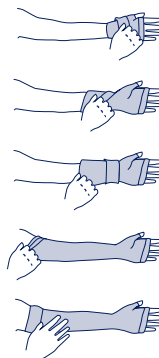Contraindications

JOBST® Bella Lite & Strong

JOBST® medical compression garments should not be worn in cases of:

- Advanced arterial insufficiency including ischemia
- Uncontrolled congestive heart failure
- Untreated septic phlebitis
- Phlegmasia coerulea dolens
- Incompatibility to fabric

Cautions JOBST® Bella Lite & Strong

- Weeping dermatosis
- Cutaneous infections
- Open wounds must be covered with an appropriate dressing before the compression garment is put on

INSTRUCTIONS FOR USE

EN

- Severely compromised skin sensibility and impaired sensitivity of the limb
- Advanced peripheral neuropathy
- Immobility (confinement to bed). Not recommended to wear during bed rest
- Conditions in which increased venous and lymphatic return is not desired.
- Appropriate hygiene and skin care should be applied
- Ill-fitting compression garments may cause constriction leading to e.g. skin necrosis, pressure ulcers and pressure related peripheral nerve damage
- The attending physician decides whether compression in general can be applied in patients suffering from the following diseases: rheumatoid arthritis, complex regional pain syndrome, CRPS (M. Sudeck), malignant lymphedema, gangrene.
- Medical compression garments should only be used after consultation with your doctor or therapist and checked regularly by them on the appropriateness of your therapy, also considering age and health status.
- Patient measurements should be taken by a trained person to determine the correct product size needed and to decide whether a made-to-measure or ready-to-wear garment is required.
- When re-ordering a garment after 6 month we recommend to take the patients measurements again as the limb circumference might have changed.
- If you notice discomfort, pain, numbness, changes in skin color, appearance of new wounds or swelling during wear, stop using the product and contact your healthcare professional.

- Any serious incident that has occurred in relation to this device should be reported to BSN medical Inc. (BSN-JOBST® GmbH within the European Union) and the competent authority of your state.

Care Instructions

JOBST® Bella Lite & Strong

Machine washable in warm water. Mild process. Use a delicate fabric detergent. Do not use fabric softener. Do not bleach. Tumble drying possible. Low temperature; exhaust temperature max. 50 °C. Do not iron. Do not dry clean. Flat drying. Drip line drying in the shade.

Special Care Instructions for JOBST®

Bella Armsleeves with Silicone band

1. Before washing the garment, turn the garment inside out.
2. Thoroughly hand wash the silicone band using a delicate fabric detergent. Concentrate on removing any lotion, powder or hair from the silicone.

The use of certain skin care products, dry peeling skin or use in hairy areas, may prevent the silicone band from keeping the garment in place.

JOBST® garments will retain their compression for 6 months under normal wear and care.

Please protect your compression garments from direct exposure to sunlight, heat, and moisture. We recommend storage at room temperature.

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