Lansinoh®

SIGNATURE PRO™
DOUBLE ELECTRIC BREAST PUMP
INSTRUCTIONS FOR USE

Important Safeguards

The Lansinoh® Signature Pro™ Double Electric Breast Pump is a personal care item and for hygienic reasons should never be shared between mothers or resold from one mother to another. Either of these actions will void the warranty.

If you are a mother who is infected with Hepatitis B, Hepatitis C or Human Immunodeficiency Virus (HIV), pumping breastmilk using the Lansinoh® Signature Pro™ Double Electric Breast Pump will not reduce or remove the risk of transmitting the virus to your baby through your breastmilk.

When using electrical products, especially when children are present, basic safety precautions should always be maintained, including the following:

Read All Instructions Before Using

Danger: To reduce the risk of electrocution:

- Always unplug this product immediately after using.
- Do not use while bathing.
- Do not place or store this product where it can fall or be pulled into a tub or sink.
- Do not place or drop pump into water or other liquid.
- Do not reach for an electrical product that has fallen into water. Unplug immediately.

Warning: To reduce the risk of burns, electrocution, fire, or injury to persons:

- This product should never be left unattended when plugged in.
- Close supervision is necessary when this product is used near children.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening or tubing.
- Do not use this product if it is broken or any parts are missing.
- Do not allow small children or pets to play with the pump unit, AC adapter, etc.
- Do not use outdoors with the AC adapter or operate where aerosol spray products are being used or where oxygen is being administered.

**Save These Instructions
along with your receipt for future reference.**

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Congratulations!

Your decision to provide your baby with breastmilk, even if you need to be apart, is one of the most precious gifts you can give your baby. The Lansinoh® Signature Pro™ is a safe, efficient and effective choice for moms who need to pump regularly.

Contents	Page
1. QUICK START GUIDE	4-5
2. PARTS DIAGRAMS	6-8
2.1 Pump Diagram	
2.2 Pump Parts	
3. CLEANING INSTRUCTIONS	9
3.1 Before First Use	
3.2 Cleaning	
3.3 Before Each Use	
4. ASSEMBLY	10-12
4.1 Assembling the Breastmilk Expression Set	
4.2 Connecting the Breastmilk Expression Set	
4.3 Connecting for Single Pumping	
5. USING THE LANSINOH SIGNATURE PRO	13-16
5.1 Power Options	
5.2 Getting Started	
5.3 2-Phase Technology	
5.5 Customization for Comfort and Efficiency	
5.6 Completing your Pumping Session	
6. PUMPING TIPS	17
7. TROUBLESHOOTING	18
8. BREASTMILK STORAGE GUIDELINES	19
9. INTRODUCING THE NATURALWAVE® NIPPLE	20
10. TECHNICAL DATA	21
11. WARRANTY	21-22

PLEASE DO NOT RETURN THIS PUMP TO THE RETAILER.

Lansinoh's dedicated Customer Service team is here to help. If you experience any issues with your Lansinoh Signature Pro™ breast pump please call 800-292-4794 or email customerservice@lansinoh.com, Monday - Friday 8am-8pm EST.

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QUICK START GUIDE

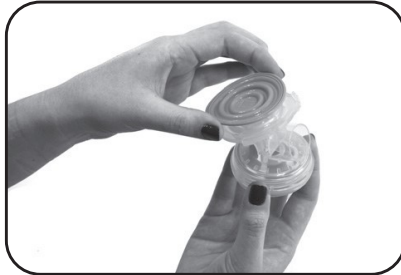
Do NOT wash or sterilize tubing

Before first use, sterilize all parts that will come into contact with milk according to instruction on pg 9, section 3.1.

Please carefully read all instructions prior to use

Assembly

See detailed instruction on pgs 10-11. Proper assembly is essential to ensure proper function and suction.

Insert purple diaphragm with arrows pointing forward

Attach clear cover and twist so tubing connector faces back

Attach flange

Insert valve

Attach bottle

Attach tubing

Secure Y connector

Secure connector

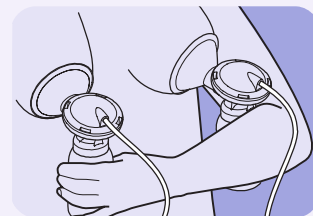
Attach to pump base

Connect pump to power using AC Adapter or six AA Alkaline Batteries (not included)

Getting Started

Hold Breastmilk Expression Set to your breasts so that the nipple is aligned to extend straight toward the nipple tunnel of the breast flange, ensuring there are no gaps for air to enter. Make sure set is tilted slightly downward to allow milk to flow into bottles.

2-Phase Technology – Using Let-Down and Expression Phases

Press the power button  to start the pump in let-down phase . Let-down phase will run automatically for two minutes.

To begin expression phase sooner, press the Let-down/Expression Toggle button  to move to the expression phase. Press the Let-down/Expression Toggle button again at any point to switch back to Let-down phase.

Customization for Comfort and Efficiency

Pumping Style Selection

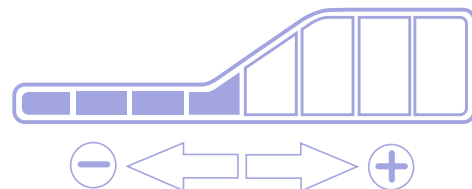
Every baby feeds differently. The Signature Pro™ offers three different pumping styles to choose from, allowing you to choose the option that is most like your baby's feeding patterns. Change the pumping style using the Style Mode Toggle button  to find the pumping style that works best for you.

NOTE: When you switch to the expression mode the pump will start in Style 2.

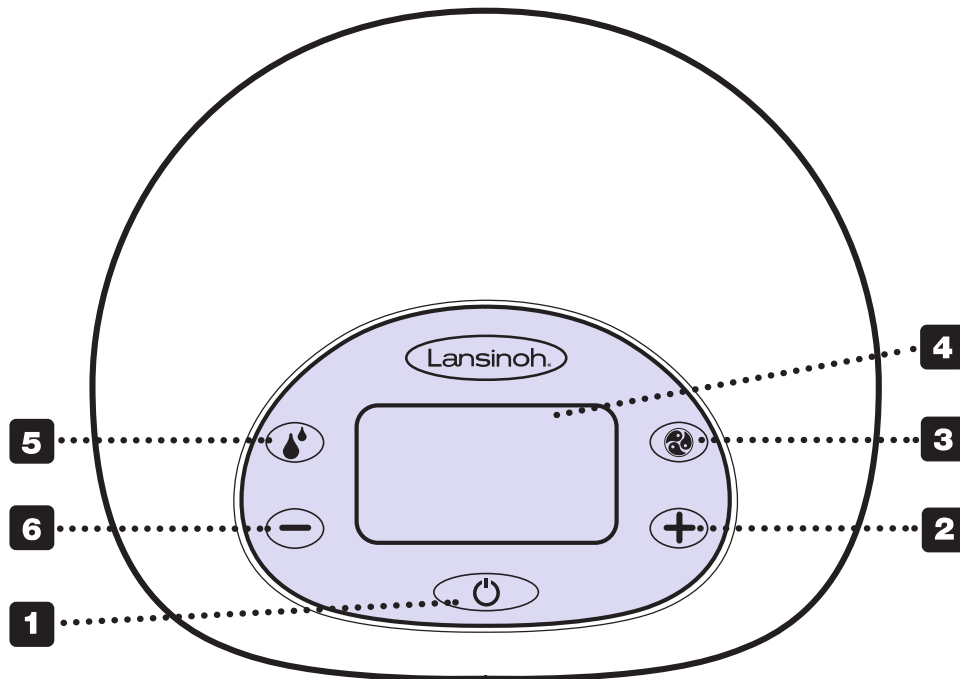
Adjusting Suction Levels

Increase or decrease the suction to the maximum level that you find comfortable. Increase suction by pressing the  button. Decrease suction by pressing the  button.

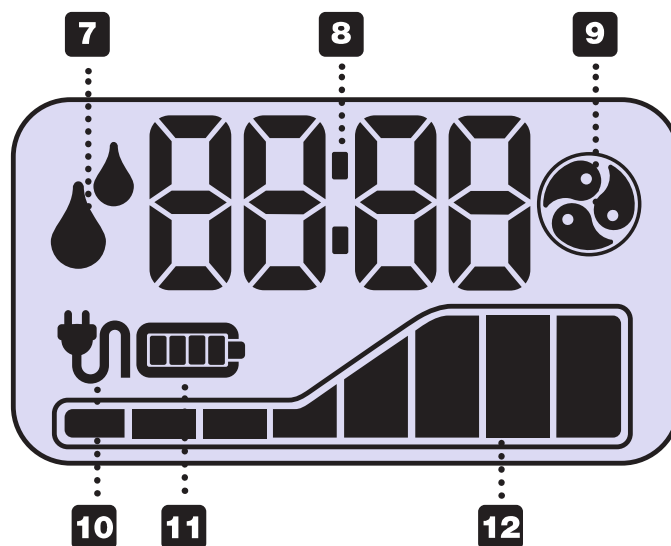
Lansinoh is here to support you. If you have questions, please do not hesitate to **contact us at (800) 292-4794**


PUMP DIAGRAM


Basic Controls, Display and Features

- | | |
|--|--|
| 1 Power, On/Off button  | 7 Let-Down/Expression Phase Indicator |
| 2 Suction Level Increase Control Button  | 8 Time Display |
| 3 Pumping Style Toggle Button  | 9 Pumping Style Indicator |
| 4 LCD Display | 10 AC Power Indicator |
| 5 Let-Down/Expression Toggle Button  | 11 Battery Power Indicator |
| 6 Suction Level Decrease Control Button  | 12 Suction Level Indicator |

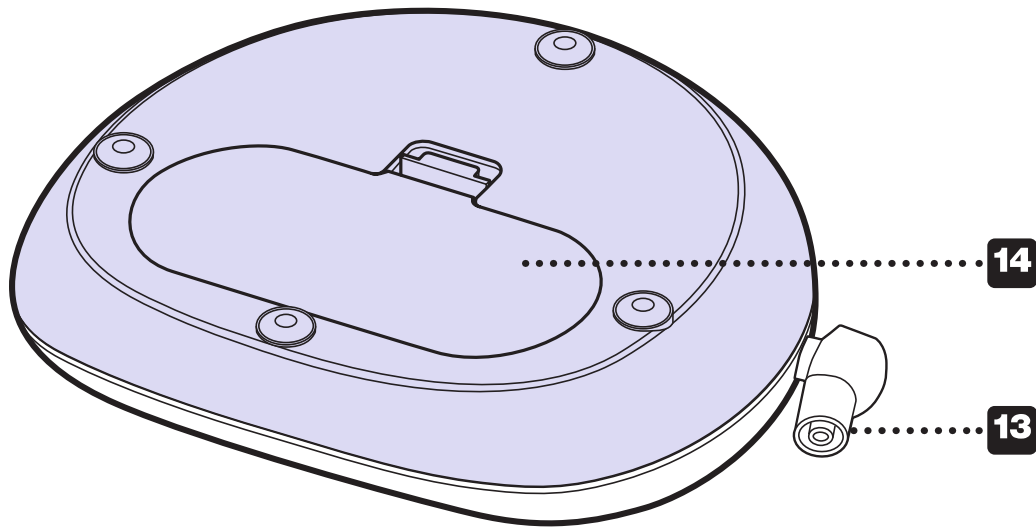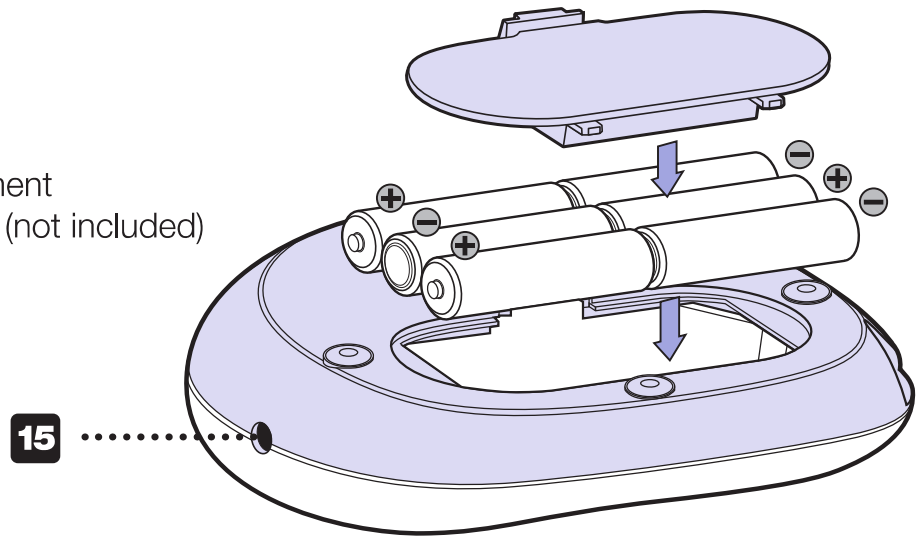
- 13** Tubing input
- 14** Battery compartment for 6 AA batteries (not included)
- 15** AC adapter input


AA Battery Compartment
(underneath pump)

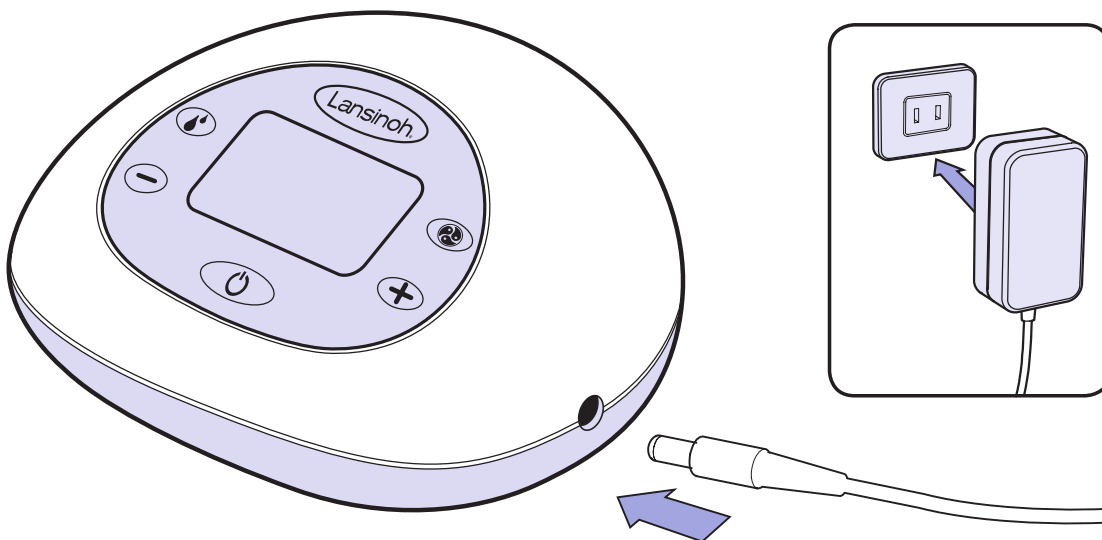AC ADAPTER

Lansinoh® Signature Pro™ Breast Pump Parts

Item	Item Number
a. Signature Pro™ Tubing	#53403
a1 Silicone Tubes 300 mm (2)	
a2 Silicone Tube 600 mm	
a3 Signature Pro™ Tubing Y Connector	#53405
a4 Signature Pro™ Tubing Strap	#50417
a5 Signature Pro™ Pump Connector	#53453
b. Breastmilk Expression Set Includes:	
b1 ComfortFit™ Breast Flange 25mm (2)	#50408
b2 Signature Pro™ Breast Flange Body (2)	#53454
b3 Signature Pro™ Diaphragm (2)	#53402
b4 Signature Pro™ Diaphragm Cap (2)	#53410
b5 White Valve (2 + 2 extra)	#53401
b6 5 oz (150 ml) Bottle (2)	#53451
c. Storage Lids	#53450
d. Signature Pro™ AC Adapter/power cord	#50409

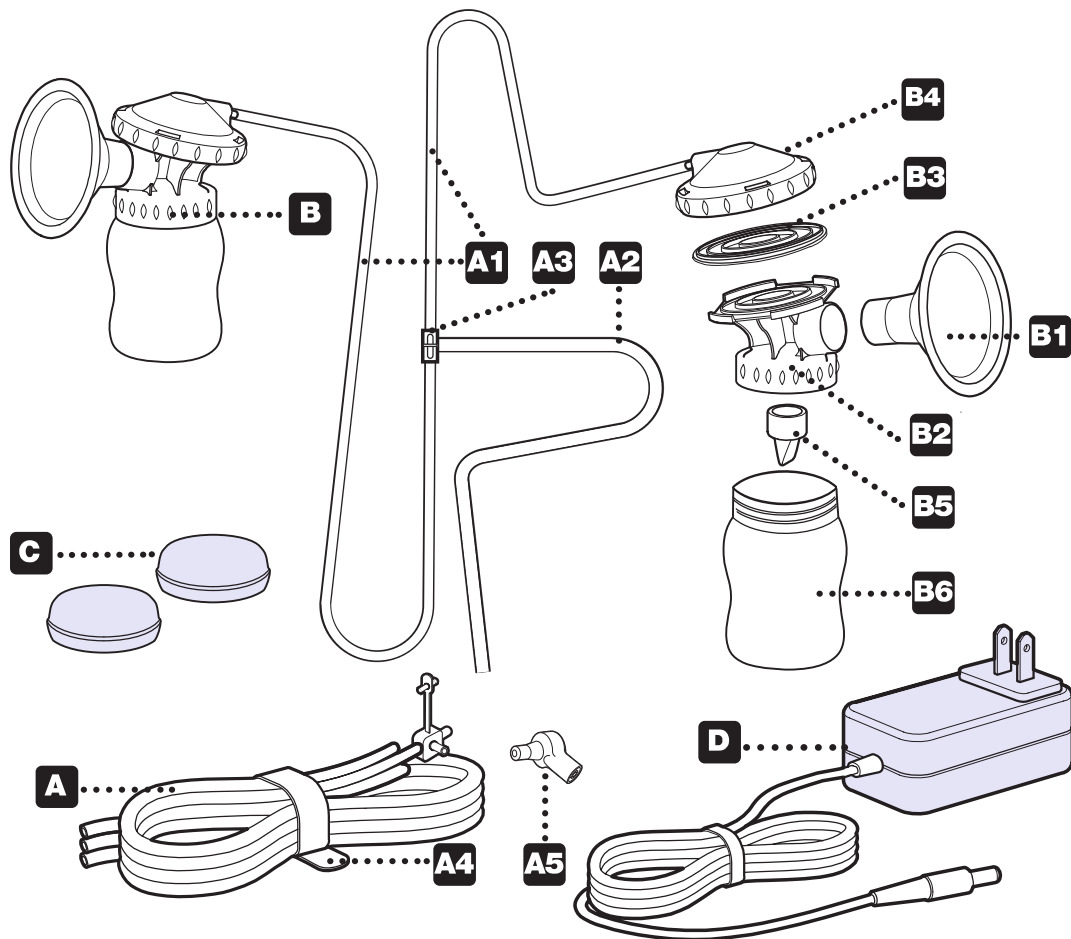Two extra white valves are included. Larger size flanges are sold separately.
 To order replacement parts or larger size flanges
 please visit www.lansinoh.com/store.

3. CLEANING INSTRUCTIONS

Unless otherwise instructed by your healthcare provider, follow these sanitizing and cleaning instructions.

Do not boil or wash tubing or connector!

3.1 Before First Use

1. Always wash and dry hands thoroughly with soap and water before handling the breastmilk expression set.
2. Disassemble the breastmilk expression set and place parts in a large pot.
3. Cover the parts in water making sure there is room for all of the parts to float freely in the pot.
4. Bring the water to a boil and allow the parts to boil for 5 minutes.
5. After boiling remove the parts from the pot with tongs and place on a clean towel to dry, covering the parts with a second towel.
6. Allow the parts to dry completely before assembly and use. Partially dried parts may affect the operation of the pump.
7. Once parts are dry they can be stored in a clean plastic bin with a lid or in a clean plastic zipper bag.

When sanitizing parts, use enough water to keep parts afloat!

3.2 Cleaning and Sanitization after Each Use

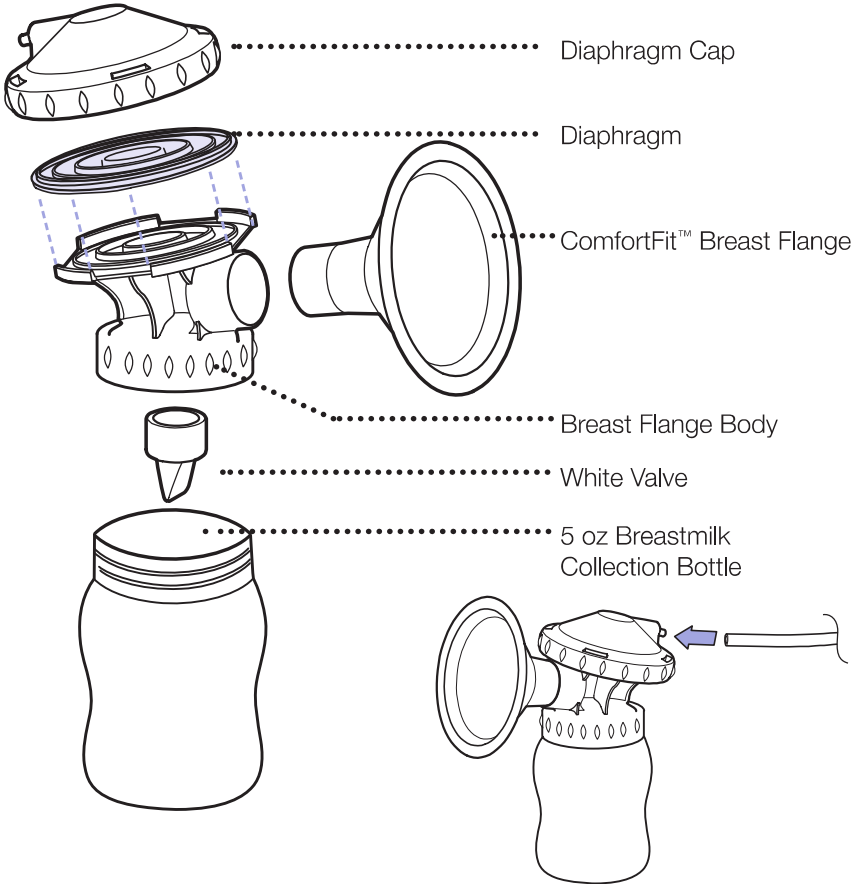
1. Disassemble the breastmilk expression set, setting aside the tubing and connectors.
2. Rinse all parts with cold water immediately after each use.
3. Hand wash all parts of the breastmilk expression set using warm soapy water. A nipple brush may be used to scrub the hard plastic pieces of the breastmilk expression set, but should not be used on the soft silicone valve or diaphragm.
4. Rinse parts with clean, warm water.
5. Place parts in a pot and cover with water making sure there is room for all of the parts to float freely in the pot.
6. Bring the water to a boil and allow the parts to boil for 5 minutes.
7. After boiling remove the parts from the pot with tongs and place on a clean towel to dry, covering the parts with a second clean towel.
8. Allow the parts to dry completely before assembly and use. Partially dried parts may affect the operation of the pump.
9. Once parts are dry they can be stored in a clean plastic bin with a lid or in a clean plastic zipper bag.

3.3 Before Each Use

1. Always wash and dry hands thoroughly with soap and water before handling the breastmilk expression set.
2. Before assembling the breastmilk expression set, be sure that all of the parts are completely dry.
3. Visually inspect the individual components for cracks, chips, tears, discoloration or deterioration. These may affect the functionality of the pump.
4. Should you have any questions or concerns about your parts please contact Lansinoh at **800-292-4794** for further guidance.
5. To purchase replacement parts for your pump please visit **Lansinoh.com/Store**

4. ASSEMBLY


4.1 Assembling the Breastmilk Expression Set

1. Place purple silicone diaphragm, arrow side up, on top of breast flange body. Ensure diaphragm fits securely between rims of breast flange body and arrows point toward breast flange.

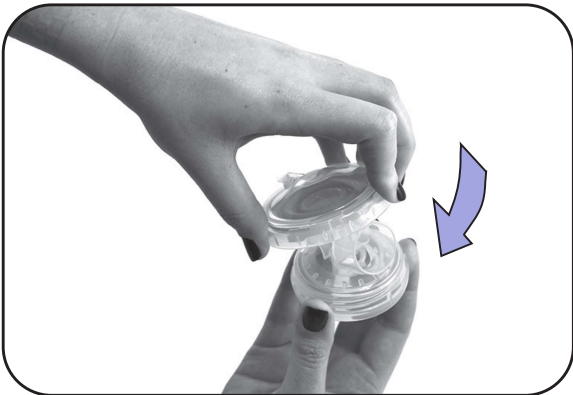
2. Place clear diaphragm cap over breast flange body and twist to secure so that the tubing portal is facing the back of the pump.

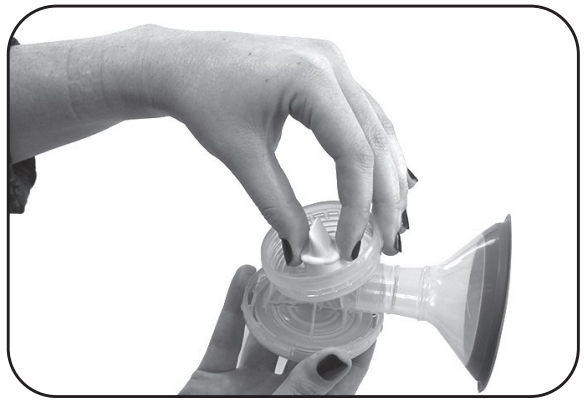3. Firmly push ComfortFit™ Flange into breast flange body.

4. Holding the base of the white valve, gently secure the valve into the bottom of the breast flange body.


NOTE: Use care when handling and cleaning the small white valve.

5. Screw collection bottle onto breast flange body.

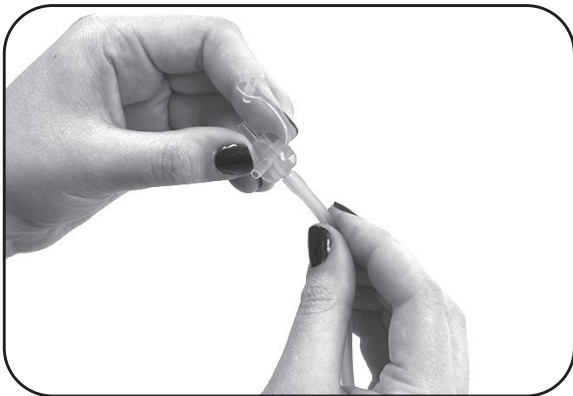6. Attach end of 300 mm tubing to back of diaphragm cap. Ensure that tubing is pushed all the way to the base of the attachment port and completely covers it.


7. Attach other end of 300 mm tubing to Tubing Y Connector.

8. Repeat assembly process for second breast flange body.

4.2 Connecting the Breastmilk Expression Set

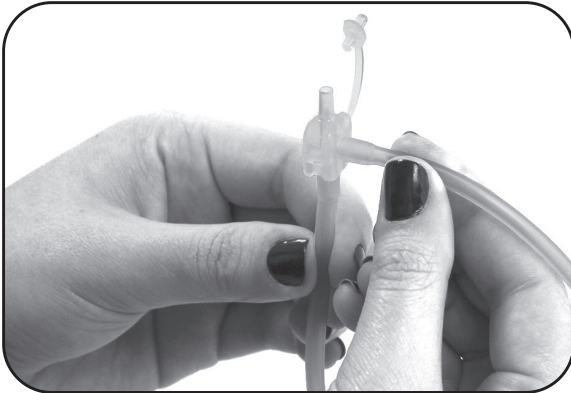Once Breastmilk Expression Set is assemble, connect to pump base:

1. Secure one end of 600 mm tubing to Tubing Y Connector.

2. Secure other end of 600 mm tubing firmly to the pump connector.

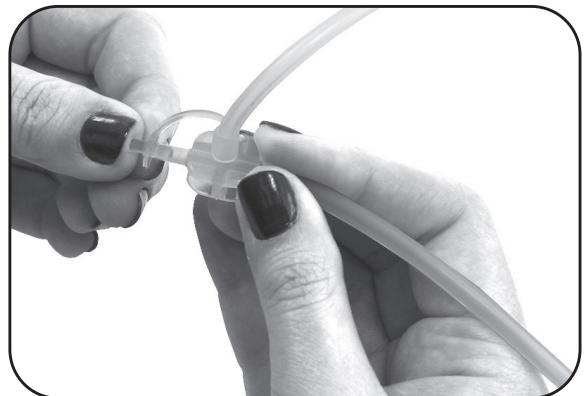
3. Insert pump connector into opening on left side of pump.

4.3 Connecting for Single Pumping

1. Remove one of the 300 mm tubing from the Tubing Y Connector.

2. Secure plug into open outlet of Tubing Y Connector.

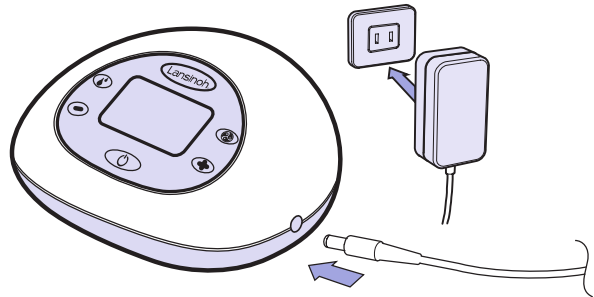
5. USING THE LANSINOH SIGNATURE PRO DOUBLE ELECTRIC BREAST PUMP

5.1 Power Options

The pump may be operated using the following power supply modes:


AC Adapter

1. Connect the AC adapter to the power receptacle located on the right side of the pump. Plug the AC adapter into an available electrical outlet.
2. Electrical safety can only be achieved when the breast pump is connected to the AC Power adapter supplied with the product. (See “Technical Data” section, page 21)

AA Alkaline Batteries (not included)

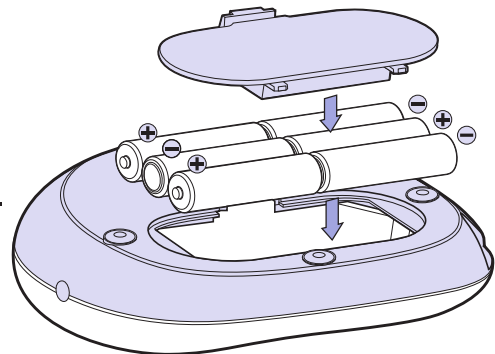1. Make sure the batteries are inserted correctly according to the polarity (+) and (-) as indicated.
2. Close the compartment by gently pressing the cover until it snaps into place.


- Remove batteries if pump is not in use for a long period of time.
- Dispose of batteries according to regulations for environmental protection and recycling.
- Replace batteries as a simultaneous set. Do not mix old and new batteries.
- Do not mix alkaline and standard (carbon-zinc) batteries.
- If batteries are in the pump, but AC adapter is in use, the pump will use the adapter to power the pump, not the batteries.

For the most economical use, we recommend using your pump with the AC adapter provided. Your AC adapter is a universal AC plug adapter type and can be used with 100 Vac to 240 Vac systems with 50Hz and 60Hz mains.

5.2 Getting Started

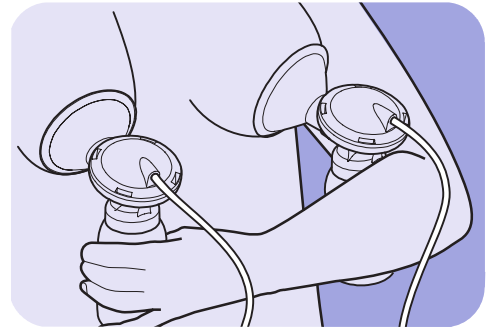If possible, wait until your breastmilk supply and breastfeeding schedule are well established. This is normally at least 3-4 weeks after birth, unless told otherwise by your doctor or breastfeeding advisor.

- 1.** Hold the Breastmilk Expression Set to your breasts so that the nipple is aligned to extend straight toward the nipple tunnel of the breast flange.
- 2.** Make sure the breast completely fills the breast flange so that a vacuum is formed and no air escapes.
- 3.** The Breastmilk Expression Set and nipple tunnel must tilt slightly downward to allow milk to flow naturally toward the bottle.

If you experience friction or discomfort on your nipples, you may need a different size flange. Larger size flanges are available. For more details, **please call 800-292-4794** or visit **www.lansinoh.com/store**.

5.3 2-Phase Technology

Using Let-Down and Expression Phases

- 1.** Press and hold the On/Off Button to turn the pump on. The Lansinoh® Signature Pro™ Double Electric Pump will start in the Let-Down Phase.
- 2.** After two minutes, the Let-Down Phase will automatically change to the Expression Phase in Cycle Mode 2.
- 3.** If you prefer to bypass the pre-set two minute Let-Down Phase or your milk lets-down sooner than two minutes, press the Let-Down/Expression Toggle Button to move to the Expression Phase.

The Style Mode Indicator will appear in the display to signify that you are in the Expression Phase.

5.4 Customization for Comfort and Efficiency

Using Customizable Pumping Styles

Choose your style

Each baby is unique – right down to the ways they feed. The Lansinoh® Signature Pro™ offers three different pumping styles to choose from – so you can select the style that best simulates your baby’s feeding style.

	Pumping Style 1	
	Pumping Style 2	
	Pumping Style 3	

- 1.** With the pump in Expression mode, select the style that is most comfortable and efficient for you by pressing the Pumping Style Toggle Button .
- 2.** Pressing the button again will select the next style.
- 3.** The mode will be indicated on the LCD Display by the Style Mode Indicator.

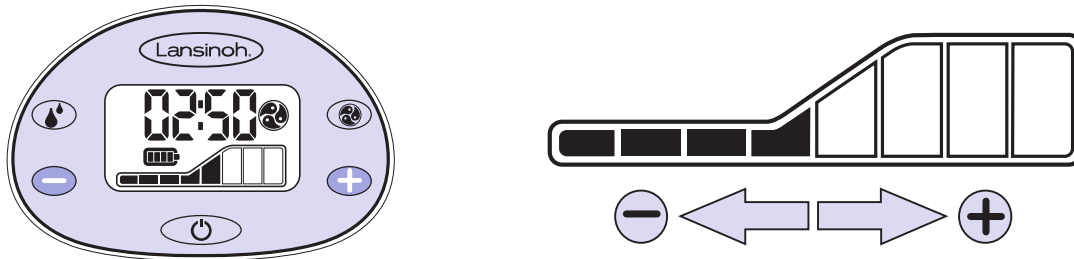
You may find that the best style for you changes over time. The customization features of the Signature Pro™ allow you to adjust with you and your baby.

5.4 Customization for Comfort and Efficiency (cont.)

Tailor for comfort

Within each pumping style, there are eight levels of suction to choose from:

1. Adjust for higher suction by pressing the  suction level increase control button.
2. Adjust for lower suction by pressing the  suction level decrease control button.
3. The current suction level will be indicated by the number of dark bars on suction level indicator.

NOTE: If at any time the suction on your breast reaches uncomfortable levels, reduce the suction level by pressing the  button on the left side of the LCD screen or disconnect the tubing from its receptacle on the electric pump.

5.5 Completing Your Pumping Session

1. When a pumping session is finished, turn the pump off with the on/off  button.
2. After turning the pump off, the display will remain on for a few seconds.
3. Make sure there is no suction on the nipple(s) before removing the breast flanges from the breast. If necessary, break the suction by inserting a finger between the breast flange and the breast. Once the suction is broken, the breast flange can be removed safely.

PLEASE NOTE: The LCD display light will automatically switch off if no control buttons are pressed for 20 seconds in order to save power. If the pump unit is left on and unattended for more than 60 minutes, the control system will automatically power off. To turn unit back on, hold the power-on button for several seconds.

6. PUMPING TIPS

Pumping Tips

The most important part of successful pumping is the ability to trigger the let-down reflex. This is when your milk begins to flow freely. Warm wet compresses applied to the breasts and gentle, circular massaging of the breasts prior to and/or during a pumping session can help encourage this natural reflex. Thinking of or looking at your baby, or a photo of your baby if you are apart, can also be helpful. Some moms can achieve more than one let-down during a pumping session. You may be able to achieve a second by returning to the Let-Down phase by pressing the Let-Down/Expression Toggle Button.

How Long To Pump

Double pumping, compared to single pumping, will reduce the time you spend expressing your milk. Breastmilk is produced on a supply and demand basis. Expressing your breastmilk regularly will stimulate your body and help you maintain your milk supply. Just like with breastfeeding, it is not necessary to time your pumping sessions. As you get more experienced with pumping, you will get a better sense of the time it takes for your body to express your milk using a breast pump. If you have any breastfeeding or breastpump related questions, consult with an International Board Certified Lactation Consultant (IBCLC), breastfeeding counselor, or other healthcare provider such as your primary physician.

Pumping and Feeding Tips

- You may pump one breast while nursing the baby on the other. This allows you to maximize pumping effectiveness using the let-down reflex stimulated by your baby's sucking.
- Always test the temperature of any liquid before feeding your baby.
- Do not let a child take a bottle to bed or self-feed for long periods.
- To prevent possible choking hazard, test strength of the nipple by pulling on the bulb portion. Discontinue use if any tear or crack appears.

7. TROUBLESHOOTING

Troubleshooting the Lansinoh® Signature Pro™ Double Electric Breast Pump

Interruption of Electrical Service

If you experience an interruption in the normal operation of the electrical and/or mechanical functioning of the pump, simply do the following:

- Unplug the AC adapter cord from the electrical outlet.
- Wait a few seconds, and then plug the AC adapter cord back into the electrical outlet. If you still experience electrical failure, try using the pump with the six AA batteries making sure the AC adapter is removed.

Discomfort While Pumping

If you feel any discomfort, you may need a different size flange.

Visit www.lansinoh.com/store or call 800-292-4794.

Decreased Suction

If it feels like the suction has decreased, try the following:

- Remove the tubing from pump and place a finger over the tubing connector.
 - If you feel suction, the pump is operating properly, but the Breastmilk Expression Set may not be correctly assembled or the white valve may be damaged.
 - Make sure the Breastmilk Expression Set has been properly assembled. Check the assembly information in the section “Assembling the Breastmilk Expression Set,” page 10-12.
 - Inspect the white valve for a tear or pinhole. This component is critical to achieving proper suction. Replace the white valve with one of the spare valves provided with the product.
 - If no suction is felt on your finger, try unplugging the adapter from the electrical outlet, wait a few seconds, and plug it back into the electrical outlet.
- When single pumping, check the tubing connector on the front of the pump to make sure it is closed for single pumping. If open, close the connector.
- If double pumping, ensure that the Breastmilk Expression Set is on both breasts before starting the pump.

Increased suction level for single pumping:

If it feels like the suction is higher than usual when single pumping, check that the plug for the pipe is clean and that the airway in the plug is not blocked.

8. BREASTMILK STORAGE GUIDELINES

Breastmilk Storage and Thawing Guidelines For Healthy Term Babies

Check with your healthcare provider or an International Board Certified Lactation Consultant for specific storage instructions. When freezing, store in the back of the freezer and do not store milk in the door as temperature can fluctuate as the door is opened and closed.

WHERE	TEMPERATURE	TIME
At Room Temperature	66-78°F (19-26°C)	4 hours (ideal) up to 6 hours (acceptable)
In a Refrigerator	<39°F (<4°C)	72 hours (ideal) up to 8 days (acceptable)
Freezer	0°F (-18°C)	3 months (ideal) up to 6 months (acceptable)
Deep Freezer	-4°F (-20°C)	6 months (ideal) up to 12 months (acceptable)

Source: La Leche League International 2014

Date all breastmilk and use the oldest breastmilk first.

Thawing

Thaw frozen milk in the refrigerator or place the containers under cool running water, gradually warming the water until the milk is thawed. Thawed milk can be kept in a refrigerator for up to 24 hours.

NOTE: DO NOT refreeze milk once it is thawed.

Gently shake thawed milk before feeding to baby to mix the layers that have separated. Breastmilk is not homogenized and the “cream” will rise to the top of the bottle. Breastmilk can vary in color and consistency depending upon the time of the day that the milk was expressed, what you have eaten, and the age of the baby at the time of pumping.

NOTE: DO NOT MICROWAVE HUMAN MILK! Microwaving alters the composition of the milk. Microwaving also may create “hot spots” in the milk that could burn the infant.

NOTE: DO NOT thaw breastmilk in hot or boiling water.

9. THE LANSINOH® NATURALWAVE® NIPPLE

Clinically proven to reduce nipple confusion in breastfed babies

The patented NaturalWave® Nipple enables baby to replicate natural sucking actions learned at the breast.

Exclusive breastfeeding is recommended for at least 3-4 weeks if possible before introducing a bottle and nipple.

Storage & Use:

Prior to use, always inspect and pull the feeding nipple in all directions to check the product. Throw away at the first signs of damage or weakness. For safety and hygiene reasons, replace the nipple every 7 weeks.

Before first use, boil all parts for 5 minutes to sterilize. Ensure all parts are completely covered in water and free floating.

Cleaning:

Before first use, place all components in boiling water for 5 minutes, allow to cool and thoroughly rinse all bottle components. Always wash all components of the product immediately after each use. The feeding nipple can be washed in hot soapy water and rinsed in clean water.

Dishwasher safe on top rack only. Safe to clean in a

microwave-safe, steam sterilizer bag and special sterilizing machine. Do not clean with solvents or harsh chemicals.

Assembly:

Once thoroughly dry, insert feeding nipple through the top of the screw cap collar. Attach the screw cap collar securely onto the bottle neck. When carrying a full bottle, make sure it is closed with the collar in place, taking care to position the feeding nipple in the center of the collar to activate the “spill-proof” function of the product.

Precautions:

Do not leave bottle in direct sunlight, heat or disinfectant (sterilizing solution), as this may weaken the nipple. Do not clean with or allow feeding nipple to come into contact with solvents or harsh chemicals. Damage could result. Fill with the required amount of liquid and tighten the screw cap before using. Do NOT overfill. Heating in a microwave is NOT recommended as it may produce uneven heating and high localized temperatures. If you choose to heat liquid this way, DO NOT SEAL BOTTLES. Swirl well before use, and always check the temperature before feeding baby. Hot or unevenly heated liquid may burn baby's mouth. Do not warm liquid by placing the bottle directly onto gas or electric rings, warming plates, in electric ovens, etc. Wash after every use. Some liquids and sterile solutions may make the silicone become opaque. This does not alter any of its properties.

WARNING!

DO NOT MICROWAVE HUMAN MILK! Microwaving alters the composition of the milk and may also create “hotspots” in the milk that could burn the infant. Allowing your child to use this product for prolonged periods separate from regular mealtimes or to go to sleep while drinking liquids, other than water, from this bottle may cause or contribute to early childhood tooth decay. Always use this product with adult supervision. Keep all components not in use out of the reach of children. Do not leave baby alone with drinking equipment due to the risk of choking. Never use feeding nipple as a pacifier or soother.

Always check liquid temperatures before feeding. Packaging materials are not part of the product. For your child's safety, please remove and discard all packaging materials prior to use but read and retain the instructions for future reference.

Warning:

Packaging materials are not part of the product. For your child's safety, please remove and discard all packaging materials prior to use but read and retain the instructions for future reference.

Latching on

Baby's lips open outward and latch onto the areola

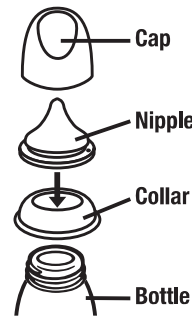
Peristaltic “wave-like” tongue movement

Baby's tongue moves in “wave-like” motion to compress nipple and extract milk

Swallowing

Back of baby's tongue rises, channeling the milk into the esophagus

Insert the nipple into top of the collar

10. TECHNICAL DATA

1. Power Supply

Battery operation: 6 AA size alkaline batteries, all connected in series

Battery operating time: Approximately 2 hours with 6 alkaline batteries

AC Adapter operation: 110Vac / 220Vac wall plug type AC adapter, galvanically isolated with 9Vdc output voltage

Pump power dissipation: 10.0 Watts maximum

Voltage regulation: Linear voltage regulator with low dropout voltage and short circuit protection

2. Control System

Control unit: Microcontroller with integral program memory

Program memory: 32K flash ROM, 2K RAM memory

HMI: LCD display with LED backlight, soft keypad

Single / double control: Single or double breast pumping possible

Selectable pumping rhythm mode: 3 pumping rhythm modes are selectable and controlled from keypad

Adjustable suction level: 8 suction levels are selectable and controlled from keypad
Suction vacuum level: Minimum 80 mmHg, maximum 220 mmHg

Technical Data

Power Requirements

AC Adapter (included)

AC Adapter – North American

Plug Configuration

110Vac/220Vac wall plug type

AC adapter, galvanically isolated

adapter with 9Vdc output voltage.

Batteries (not included)

6 AA batteries, 1.5V alkaline

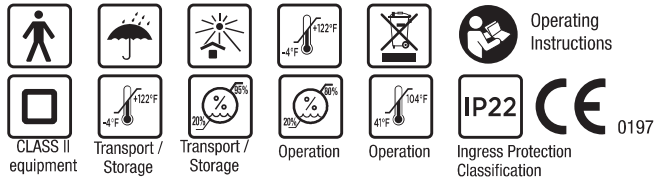
Regulatory Compliance

Bridge Power Adapter: UL (E302267) and CSA approved with North American power adapter.


Shenz Adapter: UL (E314167)

Regulatory Compliance

Type BF Applied Part

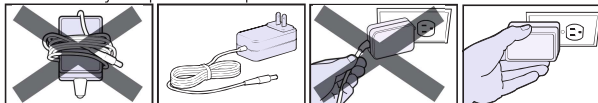To maintain your power adapter follow these instructions:

DO NOT wrap cord of power adapter around the plug body.

CORRECT

DO NOT unplug power adapter by pulling on the cord.

CORRECT

DO NOT use power adapter if damaged.

Warning: Children have STRANGLED in cords. Keep the cords out of the reach of children.

11. WARRANTY

Limited Consumer Product Warranty (United States)

This product is warranted to the original consumer purchaser (the “purchaser”) to be free from defects in material and workmanship which are not commercially acceptable for the period of one year from the date of purchase on the pump mechanism and the period of 90 days on all other product components. Warranty coverage terminates if you sell or otherwise transfer this product to another person. This warranty gives you specific legal rights and you may also have other rights, which vary, by location.

LANSINOH® MAKES NO EXPRESS WARRANTY OF ANY KIND REGARDING THE PRODUCT OTHER THAN THOSE WARRANTIES EXPLICITLY SET FORTH HEREIN. ANY IMPLIED WARRANTY, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR ANY IMPLIED WARRANTY OF FITNESS FOR A PARTICULAR PURPOSE, TO THE EXTENT PERMITTED BY LAW, SHALL BE LIMITED IN DURATION TO A PERIOD OF 90 DAYS FROM THE DATE OF PURCHASE BY THE ORIGINAL PURCHASER.

In the event the product is found by Lansinoh® to not meet the above limited warranty, as purchaser’s sole and exclusive remedy Lansinoh® will repair or, at the option of Lansinoh®, replace this product without charge for such replacement or parts or labor.

NOTE: The purchaser shall bear all expense for returning this product to Lansinoh®. This warranty does not apply to any product that has been subjected to misuse, abuse, or alteration.

Only Lansinoh® brand replacement parts are recommended for use with Lansinoh® breastfeeding products. Improper or incorrectly performed maintenance or repair, or use of non-Lansinoh® brand replacement parts may invalidate this warranty.

WARRANTY

TO THE EXTENT PERMITTED BY LAW, LANSINOH® SHALL NOT BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES INCLUDING, BUT NOT LIMITED TO, REPLACEMENT COSTS RESULTING FROM THE BREACH OF ANY WRITTEN OR IMPLIED WARRANTY.

If you wish to make a claim under this warranty, you must return the product to Lansinoh Laboratories Inc., prepaid, together with your dated bill of sale or other dated proof of purchase and a brief statement of the problem, along with a Replacement Authorization (RA) number. You may obtain the appropriate Replacement Authorization number and shipping address by calling the phone number or sending an email to the email address provided below.

Items returned without an associated Replacement Authorization Number will be returned to sender.

To obtain the address to Return Goods or a Replacement Authorization Number:

800-292-4794
customerservice@lansinoh.com

Submit warranty card online at **Lansinoh.com/warranty**

Declaration for Electromagnetic Compatibility (EMC)

For information use only

FCC Statement:

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Wireless Statement:

Wireless communications equipment such as wireless home network devices, mobile phones, cordless telephones and their base stations, walkie-talkies can affect this equipment and should be kept at least a distance 3.5m away from the equipment.

Stacking:

This device should not be used adjacent to or stacked with other equipment. And if adjacent or stacked use is necessary, this device should be observed to verify normal operation in the configuration in which it will be used.

 Manufacturer: Lansinoh Laboratories Saglik Gerecleri Tasarim San. Tic. Ltd. Sti.
Zafer SB Mah. Gunduz Sk. No:13 E.S.B Gaziemir/ Izmir/ Turkiye

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This product is a personal care product and therefore cannot be returned once opened. If you experience any problems with this product, please telephone our toll free customer service line at 800-292-4794. Our representative will be happy to assist you.

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Lansinoh® is a registered trademark of Lansinoh Laboratories, Inc.
© 2020 Lansinoh Laboratories, Inc. All rights reserved. Printed in Turkey.
DE53050IF0120

Distribuido por Lansinoh Laboratories, Inc
Alexandria, VA 22314 USA
800-292-4794
Lansinoh.com

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Fabricante: Lansinoh Laboratories Saglik Gerecleri Tasarim San. Tic. Ltd. Sti.
Distributed by Lansinoh Laboratories, Inc
Lansinoh® es una marca registrada de Lansinoh Laboratories, Inc.
© 2020 Lansinoh Laboratories, Inc. Todos los derechos reservados. Impreso en Turquía.