



Adhesives & Removers


Urofoam® Adhesive Foam Strips

“Conforms easily to the irregular contours of the body”

Urofoam® double or single-sided adhesive foam strips are soft, pliable foam strips specifically designed for the unique bonding applications of male external catheters to the skin and are suitable for use with both standard and self-adhering type catheters.

They provide a comfortable, air-tight seal that is an effective way to secure a male external catheter to the body. They are easy to apply and remove so when it comes time to change that catheter, it can be changed easily and quickly.

It is recommended that a male external catheter be changed once a day to promote healthy skin, however, with the use of Urofoam® strips, a male external catheter can be held in place for up to three days—results may vary, but it's nice to know that when you need to depend on an adhesive strip, you can depend on Urofoam®.

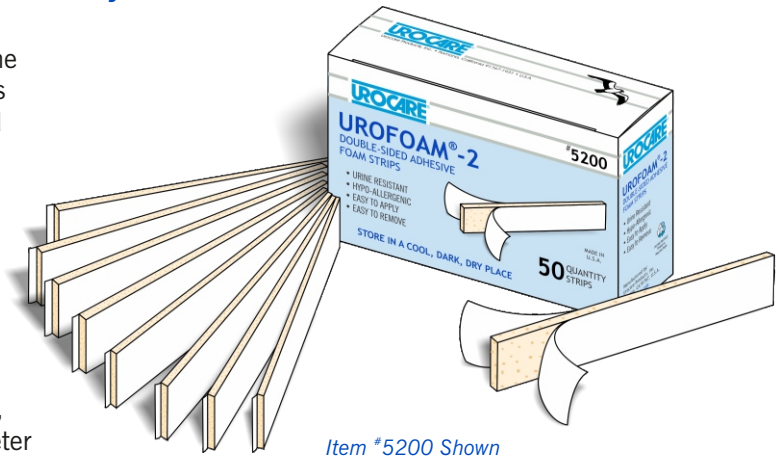

Item #5200 Shown

LEARN MORE



KEY FEATURES:

- **Ensolite Foam**—*Hypo-Allergenic.*
- **Closed-Cell Foam Structure**—*for an air-tight seal.*
- **Soft and Pliable**—*Easily conforms to irregular contours of the body.*
- **Perspiration Resistant**—*It retains its adhesive qualities even when sweating.*
- **Easy to Apply and Remove**—*Just peel 'n' Stick to apply, peel back to remove.*

If the strip is difficult to remove, Urocare® Adhesive Remover Pads (#5600) may be used to soften the adhesive and assist in the removal of the strip. They contain a **non-toxic, citrus-based** adhesive remover that may also be used to remove most types of medical adhesive or adhesive residue from skin and appliances or in preparation for use of a new strip.

Most medical grade adhesive removers will work basically in the same manner; however, consult your doctor to decide the best product for you and how to use it safely. For more information—see Urocare® Adhesive Remover Pads.

In the event of skin that has become sensitive to aggressive medical adhesives, Urocare® has developed Uro-Prep® Protective Skin Barrier Wipes (#5500) which may be used under all types of medical adhesives. They protect delicate, sensitive skin from aggressive adhesives often used in conjunction with urinary or ostomy appliances. Uro-Prep® wipes are safe, effective, and for external use only. For more information see Uro-Prep® Protective Barrier Wipes.

Sold by the box, 50 strips per box.
These are stocked items.
Non-Sterile, Latex-Free.

LATEX-FREE



Item No.	Description	Size†	Case Item No.
5200	Double-Sided	0.09" x 1" x 5.75" (0.3cm x 2.5cm x 14.6cm)	Not Available
5100	Single-Sided	0.09" x 0.88" x 5.75" (0.3cm x 2.2cm x 14.6cm)	Not Available

† Approximate dimensions in inches(") and centimeters(cm).

FIGURE A: Urofoam®-2 Method
(Foam strip under catheter)

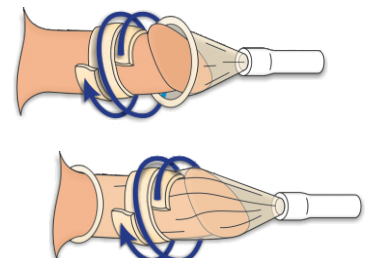


FIGURE B: Urofoam®-1 Method
(Foam strip over catheter)

