



Dremiumressure Check MATTRESSES

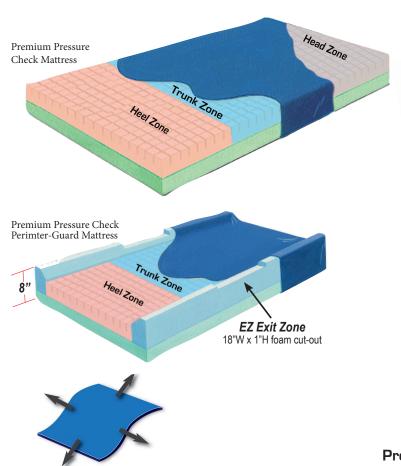

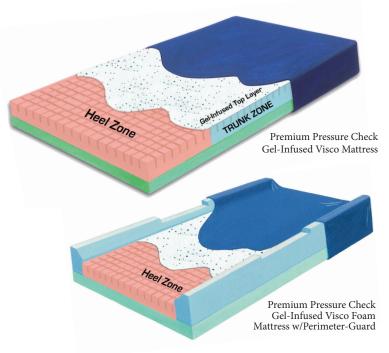

Pressure-Check Foam Mattress

Description	Size	Unit
Pressure-Check Mattress w/LSII Cover	76" x 35" x 6"	Each
Pressure-Check Mattress w/LSII Cover	80" x 35" x 6"	Each
Pressure-Check Mattress w/LSII Cover	84" x 35" x 6"	Each
Pressure-Check Mattress w/Perimeter-Guard w/LSIII Cover	76" x 35" x 6"	Each
Pressure-Check Mattress w/Perimeter-Guard w/LSIII Cover	80" x 35" x 6"	Each
Pressure-Check Mattress w/Perimeter-Guard w/LSIII Cover	84" x 35" x 6"	Each
	Pressure-Check Mattress w/LSII Cover Pressure-Check Mattress w/LSII Cover Pressure-Check Mattress w/LSII Cover Pressure-Check Mattress w/Perimeter-Guard w/LSIII Cover Pressure-Check Mattress w/Perimeter-Guard w/LSIII Cover	Pressure-Check Mattress w/LSII Cover 76" x 35" x 6" Pressure-Check Mattress w/LSII Cover 80" x 35" x 6" Pressure-Check Mattress w/LSII Cover 84" x 35" x 6" Pressure-Check Mattress w/Perimeter-Guard w/LSIII Cover 76" x 35" x 6" Pressure-Check Mattress w/Perimeter-Guard w/LSIII Cover 80" x 35" x 6"

Pressure-Check Gel-Infused Visco Foam Mattress

558123	Gel-Infused Visco Foam Mattress	76" x 35" x 6"	Each
558125	Gel-Infused Visco Foam Mattress	80" x 35" x 6"	Each
558127	Gel-Infused Visco Foam Mattress	84" x 35" x 6"	Each
558133	Gel-Infused Visco Foam Mattress w/Perimeter-Guard w/LSIII Cover	76" x 35" x 6"	Each
558135	Gel-Infused Visco Foam Mattress w/Perimeter-Guard w/LSIII Cover	80" x 35" x 6"	Each
558137	Gel-Infused Visco Foam Mattress w/Perimeter-Guard w/LSIII Cover	84" x 35" x 6"	Each

Pressure-Check™ Foam

Pressure-CheckTM Foam has 400 support pods that adjust individually for optimal weight distribution and effective pressure relief. The 400 individual "Support Pods" separate the mattress into three specific pressure redistribution zones; Head, Heel and Trunk.



LSIII Four-Way Stretch

High-Density Support Foam

The bottom layer provides uniform compression-resistant support for weight distribution, comfort and extended mattress life.

Important Fire Code Information

All Skil-Care mattresses comply fire code 16 CFR 1633.

Cover

The NEW Low-Shear III cover provides an industry top choice of 4-Way Stretch fabric that provides less moisture, more comfort, and cooler feel than a standard 2-Way stretch fabric. The waterfall zipper prevents liquid penetration.

Application

The 6" Mattress is specifically designed to provide an industry leading pressure redistribution and pressure relief environment for the resident. The "Support Pod" technology optimally addresses the key issues of weight distribution, pressure relief, anti-shearing, and provides a positive thermo-climate by eliminating heat buildup.

Cleaning Instructions

The cover can be wiped clean with mild detergents or anti-bacterial solutions. Do not use bleach or chlorine products. Do not laundry.

Caution

The mattress should not be flipped over, the bottom section of the mattress cover is made with slip-resistant fabric. The foam is encapsulated with a fiber-glass fire reduction covering and should not be removed.