


B[®] BIONIX[®]

GOING-BEYOND 

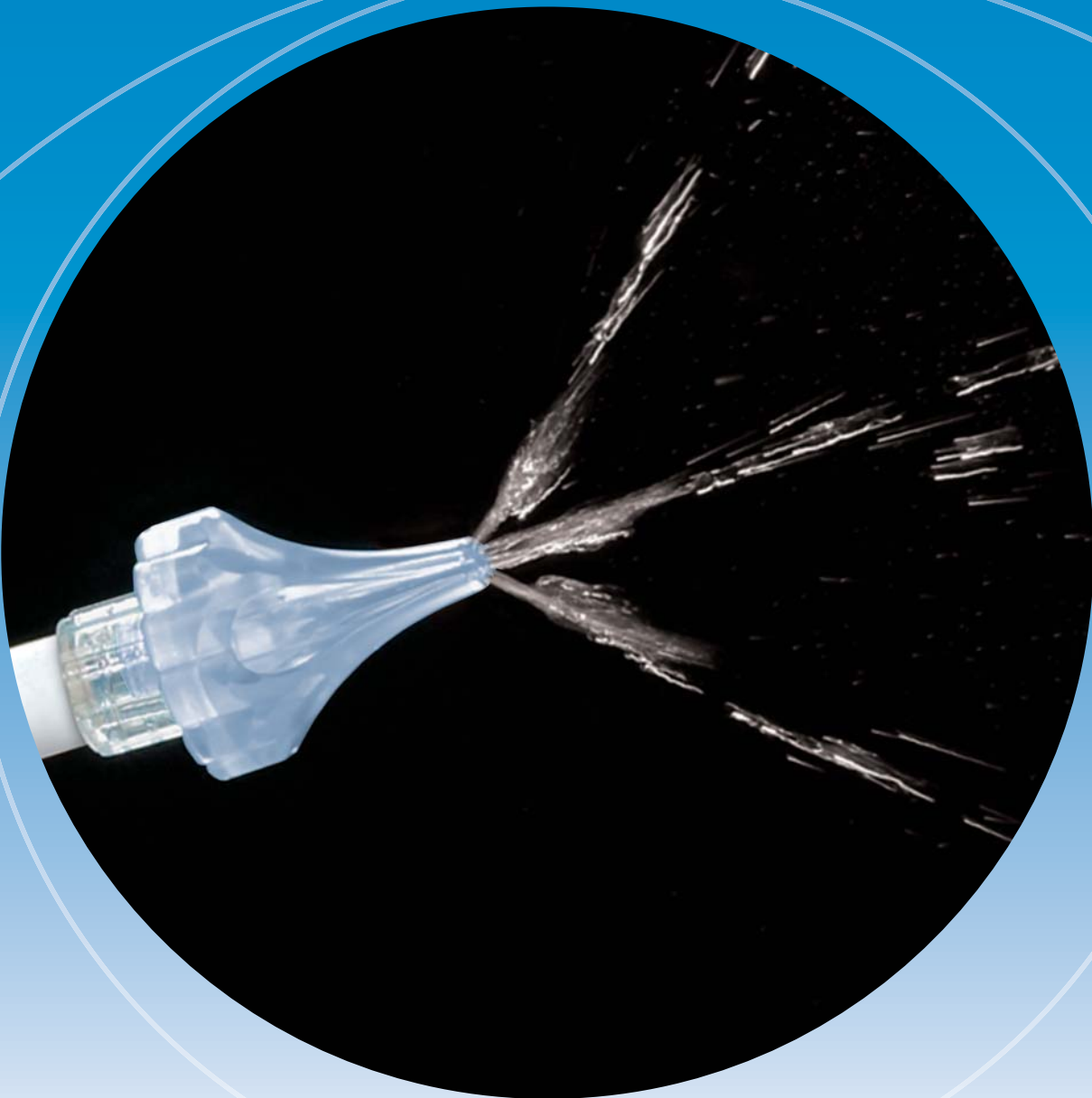

OtoClear[®] Ear Irrigation

Bionix OtoClear Ear Irrigation helps healthcare workers provide safe and comfortable ear irrigation without the backsplash, at half the cost of a self-constructed device.



“

“Our office is very happy with the soft (OtoClear) tips. They are easier to maneuver in the ear canal and more comfortable for the patient.”

”

- Laura Griffith, MA

OtoClear® Ear Irrigation with Gentle Touch™ Design

The Concept


Problem: Traditional ear irrigation shoots water directly into the ear canal, which can be painful, dangerous and messy.

Solution: OtoClear Ear Irrigation Tips, with Gentle Touch Design, directs water to the ear canal wall, with three divergent streams, instead of directly at the Tympanic Membrane.

Results: A safe, clean ear irrigation tip at half the cost of a self-constructed device.




How to Use

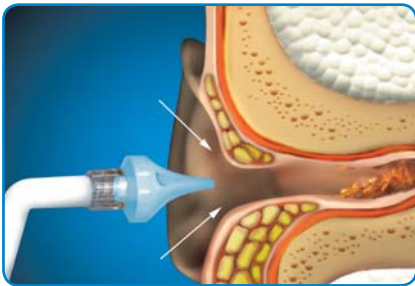

Insert *OtoClear* tip fully into the ear canal. The flared design of the *OtoClear* tip protects against over-insertion.



Begin irrigation. Note: The *OtoClear* tip directs the water towards the walls of the ear canal creating a turbulent lavage that effectively dislodges the cerumen. The exit portals control back splash and deliver the effluent into the basin.

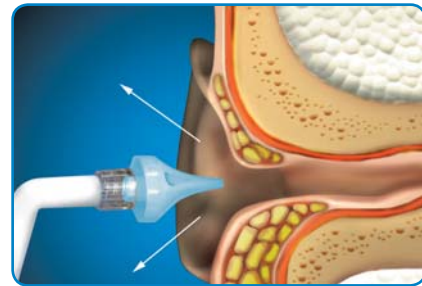

When procedure is complete, remove and dispose of the single-use *OtoClear* tip.

OtoClear Ear Lavage System Video



OtoClear Spray Wash Kit Video



Product Offerings



#7200
40 single-patient
tips per box



#7245
Portable Waterpik®



#3700
Ear Irrigation
Basins (box of 3)



#7280
Ear Lavage System

- 1 Portable Waterpik®
- 40 OtoClear Tips
- 1 Ear Basin



#7290
OtoClear Spray Wash Kit

- 1 Spray Wash Bottle
- 20 OtoClear Tips
- 1 Ear Basin

5154 Enterprise Blvd. • Toledo, Ohio 43612
Phone: 800.551.6810 • Fax: 800.455.5678
www.BionixMed.com

©2016 Bionix Development Corp. Rev. 1/27/16
*OtoClear® Tips - U.S. Patent #6,706,023. Waterpik® is a registered trademark of Water Pik, Inc.
The OtoClear® Ear Irrigation Tip is not affiliated with Water Pik Technologies, Inc.

B **BIONIX**[®]
GOING **BEYOND** 