

Taking performance to the max

AirLife™ brand Misty Max 10™ disposable nebulizer



AirLife brand Misty Max 10 disposable nebulizer

An advanced jet design for:

- Fast treatment times
- High percentage of respirable particles
- Low residual volume
- Consistent and reliable performance

Top 10 Features

1 Full 10 cc capacity provides medication dosing flexibility.



2 Unique, integrated (jet/baffle) design increases consistency and reliability in medication delivery.

3 Tinted snap-in jet remains in place yet allows for removal during cleaning.



4 Single twist closure of top and bottom gives confirmation of proper assembly.

5 Secure top and bottom seals prevent leakage.



6 Anti-spill design helps prevent medication waste.

7 Angled operation accommodates various patient positions.

8 Jet technology delivers efficient, high output rates, resulting in fast treatment times.

9 Clear base with visible graduations facilitates medication loading and visualization of aerosol delivery.

10 Easy to assemble rigid tubing connector and crush-resistant tubing.




Delivering fast and effective treatment, more efficiently. Accommodating all types of patients in all types of situations. These were the forces behind the development of the Misty Max 10 nebulizer. Featuring a new type of jet design, the Misty Max 10 nebulizer can speed up treatment and help ensure consistent medication delivery while offering you more flexibility in treatment options.

Nebulizer is available preassembled and fits virtually all common masks and adapters. All oxygen tubing is crush-resistant. Recommended for single patient use. The one-piece jet design and single-thread top connection make the unit easy to disassemble and reassemble.

Technical information (@ 8 L/min flow)*

MMAD	1.61 μm
GSD	2.18
Respirable fraction % < 4.7 μm (by mass)	85.2
Output rate	0.325 g/min
Treatment time	8 min, 9 sec
Residual volume	0.92 g
Fill capacity	10 mL

* Please see the Misty Max 10 nebulizer technical data sheet for additional information, including performance characteristics at 5 L/min and 10 L/min flow.

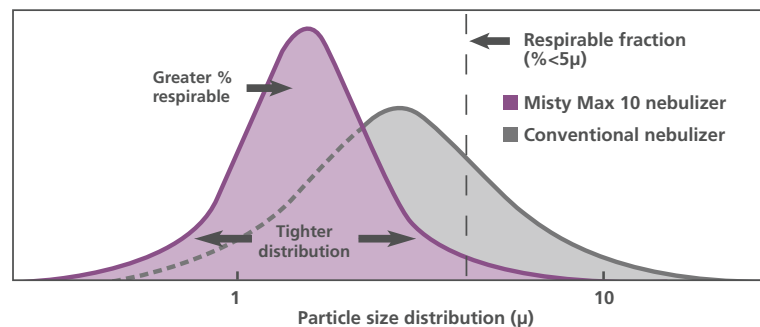

Particle size is an important characteristic of nebulizer performance. The Misty Max 10 nebulizer provides fine particles with a tighter distribution.

Ordering information


Cat. no.	Description
Nebulizers without mask, 50 per package	
002430	Nebulizer only
002431	With 7' oxygen tubing
002434	With baffled tee adapter, 7' oxygen tubing and mouthpiece
002435	With tee adapter, 7' oxygen tubing and mouthpiece
002437	With baffled tee adapter and mouthpiece
002438	With baffled tee adapter, 7' oxygen tubing, mouthpiece and 6" flextube
002439	With baffled tee adapter, 14' oxygen tubing, mouthpiece and 6" flextube
002442	With baffled tee adapter and intermittent finger control
002445	With baffled tee, 7' U/Connect-It® oxygen tubing and mouthpiece
002446	With baffled tee, 7' U/Connect-It oxygen tubing, mouthpiece and 6" flextube
002447	With 22 mm valved tee adapter (22 mm I.D. x 22 mm O.D.) and 7' oxygen tubing
Nebulizer without mask, 25 per package	
002450	With patient "Y" mouthpiece, 7' oxygen tubing, 6" flextube and bacteria filter
Nebulizers with adult mask, 50 per package	
002432	With adult aerosol mask
002433	With adult aerosol mask and 7' oxygen tubing
Nebulizers with pediatric mask, 50 per package	
002443	With pediatric aerosol mask
002444	With pediatric aerosol mask and 7' oxygen tubing
002440**	With pediatric aerosol elephant mask
002455**	With pediatric aerosol elephant mask, flextube adapter and 7' oxygen tubing
Masks only, 50 per case	
001271**	Infant turtle aerosol mask
001261	Pediatric aerosol mask, short style
001262	Pediatric aerosol mask, under-the-chin style
001276**	Pediatric elephant aerosol mask with flextube
001266**	Pediatric dragon aerosol mask
001206	Adult aerosol mask

Nebulizer (shown with dragon mask 001266—order separately)



For more information about the Misty Max 10 nebulizer or the complete line of *AirLife* products and services, please contact your Respiratory Care Products and Services sales representative, call customer service at **800.964.5227** or visit us at **carefusion.com/respiratory**.

 **WARNING**—U.S. Federal Law restricts this device to sale by or on the order of a physician.

 CareFusion
22745 Savi Ranch Parkway
Yorba Linda, CA 92887

800.231.2466 toll-free
714.283.2228 tel
714.283.8493 fax

 CareFusion France 205 S.A.S
5 rue Charles Lindbergh
F-44110 Châteaubriant, France

+33 240 813 804 tel
+33 0 240 815 634 fax


0086

CareFusion
Yorba Linda, CA

carefusion.com

**KidsMED, Inc.

© 2011 CareFusion Corporation or one of its subsidiaries. All rights reserved. AirLife and Misty Max 10 are trademarks or registered trademarks of CareFusion Corporation or one of its subsidiaries. KidsMED is a registered trademark of KidsMED Inc. L3374 Rev. A RC2405 (0211/1000)

