



M50 STATIONARY OXYGEN CONCENTRATOR



GCE[®]
MEDICAL

GCE[®]
Healthcare

Ohio Medical

M50 Oxygen Concentrator

The lightweight, ergonomically designed M50 oxygen concentrator makes life simpler for patients that require Long Term Oxygen Therapy (LTOT). The M50 delivers high concentration of oxygen with flow rates of up to 5 litres per minute. Featuring a compact design, the M50 is easy to move, store and maintain - with just one minute needed to replace any part to be maintained by the user. It also comes with an integrated oxygen sensing device for monitoring oxygen purity levels, and a low oxygen purity alarm to alert the patient if there is low oxygen supply.


- **Lightweight** - a mere 33lbs
- **Quiet operation**
- **Sleek, compact design** for easy handling
- **Adjustable flow rate** up to 5 litres per minute
- **Fire safe oxygen outlet port** for added safety
- **Low oxygen purity alarm**
- **Compatible** with various bubble humidifiers
- **Three year warranty**



Flowmeter



Two stage paper filter

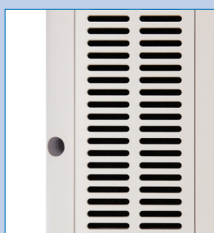

Air filter (Located behind the cover)



Technical Data

Delivery Rate	0.5 - 5 L/min
Flow Mode	Continuous
Oxygen Concentration	87 - 96%
Weight	33 lbs
Dimensions	600 H x 390 W x 230 D mm
Outlet Pressure	0.05+/-105 MPa
Power requirements	110 - 120v/60Hz
Sound Level	42 dBA
Operating Environment	10 - 37 (°C)
Equipment Class	Class IIa
Humidity Range	15% - 95%
Oxygen Purity Alarm	Low oxygen: 82%
Additional Functions	Metal oxygen outlet for fireproofing

Part Numbers

M50 Oxygen Concentrator with USA plug	304051011
Bacteria Filter	14090328
Air filter	501150007
Two Stage Paper filter	501010261
M50 air intake vent cover	501010407
M50 paper filter access panel	411340026
M50 paper filter housing	501010191