



In Wound



On Wound



Around Wound

SILICONE SAP DRESSINGS


Zetuvit® Plus Silicone Border

Zetuvit® Plus Silicone Border is a self-adhesive dressing that is suitable for a wide range of acute and chronic exuding wounds. The absorbent core traps fluid with a unique combination of loose cellulose and superabsorbent polymer (SAP). The silicone coated contact layer and border allows for repeated lifting and repositioning while enabling atraumatic removal.


Order Number	Dressing Size	Pad Size	Selling UOM	Dressings per Box	HCPCS Code	Boxes per Case
413900	3" x 3" / 8cm x 8cm	1.4" x 1.4" / 3.5cm x 3.5cm	Box	10	A6196	6
413901	4" x 4" / 10cm x 10cm	2" x 2" / 5cm x 5cm	Box	10	A6196	6
413902	5" x 5" / 12.5cm x 12.5cm	2.75" x 2.75" / 7.5cm x 7.5cm	Box	10	A6196	6
413903	6" x 6" / 15cm x 15cm	3.7" x 3.7" / 9.5cm x 9.5cm	Box	10	A6196	6
413924	7" x 7" / 17.5cm x 17.5cm	4.7" x 4.7" / 12.1cm x 12.1cm	Box	10	A6197	6
413904	6" x 10" / 16cm x 26cm	4.1" x 8.1" / 10.5cm x 20.5cm	Box	10	A6197	6
413905	8" x 8" / 20cm x 20cm	5.7" x 5.7" / 14.5cm x 14.5cm	Box	10	A6197	6
413906	Sm Sacrum 7" x 7" / 18cm x 18cm	5" x 5" / 12.6cm x 12.6cm	Box	10	A6197	6
413907	Lg Sacrum 9" x 9" / 23cm x 23cm	6.9" x 6.9" / 17.6cm x 17.6cm	Box	10	A6197	6
413908	Oval 5" x 6" / 13cm x 15.5cm	3" x 4" / 7.5cm x 10cm	Box	10	A6196	6
413909	Oval 4.7" x 9" / 12cm x 23cm	2.6" x 6.9" / 6.5cm x 17.5cm	Box	10	A6197	6
413911	Heel 10" x 10" / 25cm x 25cm	7.3" x 6.8" / 18.6cm x 17.2cm	Box	10	A6197	6




RECTANGLE

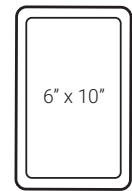


SQUARE

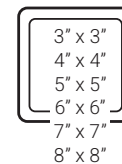


SACRUM

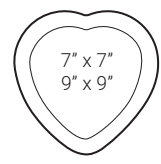

OVAL

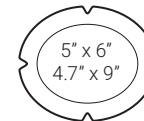

HEEL

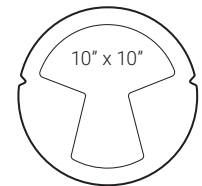


Zetuvit® Plus Silicone (Non-Bordered)

Zetuvit® Plus Silicone is an adhesive, non-border version of Zetuvit® Plus Silicone Border.


Order Number	Dressing Size	Selling UOM	Dressings per Box	HCPCS Code	Boxes per Case
413811	3" x 3" / 8cm x 8cm	Box	10	A6196	12
413812	5" x 5" / 12.5cm x 12.5cm	Box	10	A6197	6
413813	4" x 8" / 10cm x 20cm	Box	10	A6197	6
413814	8" x 8" / 20cm x 20cm	Box	10	A6198	6
413815	8" x 10" / 20cm x 25cm	Box	10	A6198	6



RECTANGLE

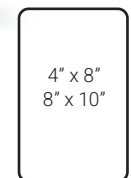

SQUARE

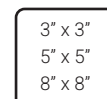

SUPERABSORBENT DRESSINGS


Zetuvit® Plus Superabsorber

Zetuvit® Plus is a non-bordered superabsorbent dressing that does not contain an adhesive layer. The diffusion layer offers quick and uniform distribution of exudate into a superabsorbent core for optimized fluid handling. The dressing can be used with compression wraps or with fixation tapes and bandages.


Order Number	Dressing Size	Selling UOM	Dressings per Box	HCPCS Code	Boxes per Case
413729	2.5" x 2.5" / 6cm x 6cm	Box	10	A6196	12
413600	2.6" x 6.9" / 6.5cm x 17.5cm, "Oval"	Box	10	A6196	6
413110	4" x 4" / 10cm x 10cm	Box	10	A6196	8
413111	4" x 8" / 10cm x 20cm	Box	10	A6197	8
413112	6" x 8" / 15cm x 20cm	Box	10	A6197	8
413601	6.9" x 6.9" / 17.6cm x 17.6cm, "Heart"	Box	10	A6197	6
413602	7.3" x 6.8" / 18.6cm x 17.2cm, "Multisite"	Box	10	A6197	6
413113	8" x 10" / 20cm x 25cm	Box	10	A6198	6
413715	8" x 16" / 20cm x 40cm	Box	10	A6198	5
413739	8" x 24" / 20cm x 60cm	Box	10	A6198	2



RECTANGLE

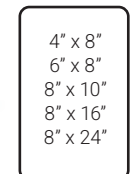

SQUARE

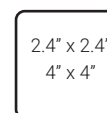

MULTISITE



HEART



OVAL

