



PB6034




PB7035



PBHBMATIS



PBSATINFOAM3680

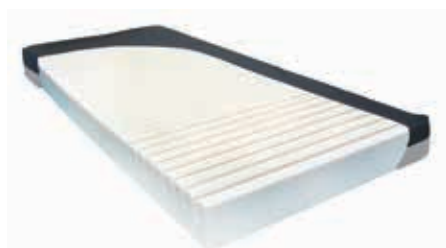

ARUBA-200

Semi-Electric Bed Packages

PBSMB	Semi Electric Bed Frame Only
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PBSMB-FR	Semi Electric Bed w/Full Rails (No Support Surface)
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB6034	ProBasics Full Bed Rails
PBSMB-HR	Semi Electric Bed w/Half Rails (No Support Surface)
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB7035	ProBasics Half Bed Rails
PBSMB-FRFCPKG	Semi Electric Bed w/Full Rails & Fibercore Mattress
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB6034	ProBasics Full Bed Rails
PBHBMATFB	ProBasics Fibercore Mattress
PBSMB-HRFCPKG	Semi Electric Bed w/Half Rails & Fibercore Mattress
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB7035	ProBasics Half Bed Rails
PBHBMATFB	ProBasics Fibercore Mattress
PBSMB-FRISPKG	Semi Electric Bed w/Full Rails & Innerspring Mattress
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB6034	ProBasics Full Bed Rails
PBHBMATIS	ProBasics Innerspring Mattress
PBSMB-HRISPKG	Semi Electric Bed w/Half Rails & Innerspring Mattress
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB7035	ProBasics Half Bed Rails
PBHBMATIS	ProBasics Innerspring Mattress
PBSMB-FRGP1PKG	Semi Electric Bed w/Full Rails & Group 1 Boxed Mattress
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB6034	ProBasics Full Bed Rails
PBSATINFOAM3680	ProBasics Group 1 Compressed Foam Mattress - Boxed
PBSMB-HRGP1PKG	Semi Electric Bed w/Half Rails & Group 1 Boxed Mattress
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB7035	ProBasics Half Bed Rails
PBSATINFOAM3680	ProBasics Group 1 Compressed Foam Mattress - Boxed
PBSMB-FRARBPKG	Semi Electric Bed w/Full Rails & Aruba 200 Mattress
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB6034	ProBasics Full Bed Rails
ARUBA-200	ARUBA Group 1 Sculptured Foam Mattress - Boxed
PBSMB-HRARBPKG	Semi Electric Bed w/Half Rails & Aruba 200 Mattress
HBSMBEM	ProBasics Semi-Electric Bed Ends w/Casters
HBSMSS	ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
PB7035	ProBasics Half Bed Rails
ARUBA-200	ARUBA Group 1 Sculptured Foam Mattress - Boxed



Lightweight Semi-Electric Bed

The **ProBasics Lightweight Semi-Electric Bed** is easy for one person to transport - with a foot section that weighs just 42 lbs. The robust, lightweight frame and snap on or off low-voltage DC motor box ensure quick, efficient set-ups and easy routine maintenance.

Your clients will take comfort that in the event of a power failure with a 9V battery back-up that allows movement to the lowest bed position. Plus, the smart, recessed design of the head and foot spring assemblies help reduce pinching and injuries.

FEATURES

- Easy to transport – foot section weighs just 42 lbs.
- Quick set-up & easy maintenance with motor & electronics that can be removed for service or cleaning
- Head and foot sections can be raised independently or simultaneously
- Large button, simple-to-use hand pendant
- Manual crank at the foot adjusts the frame height
- Scratch-resistant bed ends keep the bed looking new
- 9V battery back-up lets users go to the lowest bed position during power failures (batteries not included)
- Headspring support and foot assembly are recessed to help reduce pinching and injuries
- Trendelenberg & reverse Trendelenberg position capable
- Beds include 2 locking and 2 non-locking casters
- Quiet, efficient ETL-approved DC low voltage single motor box
- HCPCS Code: E0295***

SHIPPING SPECIFICATIONS

Component:	Carton Size:	Shipping Weight:
Foot/Head Springs, Motor Box pendant:	45.5" x 10.5" x 37"	104 lbs.
Head/Foot Bed ends, Casters:	38" x 7.5" x 2.25"	56 lbs.

ACCESSORIES

Works with ProBasics full/half rails, patient trapeze, bed extender and other components.

INDICATIONS FOR USE

Designed and developed for patients with permanent, long-term disabilities and conditions, the ProBasics Electric Hospital Beds are intended to provide assistance to a patient and/or their caregiver for proper positioning of patients who:

1. Require the elevation of the head of the bed 30 or more degrees most of the time due to, but not limited to, congestive heart failure, chronic pulmonary disease, or problems with aspiration.
2. Require special attachments, such as traction equipment, that can only be attached to a hospital bed.
3. Require an adjustable bed to permit transfers to chair, wheelchair, or standing position. This includes, but is not limited to:
 - Severe arthritis and other injuries to lower extremities
 - Spinal cord injuries (including quadriplegia and paraplegia), muscular dystrophies, multiple limb amputations, and stroke, where the patient is able to transfer from a bed to a wheelchair with or without help
4. Require frequent changes in body position or might require an immediate need for a change in body position to avert a life threatening situation.

PRODUCT SPECIFICATIONS

Surface Height:	15" min. - 20" max. (without mattress)
Sleep Surface:	80" L x 36" W
Overall Bed:	87" L x 36" W
Weight Capacity:	450 lb* / Patient Weight 350 lb
Limited Warranty:	5 year on frame** 2 year on motor - mechanical/electrical

Item #	Description	UOM
PBSMB	ProBasics Semi-Electric Bed	each

* - Total distributed weight capacity including user, mattress, bedding and any other items on the bed surface. Please refer to the Owner's Manual for complete details.

** - Limited lifetime warranty is offered on welds, please see Owner's Manual for complete warranty details.

*** - Providers are responsible for determining appropriate billing codes for Medicare Program claims. This information is only included for your convenience and is not intended, nor should be considered as Medicare billing advice.