

Needleless IV Access Devices

Blunt Cannula Split-Septum

A simple split-septum penetrated by a blunt cannula to become the fluid path



Interlink® Components
by BD



Interlink® Site
by Baxter
Flow Rate** 10.3 L/hr
Flow Rate*** 16.7 L/hr





VersaSafe®
by Alaris



SAFELINE®
by B. Braun



LifeShield® Prepierced
by Hospira

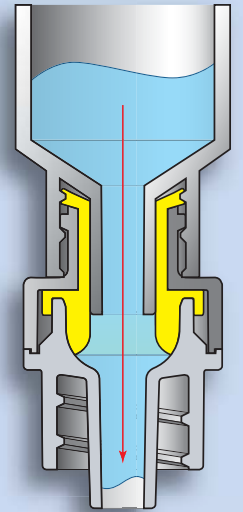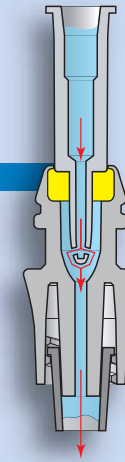

BD Q-Syte™ Luer Access Split-Septum

Flow Rate* 32.0 L/hr

A simple split-septum penetrated by a standard male luer (syringe tip or administration set) to become the fluid path

Luer Access Mechanical Valve

An internal mechanism actuated with a standard male luer (syringe tip or administration set) to open a fluid path



Clearlink®
by Baxter
Flow Rate*
7.3 L/hr



CLAVE®
by Hospira/ICU Medical
Flow Rate*
8.0 L/hr




MicroCLAVE®
by Hospira/ICU Medical
Flow Rate*
5.4 L/hr



InVision-Plus®
by RyMed
Flow Rate*
7.9 L/hr



Smartsite®
by Alaris
Flow Rate*
8.9 L/hr

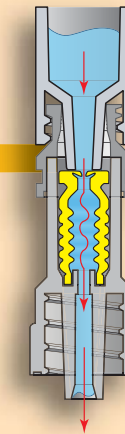

SPLIT-SEPTUM

- Simple
- Clear
- High Flow

Luer Access Mechanical Valve with Positive Displacement Feature

A luer access mechanical valve with the addition of a fluid reservoir that creates positive displacement upon disconnection



Posiflow™
by BD
Flow Rate*
3.2 L/hr



MaxPlus™
by Maximus
Flow Rate*
4.4 L/hr




CLC2000®
by Hospira/ICU Medical
Flow Rate*
14.5 L/hr



ULTRASITE®
by B. Braun
Flow Rate*
12.8 L/hr



SmartSite® Plus
by Alaris
Flow Rate*
11.9 L/hr

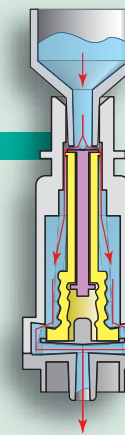

*Test Method; ISO 10555-5 Annex B (max flow @ 525 ± 25ml/min)
**Test Method; ISO 10555-5 Annex B using BD Blunt Plastic Cannula (max flow @ 525 ± 25ml/min)
***Test Method; ISO 10555-5 Annex B using Interlink Lever Lock Cannula (max flow @ 525 ± 25ml/min)

SmartSite and VersaSafe are trademarks of Alaris Medical Systems. ULTRASITE and SAFELINE are trademarks of B. Braun Medical Inc. Interlink and Clearlink are trademarks of Baxter International Inc. CLAVE, MicroCLAVE and CLC2000 are trademarks of ICU Medical, Inc. MaxPlus is a trademark of Maximus Medical Products, Inc. InVision-Plus is a trademark of RyMed Technologies, Inc. BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. ©2007 BD.

