

2020-2021

OSTOMY

Product Guide

Shaping Lives with Complete
Ostomy Products, Care and Support



ConvaTec

Ostomy Care



Together, We Can Help SHAPE LIVES

Innovation begins by listening
and working together.

ConvaTec is a global medical products and technologies company focused on therapies for the management of chronic conditions. Our portfolio of ostomy products, me+™ patient support program and collaboration with health care professionals help to improve the lives of the people we touch.



TABLE OF CONTENTS

me+™ SUPPORT	1	ACUTE KITS	59
OSTOMY OVERVIEW	5	PEDIATRIC PRODUCTS	69
• Introduction to Ostomy Care.....	5	ACCESSORIES	75
• Choosing the Right Ostomy System.....	7	WOUND & FISTULA POUCHES	87
GUIDE TO OSTOMY SYSTEMS	9	DISCONTINUED PRODUCTS	90
• A Guide to Adhesives	9	OSTOMY RESOURCES & SUPPORT	95
• A Guide to Skin Barriers.....	11	MEDICARE UTILIZATION CHART	96
• Moldable Skin Barriers.....	13	PRODUCT CODE INDEX	97
• Accordion Skin Barriers.....	15	REFERENCES	102
• The Purpose of Convexity.....	17		
• Ostomy Pouch Options.....	21		
• Common Pouch Features.....	23		
TWO-PIECE SYSTEMS	25		
• Natura™, Natura™+, and Sur-FIT® Natura™ Systems.....	27		
• Esteem Synergy® Systems.....	41		
ONE-PIECE SYSTEMS	45		
• Esteem™+ Pouches	46		
• ActiveLife® Pouches.....	55		

THERE IS A NEED FOR ONGOING PATIENT SUPPORT

The **2nd most common** reason for hospital readmissions is ileostomy and other enterostomy surgery¹

40%

of readmissions following colorectal surgery *may be preventable*²

30%

of ostomy patients exhibited peristomal skin complications within 90 days of surgery³

49%

of respondents reported a need for *adequate information and communication post-op*⁴

Most common ostomy self-care challenges are related to *leakage, skin issues, frequency of changes and difficulty caring for an ostomy*⁵



THE **me+**TM PROGRAM SUPPORTS PATIENTS

The me+TM program is a comprehensive support program designed to complement clinical care by providing patients with support, inspiration, insights and ideas throughout their entire ostomy journey.

me+TM program members can:

- ✓ Request free ostomy samples
- ✓ Speak to a WOC nurse or product specialist
- ✓ Connect to our me+ community
- ✓ Receive ongoing lifestyle education, support and much more

Health Care Professionals

TO ENROLL YOUR PATIENTS in the me+ program

call us at 877-585-0470 or email concierge@convatec.com

Not a Health Care Professional?

CONTACT US TODAY to join the me+ program

call 1-800-422-8811 or email CIC@convatec.com



Our goal is to *partner with you and your patients to deliver best-in-class care* throughout the patient journey:

✓ **Pre-Surgery / Surgery**

- Access online resources at mepluscare.com
- Request me+™ patient Starter Kit and education materials
- Enroll in me+ program
- Learn about the me+ recovery series, developed to aid rehabilitation after surgery

✓ **Weeks Immediately After Surgery**

- Consult with me+ WOC nurses and product specialists
- Personalized Starter Kit
- me+ welcome email and support
- me+ recovery phase 1 exercise booklet*

✓ **First Year & Beyond**


- Explore online resources and tools
- Get inspired by the me+ community
- Sample new products
- Continue me+ recovery series exercises*

Throughout their journey, your patients can call **1-800-422-8811** for live support from me+™ WOC nurses and product specialists

A me+™ Concierge team member can assist Healthcare professionals with:

- ✓ Personalized support
- ✓ Enrolling your patients into the me+™ program
- ✓ Providing a me+ Starter Kit and free product samples for your patients

TO LEARN MORE contact your me+ concierge directly at 877-585-0470

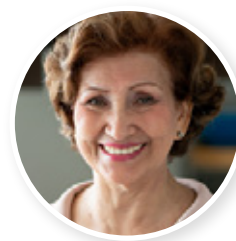

"You guys have been a rock for me. Your expert insight and help picking the right product for my needs has been crucial at setting my mind at ease, and feeling that even with this significant life change that it's all going to be alright." - Sofia†

*Model portrayal for illustrative purposes only.

*The information in the series is not intended to provide medical advice. Patients and users should speak with their health care professional before engaging in physical activity following ostomy surgery. In the event they experience pain or discomfort, they should stop the activity immediately and contact their health care professional.

An Introduction to OSTOMY CARE

There are 3 primary types of ostomy surgeries. Understanding the fundamental difference is important to select the right product and provide the right support.



COLOSTOMY (large intestine)

A colostomy is a surgical procedure that creates an opening on the abdomen through which the large intestine is redirected to bypass diseased or damaged parts of the large intestine. The large intestine is brought to the surface of the skin, and stitched to the abdomen to form a new exit for the stool from the bowel.

Things to know about a Colostomy

- ✓ Primarily formed due to colon or rectal cancer or trauma
- ✓ Can be temporary or permanent
- ✓ Stool tends to be semi-formed to formed and predictable 1-2 times per day

ILEOSTOMY (small intestine)

An ileostomy is a surgical procedure that creates an opening on the abdomen through which the end of your small intestine is redirected to bypass diseased or damaged parts of the digestive system. The small intestine is brought to the surface of the skin, and stitched to your abdomen to form a new exit for waste from your bowels.

Things to know about an Ileostomy

- ✓ Primarily formed due to inflammatory bowel disease (Crohn's, Ulcerative Colitis)
- ✓ Can be temporary or permanent
- ✓ Stool is liquid or pasty and constant throughout the day
- ✓ The small intestine contains digestive enzymes and acids, and the output may cause skin irritation if it comes in contact with the skin

UROSTOMY (urinary tract)

A urostomy is a surgical procedure that creates an opening on the abdomen through which urine is redirected away from a bladder that is damaged or diseased. After the bladder is removed, a piece of small intestine (ileum) is separated from the rest of the bowel. This piece is used as a conduit for urine to pass out of the body. The ureters are attached to this new passageway and one end is brought up to an opening made on your abdomen.

Things to know about a Urostomy

- ✓ Primarily formed due to cancer or inflammation of the bladder
- ✓ Is always permanent
- ✓ Urine flows through the stoma as it is produced which means output is constant
- ✓ Urine may cause skin irritation if it comes in contact with the skin

Choosing the Right OSTOMY SYSTEM

The ostomy pouch collects stoma output. The pouch comes as part of a two-piece or one-piece system.



2 PIECE SYSTEM

The two-piece system includes a skin barrier and a separate pouch that securely join together. This option allows the pouch to be removed to empty or change, without removing the skin barrier as frequently.



Designed for patients

- ✓ who want to change their pouch more frequently than their skin barrier
- ✓ with an ileostomy who need to empty/change the pouch frequently
- ✓ who exercise, travel and play sports

1 PIECE SYSTEM

The skin barrier and pouch are provided as an integrated solution that cannot be separated. It offers one simple and convenient step to attach or remove, providing a fresh start every time. This option is thin and flexible, and designed for more frequent system changes.



Designed for patients

- ✓ post-operatively – removes any added pressure on the abdomen caused when connecting a skin barrier and pouch
- ✓ who prefer to change their ostomy system frequently
- ✓ who prefer a flexible and low profile ostomy system

A Guide to ADHESIVES

The adhesive is an important component of the skin barrier. Our product portfolio features a range of adhesive technologies to enable you to select an appropriate solution that meets the needs of your patients.



STOMAHESIVE® (standard wear)

Designed for people with mostly solid output, our unique Stomahesive® technology is a hydrocolloid adhesive that gently adheres to the skin, protecting it from the damaging effects of stoma effluent.

We also offer a Modified Stomahesive® option, which features a touch of mineral oil that makes the barrier more flexible during wear, and removal easier on the skin.

Designed for patients

- ✓ with semi-formed/formed output
- ✓ who have sensitive skin
- ✓ with frequent changes
- ✓ who prefer a standard weartime

DURAHESIVE® (extended wear)

Designed for people with mostly liquid output, our unique Durahesive® technology is a hydrocolloid adhesive that swells up upon hydration to protect the stoma. This special effect called turtlenecking, helps create a more secure seal without harming the skin.

We also offer Durahesive® Plus, which features an enhanced tack, strength and cohesion designed for extended weartime or aggressive output.

Designed for patients

- ✓ with semi-liquid/liquid output, including urine
- ✓ who want or need longer weartime

ESTEEM™ + Soft & Flex CONVEX

Adhesive designed with a high initial tack, to aid in the initial attachment to the skin. The hydrocolloid barrier features a thin tapered edge to help prevent the edge catching on clothing, and a thicker center to help absorption around the stoma.⁶

Designed for patients

- ✓ with both solid and liquid output
- ✓ who prefer a standard weartime

TAPE COLLARS
are available in two forms:

Acrylic: Tape collar designed with good tack
Hydrocolloid: Designed for those with sensitivity to other tape collars

A Guide to SKIN BARRIERS

The skin barrier is the piece of the pouching system that sticks to the body. It should fit snugly around the stoma to help prevent leakage and have an opening customized to stoma size and shape.



CUT-TO-FIT

Cut-to-Fit barriers can be cut to match the shape and size of individual stomas. Using scissors, one can cut a circle, oval or any shape needed to fit the stoma.




Designed for patients

- ✓ with a stoma size in-between pre-cut sizes
- ✓ with an irregular stoma shape

PRE-CUT

Pre-cut barriers have a uniform, circular stoma hole and are ready to apply right out of the box. They are available for the most common stoma sizes.

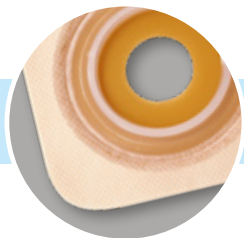


Designed for patients

- ✓ with an evenly round stoma that fits the available pre-cut size
- ✓ with dexterity issues who find cutting challenging

MOLDABLE

ConvaTec Moldable Technology™ skin barrier with Rebounding Memory Technology™ creates an elastic-like seal that rebounds to fit most shapes and sizes, and then “turtlenecks” upward, to enhance the seal and protect the skin.⁷

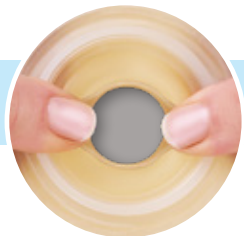

Designed for patients

- ✓ who have large, round, oval or irregular shaped stomas
- ✓ with dexterity issues who find cutting challenging

MOLDABLE

Skin Barriers



ConvaTec Moldable Technology™ Skin Barriers are designed to hug the stoma opening for a more custom and secure fit.

CLINICALLY PROVEN


to prevent and improve peristomal skin complications⁷

Up to 55% of ostomates have skin issues⁸

- 95% of people who started on ConvaTec Moldable Technology kept their skin healthy⁷
- 86% of people who switched to ConvaTec Moldable Technology improved their skin condition⁷

Features and Benefits


- ✓ Reduced application steps—fewer steps than traditional cut-to-fit skin barriers
- ✓ No cutting needed
- ✓ Memory – flexible materials gently expand and contract with the stoma
- ✓ Turtlenecking – the adhesive gently swells when in contact with output & stool to help protect the skin and secure the seal.

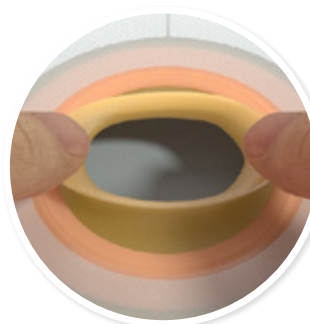

Roll, mold, and hold to the stoma's shape and size



Turtlenecks around the stoma to protect against leaks

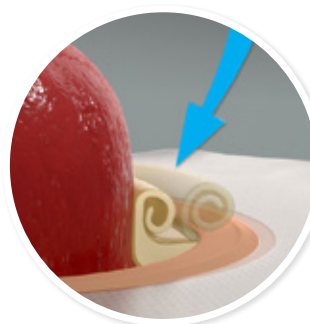

Rebounding Memory Technology™ - rolls back to achieve a custom fit

ACCORDION

Skin Barriers

The Accordion Flange is an option on a range of Natura® skin barriers. The accordion flange lifts to provide generous finger room so coupling is easy and comfortable.



PRODUCT FEATURES

- 2X more finger space than any other lifting flange⁹
- Simple to measure and cut or mold
- Range of flat, convex, Cut-to-Fit and Moldable options
- Compatible with SUR-FIT® Natura™ and Natura™ + pouch option

Designed for patients

- ✓ with dexterity issues as it provides finger room so coupling is easy and comfortable
- ✓ designed to prevent pressure or trauma if abdominal tenderness is present due to post-operative stoma surgery/incisions

HOW IT WORKS

- Designed to help minimize abdominal pressure and pain when attaching the pouch.
- Connects with an audible 'click' and physical 'snap' to confirm the pouch is firmly attached.



Lift 'UP' flange to easily attach the pouch away from the body



Low profile in the 'DOWN' position so it's not visible under clothing

The Purpose of CONVEXITY



WHAT IS CONVEXITY?

All bodies are different. Convex wafers are designed to help manage irregular stomas and abdomens to help reduce leakage and simplify stoma management.

- Convexity is the outward curving of the baseplate or skin barrier.^{10,11}
- Convexity applies a uniform pressure to the immediate peristomal skin, improving the projection of the stoma for drainage.¹²
- Convexity helps maintain a secure seal between the pouch and the skin.¹³

GROWING CONVEXITY NEED

Data shows the need for convex skin barriers growing year on year.¹⁴ This is due to a number of factors including:


- Changing surgical techniques/outcomes of stoma shape & profile
- Increased obesity rates

WHEN TO USE CONVEXITY




Flush Stomas

Located at the same level as the abdomen skin surface, and does not protrude¹⁵

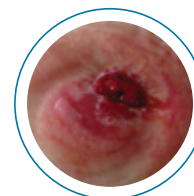


Retracted Stomas

Stoma pulled into the abdomen below skin level or having little or no spout¹⁶

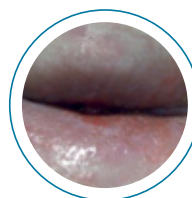

Skin Folds

A stoma which is deep in a skin fold/crease¹³

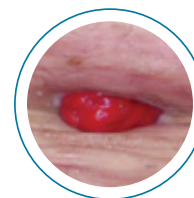

Soft or Flaccid Abdomens

An abdomen requiring firm level of support to stabilize the pouching system¹⁵

The Characteristics of CONVEXITY



RIGIDITY

Firm/Rigid Convexity

Typically used on soft or flaccid abdomens.


Soft Convexity

Typically used on firm or toned abdomens.

DEPTH

Height of the convex plateau.

- ✓ Deep
- ✓ Shallow

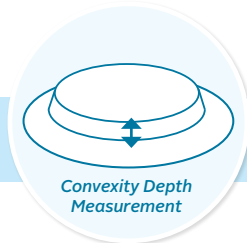


PLATEAU WIDTH

Distance across the raised flat surface from edge-to-edge

Why Plateau Width Matters

Successful convexity must press down on the peristomal skin enough to make a good seal around the base of the stoma and help with protrusion.¹³

- ✓ **Convexity Pressure:** A narrow plateau width may be used to create pressure directly around the stoma with lesser pressure as the plateau width gets wider.
- ✓ **Plateau Width:** A wider plateau width may be used to flatten out skin creases away from the stoma.

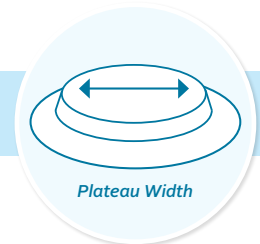

Choosing the right convexity solution for different body profiles can be challenging. Please feel free to consult with one of our me+™ Specialists or WOC Nurses for further guidance.

Ostomy Pouch OPTIONS

ConvaTec offers a wide range of pouching options across both 2 Piece and 1 Piece Systems.



CLOSED

A closed pouch doesn't require draining. Once they are full, you simply remove and replace with a new pouch.

Additionally, ideal during physical activity, intimate moments, swimming, and travel.



Designed for patients

- ✓ who have thick, formed stool that is predictable (1-2 times per day)

DRAINABLE

Drainable pouches allow you to drain the output periodically throughout the day so that you do not have to change pouches as frequently.




Designed for patients

- ✓ who have loose, liquid stool or semi-formed stool that can be unpredictable or constant throughout the day

UROSTOMY

A Urostomy pouch has a faucet tap, soft tap, or fold over tap at the bottom that allows for draining urine as frequently as needed.

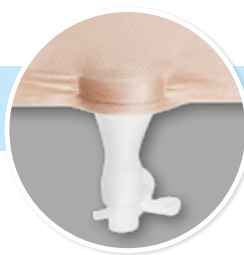

Designed for patients

- ✓ with a Urostomy

Common POUCH FEATURES



POUCH COLOR



Transparent (*Clear front of pouch*)

Clear pouch front for easy output assessment and determining when a pouch needs to be changed or emptied. Soft and quiet material on the body side of the pouch.

Opaque (*Tan-colored pouch*)

Soft and quiet material cover both the front and the body side of the pouch for more discretion and comfort.



Easy-View Window

Some options come with an easy-view window for accurate positioning and observation.

FILTER or NO FILTER



Filters help release and deodorize gas and keep the closed or drainable pouch flat. Note that the filter must be covered with a filter cover while swimming, bathing or showering.

DRAINABLE POUCH CLOSURE



InvisiClose™ Clipless Tail Closure

Self closing drainable pouch closure which provides an easy closure option that does not require a separate clip for added security.



Lock-it Pocket™

Folds inside out to hide away the secured drainage tail—providing the flexibility of a drainable pouch, with the feel and comfort of a closed-end pouch.

Traditional drainable pouches are closed with a separate clip.

UROSTOMY POUCH DRAINAGE



The Accuseal™ Tap with valve is a faucet type tap designed for easy manipulation and emptying. Other drainage options include a Soft Tap or Fold Over Tap.

All Urostomy pouches include an anti-reflux valve which is designed to help prevent backflow.

TWO-PIECE Systems


Two-piece systems include a **skin barrier and a separate pouch** that securely join together. There are two brands that make up ConvaTec's Two-Piece Ostomy Systems: Natura™ + and Esteem Synergy® +.



NATURA™ 



The Natura™ brand (inclusive of Natura™, Sur-Fit Natura™, and Natura™+) is an extensive range of skin barriers and pouches that connect together with a plastic ring called a flange in a simple snap-off, snap-on design.

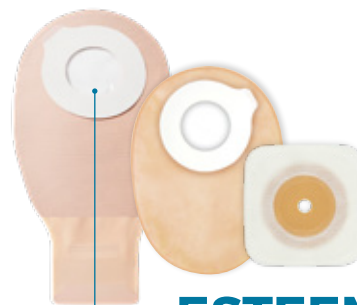







ESTEEM SYNERGY® 






The Esteem Synergy®+ brand are skin barriers and pouches that connect together with Adhesive Coupling Technology™. This is a flat, flexible adhesive foam ring on the pouch that sticks to a flat 'landing zone' on the skin barrier.

Color Key & Sizing Guide

COLOR MATCH						
FLANGE SIZE	1 1/4" (32mm)	1 1/2" (38mm)	1 3/4" (45mm)	2 1/4" (57mm)	2 3/4" (70mm)	4" (100mm)

Color Key & Sizing Guide

COLOR MATCH			
COUPLING SIZE	Small Stoma Opening	Medium Stoma Opening	Large Stoma Opening



When ordering a Two-Piece System, **make sure you select the same flange/coupling size for the skin barrier and pouch.** Use the colored dots in the table to help.

NATURA™ DURAHESIVE® SKIN BARRIERS

Compatible with Natura™ + and SUR-FIT® Natura™ pouches

 Recommended Product



 SUR-FIT® NATURA™ SKIN BARRIER

Cut-to-Fit

- Durahesive® (Extended Wear)
- Quantity = 10 (*Box of 5)

HCPCS A4409

(413157 and 401905 are HCPCS A4410)

Flange Size	Cutting Range	Acrylic Tape Collar WHITE	Acrylic Tape Collar TAN	Without Tape Collar
● 1 1/4" (32mm)	1/2 - 3/4"	413159	413164	413153
● 1 1/2" (38mm)	1/2 - 1"	413160	413165	413154
● 1 3/4" (45mm)	1/2 - 1 1/4"	413161	413166	413155
● 2 1/4" (57mm)	1/2 - 1 3/4"	413162	413167	413156
● 2 3/4" (70mm)	1/2 - 2 1/4"	413163	413168	413157
● 4" (100mm)	1/2 - 3 1/2"	N/A	N/A	401905*



 NATURA™ SKIN BARRIER

Cut-to-Fit | Accordion Flange

- Durahesive® (Extended Wear)
- Quantity = 10

HCPCS A4409

Flange Size	Cutting Range	Acrylic Tape Collar WHITE
● 1 3/4" (45mm)	1/2 - 13/16" (13mm-21mm)	421454
● 2 1/4" (57mm)	1/2 - 1 5/16" (13mm-33mm)	421458
● 2 3/4" (70mm)	1/2 - 1 3/4" (13mm-45mm)	421462



 NATURA™ SKIN BARRIER

Moldable | Accordion Flange

- Durahesive® (Extended Wear)
- Quantity = 10

HCPCS A4409

Flange Size	Stoma Size	Hydrocolloid Tape Collar TAN
● 2 1/4" (57mm)	1/2 - 7/8" (13-22mm)	421039
● 2 1/4" (57mm)	7/8 - 1 1/4" (22-33mm)	421040
● 2 3/4" (70mm)	1 1/4 - 1 3/4" (33-45mm)	421041

 Recommended Product

 NATURA™ SKIN BARRIER

Cut-to-Fit | Convex | Accordion Flange

- Durahesive® (Extended Wear)
- Convex Wafer Depth = 6.65mm
- Quantity = 10

HCPCS A4407

Flange Size	Cutting Range	Acrylic Tape Collar TAN	Hydrocolloid Tape Collar TAN
● 1 3/4" (45mm)	1/2 - 7/8" (13-22mm)	421638	421639
● 2 1/4" (57mm)	1/2 - 1 3/8" (13-35mm)	421640	421641
● 2 3/4" (70mm)	1/2 - 1 7/8" (13-48mm)	421642	421643


 NATURA™ SKIN BARRIER

Moldable

- Durahesive® (Extended Wear)
- Quantity = 10

HCPCS A4409

Flange Size	Stoma Size	Acrylic Tape Collar WHITE	Hydrocolloid Tape Collar WHITE
● 1 3/4" (45mm)	1/2 - 7/8" (13-22mm)	413417	411800
● 1 3/4" (45mm)	7/8 - 1 1/4" (22-33mm)	413418	411802
● 2 1/4" (57mm)	1 1/4 - 1 3/4" (33-45mm)	413419	411804
● 2 3/4" (70mm)	1 1/4 - 1 3/4" by 1 3/4 - 2 1/8" (33x45mm by 45x56mm) Oval shape	413420	411806


 NATURA™ SKIN BARRIER

Moldable | Convex

- Durahesive® (Extended Wear)
- Quantity = 10

HCPCS A4407

Flange Size	Stoma Size	Acrylic Tape Collar WHITE	Convex Wafer Depth
● 1 3/4" (45mm)	1/2 - 7/8" (13-22mm)	404592	6.35mm
● 1 3/4" (45mm)	7/8 - 1 1/4" (22-33mm)	404593	7mm
● 2 1/4" (57mm)	1 1/4 - 1 3/4" (33-45mm)	404594	7mm

SUR-FIT® NATURA™ SKIN BARRIER



Pre-cut | Convex

- Durahesive® (Extended Wear)
- Convex Wafer Depth = 6.65mm
- Quantity = 10

HCPCS A4407

Flange Size	Stoma Size	Acrylic Tape Collar WHITE
● 1 3/4" (45mm)	1/2" (13 mm)	413177
● 1 3/4" (45mm)	5/8" (16 mm)	413178
● 1 3/4" (45mm)	3/4" (19 mm)	413179
● 1 3/4" (45mm)	7/8" (22 mm)	413180
● 1 3/4" (45mm)	1" (25 mm)	413181
● 1 3/4" (45mm)	1 1/8" (28 mm)	413182
● 1 3/4" (45mm)	1 1/4" (32 mm)	413183
● 1 3/4" (45mm)	1 3/8" (35 mm)	413184
● 2 1/4" (57mm)	1 1/2" (38 mm)	413185
● 2 1/4" (57mm)	1 5/8" (41 mm)	413186
● 2 1/4" (57mm)	1 3/4" (45 mm)	413187
● 2 1/4" (57mm)	2" (50 mm)	413188

NATURA™ STOMAHESIVE® SKIN BARRIERS

Compatible with Natura™ + and SUR-FIT® Natura™ pouches

 Recommended Product









 SUR-FIT® NATURA™ SKIN BARRIER


Cut-to-Fit

- Stomahesive® (Standard Wear)
- Quantity = 10 (*Box of 5)

HCPCS A4414

(401577 and 401906 are HCPCS A4415)

Flange Size	Cutting Range	Acrylic Tape Collar WHITE	Acrylic Tape Collar TAN	Without Tape Collar
 1 1/4" (32mm)	1/2 - 3/4"	125257	125262	401573
 1 1/2" (38mm)	1/2 - 1"	125258	125263	401574
 1 3/4" (45mm)	1/2 - 1 1/4"	125259	125264	401575
 2 1/4" (57mm)	1/2 - 1 3/4"	125260	125265	401576
 2 3/4" (70mm)	1/2 - 2 1/4"	125261	125266	401577
 4" (100mm)	1/2 - 3 1/2"	N/A	N/A	401906*

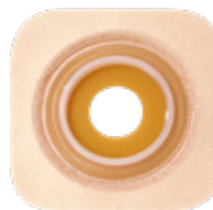












 SUR-FIT® NATURA™ SKIN BARRIER

Pre-cut

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Flange Size	Stoma Size	Acrylic Tape Collar TAN
 1 3/4" (45mm)	1/2" (13 mm)	125267
 1 3/4" (45mm)	5/8" (16 mm)	125268
 1 3/4" (45mm)	3/4" (19 mm)	125269
 1 3/4" (45mm)	7/8" (22 mm)	125270
 1 3/4" (45mm)	1" (25 mm)	125271
 1 3/4" (45mm)	1 1/8" (28 mm)	125272
 1 3/4" (45mm)	1 1/4" (32 mm)	125273
 1 3/4" (45mm)	1 3/8" (35 mm)	125274
 2 1/4" (57mm)	1 1/2" (38 mm)	125275
 2 1/4" (57mm)	1 5/8" (41 mm)	125276
 2 1/4" (57mm)	1 3/4" (45 mm)	125277






 NATURA™ SKIN BARRIER

Cut-to-Fit | Accordion Flange

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Flange Size	Cutting Range	Acrylic Tape Collar WHITE
 1 3/4" (45mm)	1/2 - 13/16" (13-21mm)	421456
 2 1/4" (57mm)	1/2 - 1 5/16" (13-33mm)	421460
 2 3/4" (70mm)	1/2 - 1 3/4" (13-45mm)	421464

 Recommended Product

 NATURA™ SKIN BARRIER
Moldable | Accordion Flange

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Flange Size	Stoma Size	Hydrocolloid Tape Collar TAN
● 2 1/4" (57mm)	1/2 - 7/8" (13-22mm)	421033
● 2 1/4" (57mm)	7/8 - 1 1/4" (22-33mm)	421034
● 2 3/4" (70mm)	1 1/4 - 1 3/4" (33-45mm)	421035


 NATURA™ SKIN BARRIER
Moldable

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Flange Size	Stoma Size	Acrylic Tape Collar WHITE	Hydrocolloid Tape Collar WHITE
● 1 3/4" (45mm)	1/2 - 7/8" (13-22mm)	413421	411801
● 1 3/4" (45mm)	7/8 - 1 1/4" (22-33mm)	413422	411803
● 2 1/4" (57mm)	1 1/4 - 1 3/4" (33-45mm)	413423	411805
● 2 3/4" (70mm)	1 1/4 - 1 3/4" by 1 3/4 - 2 1/8" (33x45mm by 45x56mm) Oval shape	413424	411807

NATURA™ + POUCHES

Compatible with SUR-FIT® Natura™ and Natura™ skin barriers

 Recommended Product

 NATURA™+ POUCH
Drainable | InvisiClose™ Tail Closure

- Pouch Length = 12" (*pouch length = 14")
- Quantity = 10

HCPCS A5063 | No Filter

	Flange Size	Transparent	Opaque
●	1 3/4" (45mm)	416418	421738
●	2 1/4" (57mm)	416421	421739
●	2 3/4" (70mm)	416424	421740
●	4" (100mm)	416472*	N/A

HCPCS A4425 | Filter (Filter covers included)

	Flange Size	Transparent	Opaque w/ Easy-View Window	Opaque
●	1 1/2" (38mm)	416414	421746	416415
●	1 3/4" (45mm)	416416	421747	416417
●	2 1/4" (57mm)	416419	421748	416420
●	2 3/4" (70mm)	416422	421749	416423

SUR-FIT® NATURA™ POUCH

Drainable | InvisiClose™ Tail Closure

- Pouch Length = 12"
- Quantity = See table below



HCPCS A5063 | No Filter

	Flange Size	Transparent Qty = 10	Opaque Qty = 20
●	1 3/4" (45mm)	413312	413334
●	2 1/4" (57mm)	413313	N/A
●	2 3/4" (70mm)	413314	N/A

HCPCS A4425 | Filter (Filter covers included)

	Flange Size	Transparent Qty = 10	Opaque Qty = 20
●	1 3/4" (45mm)	411360	411310
●	2 1/4" (57mm)	411361	411311
●	2 3/4" (70mm)	411362	N/A

SUR-FIT® NATURA™ POUCH

Drainable | Clipped Closure

- Pouch Length = 10"
- 2 Tail Closures per box
- Quantity = See table below



HCPCS A5063 | No Filter

	Flange Size	Opaque Qty = 10	Opaque Qty = 20
●	1 1/4" (32mm)	401505	404030
●	1 1/2" (38mm)	401506	404031
●	1 3/4" (45mm)	401507	404032
●	2 1/4" (57mm)	401508	404033
●	2 3/4" (70mm)	401509	404034

SUR-FIT® NATURA™ POUCH

Drainable | Clipped Closure

- Pouch Length = 12" (*pouch length = 14"; **pouch length = 10")
- 2 Tail Closures per box
- Quantity = See table below



HCPCS A5063 | No Filter

	Flange Size	Transparent Qty = 10	Opaque Qty = 10	Opaque Qty = 20	Opaque w/ 1 comfort panel Qty = 10	Opaque w/ 1 comfort panel Qty = 20
●	1 1/4" (32mm)	401510	401500	404025	401932	404014
●	1 1/2" (38mm)	401511	401501	404026	401933	404015
●	1 3/4" (45mm)	401512	401502	404027	401934	404016
●	2 1/4" (57mm)	401513	401503	404028	401935	404017
●	2 3/4" (70mm)	401514	401504	404029	401936	404018
●	4" (100mm)	401515*	N/A	N/A	N/A	N/A

HCPCS A4425 | Filter (Filter covers included)

	Flange Size	Transparent Qty = 10	Opaque Qty = 20	Opaque Qty = 20**
●	1 1/4" (32mm)	N/A	411489	411288**
●	1 1/2" (38mm)	411264	411490	411289**
●	1 3/4" (45mm)	411265	411491	411290**
●	2 1/4" (57mm)	411266	411492	411291**
●	2 3/4" (70mm)	411267	411493	411292**

 Recommended Product

SUR-FIT® NATURA™ POUCH

Drainable | High-Output

- Anti-Reflux Valve
- Outlet with spout and cap for high output drainage
- Pouch Length = 14"
- 1 Tail Closure, 1 Night Drainage Adapter, and 5 Replacement Filters per box
- Quantity = 5



HCPCS A4413 | Filter

	Flange Size	Transparent	Opaque
●	1 3/4" (45mm)	420695	401557
●	2 1/4" (57mm)	420696	401558
●	2 3/4" (70mm)	420697	401559


 NATURA™+ POUCH

Closed-End

- Pouch Length = 8"
- Quantity = 30

HCPCS A5054 | No Filter

	Flange Size	Opaque w/ Easy-View Window
●	1 1/4" (32mm)	421677
●	1 1/2" (38mm)	421679
●	1 3/4" (45mm)	421681
●	2 1/4" (57mm)	421683
●	2 3/4" (70mm)	421893

HCPCS A4419 | Filter (Filter covers included)

	Flange Size	Opaque w/ Easy-View Window	Opaque
●	1 1/4" (32mm)	421676	N/A
●	1 1/2" (38mm)	421678	416403
●	1 3/4" (45mm)	421680	416406
●	2 1/4" (57mm)	421682	416409
●	2 3/4" (70mm)	421798	416412

NATURA™+ POUCH

Closed-End

- Pouch Length = 6"
- Quantity = 30



HCPCS A4419 | Filter (Filter covers included)

	Flange Size	Opaque w/ Easy-View Window
●	1 3/4" (45mm)	421796
●	2 1/4" (57mm)	421797

SUR-FIT® NATURA™ POUCH

Closed-End

- Pouch Length = 8"
- Quantity = See table below



HCPCS A5054 | No Filter

	Flange Size	Opaque Qty = 30	Opaque Qty = 60
●	1 3/4" (45mm)	401522	413170
●	2 1/4" (57mm)	401523	413171
●	2 3/4" (70mm)	401524	N/A

HCPCS A4419 | Filter (Filter covers included)

	Flange Size	Opaque Qty = 30	Opaque Qty = 60
●	1 3/4" (45mm)	401526	413174
●	2 1/4" (57mm)	401527	413175
●	2 3/4" (70mm)	401528	N/A



SUR-FIT® NATURA™ POUCH

Closed-End | Flange Cap

- Quantity = 25



HCPCS A5055 | Filter (Filter covers included)

	Flange Size	Opaque w/ 1 comfort panel
●	1 3/4" (45mm)	401909



SUR-FIT® NATURA™ POUCH

Urostomy | Accuseal™ Tap with Valve

- Anti-Reflux Valve
- Pouch Length = 10" (*pouch length = 9")
- 2 Accuseal Night Drainage Adapters per box
- Quantity = 10

HCPCS A4432

	Flange Size	Transparent	Opaque	Opaque*
●	1 1/4" (32mm)	401542	401551	401547*
●	1 1/2" (38mm)	401543	401552	401548*
●	1 3/4" (45mm)	401544	401553	401549*
●	2 1/4" (57mm)	401545	401554	401550*
●	2 3/4" (70mm)	401546	401555	N/A

NATURA™+ POUCH

Urostomy | Soft Tap

- Anti-Reflux Valve
- Pouch Length = 10" (*pouch length = 9")
- 2 Soft Tap Drainage Adapters per box
- Quantity = 10



HCPCS A5073

	Flange Size	Transparent	Transparent*
●	1 1/4" (32mm)	N/A	421074*
●	1 3/4" (45mm)	413437	N/A
●	2 1/4" (57mm)	413438	N/A
●	2 3/4" (70mm)	413439	N/A

SUR-FIT® NATURA™ POUCH

Urostomy | Fold Over Tab

- Anti-Reflux Valve
- Pouch Length = 10" (*pouch length = 9")
- 2 Universal Night Drainage Adapters per box
- Quantity = 10



HCPCS A5073

	Flange Size	Transparent	Transparent*
●	1 1/4" (32mm)	401533	401538*
●	1 1/2" (38mm)	401534	401539*
●	1 3/4" (45mm)	401535	401540*
●	2 1/4" (57mm)	401536	401541*
●	2 3/4" (70mm)	401537	N/A

ESTEEM SYNERGY® SKIN BARRIERS

Compatible with Esteem Synergy®+ Pouches

ESTEEM SYNERGY® SKIN BARRIER



Cut-to-Fit

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Coupling Size	Cutting Range	Acrylic Tape Collar WHITE
● Small Stoma Opening	up to 1 3/8" (35mm)	405456
● Medium Stoma Opening	up to 1 7/8" (48mm)	405457
● Large Stoma Opening	up to 2 3/8" (61mm)	405458

HCPCS A4415 (*HCPCS is A4414)

Coupling Size	Cutting Range	Without Tape Collar
● Small Stoma Opening	up to 1 3/8" (35mm)	405467*
● Medium Stoma Opening	up to 1 7/8" (48mm)	405468
● Large Stoma Opening	up to 2 3/8" (61mm)	403948

ESTEEM SYNERGY® SKIN BARRIER



Moldable

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Coupling Size	Stoma Size	Hydrocolloid Tape Collar WHITE
● Medium Stoma Opening	1/2 - 7/8" (13-22mm)	411657
● Medium Stoma Opening	7/8 - 1 1/4" (22-33mm)	411658
● Large Stoma Opening	1 1/4 - 1 3/4" (33-45mm)	411659

ESTEEM SYNERGY® SKIN BARRIER



Moldable | Convex

- Durahesive® (Extended Wear)
- Convex Wafer Depth = 6.65mm
- Quantity = 10

HCPCS A4407

Coupling Size	Stoma Size	Acrylic Tape Collar WHITE
● Small Stoma Opening	1/2 - 7/8" (13-22mm)	409268
● Medium Stoma Opening	7/8 - 1 1/4" (22-33mm)	409269
● Large Stoma Opening	1 1/4 - 1 3/4" (33-45mm)	409270

ESTEEM SYNERGY® SKIN BARRIER



Pre-Cut

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Coupling Size	Stoma Size	Acrylic Tape Collar WHITE
● Small Stoma Opening	5/8" (16 mm)	405473
● Small Stoma Opening	3/4" (19 mm)	405474
● Small Stoma Opening	7/8" (22 mm)	405475
● Small Stoma Opening	1" (25 mm)	405476
● Small Stoma Opening	1 1/8" (28 mm)	405477
● Small Stoma Opening	1 1/4" (32 mm)	405478
● Small Stoma Opening	1 3/8" (35 mm)	405479
● Medium Stoma Opening	1 1/2" (38 mm)	405480
● Medium Stoma Opening	1 5/8" (41 mm)	405481
● Medium Stoma Opening	1 3/4" (45 mm)	405482

ESTEEM SYNERGY® POUCHES

Compatible with Esteem Synergy® Skin Barriers

ESTEEM SYNERGY® + POUCH

Drainable | InvisiClose™ Tail Closure

- Pouch Length = 12"
- Quantity = 10



HCPCS A5063 | No Filter

Coupling Size	Transparent
● Large Stoma Opening	416903

HCPCS A4425 | Filter (Filter covers included)

Coupling Size	Transparent	Opaque
● Small Stoma Opening	416795	416796
● Medium Stoma Opening	416798	416799
● Large Stoma Opening	416901	416902

ESTEEM SYNERGY® + POUCH

Closed-End

- Pouch Length = 8"
- Quantity = 30



HCPCS A4419 | Filter (Filter covers included)

Coupling Size	Transparent	Opaque
● Small Stoma Opening	416774	416775
● Medium Stoma Opening	416778	416779
● Large Stoma Opening	416782	416783

ESTEEM SYNERGY® + POUCH

Closed-End

- Pouch Length = 6"
- Quantity = 30



HCPCS A4419 | Filter (Filter covers included)

Coupling Size	Transparent	Opaque
● Small Stoma Opening	416776	416777
● Medium Stoma Opening	N/A	416781

ESTEEM SYNERGY® POUCH

Urostomy | Accuseal™ Tap with Valve

- Anti-Reflux Valve
- Pouch Length = 10" (*pouch length = 9")
- 2 Accuseal Night Drainage Adapters per box
- Quantity = 10



HCPCS A4432

Coupling Size	Transparent	Transparent*	Opaque
● Small Stoma Opening	405451	405446*	405448
● Medium Stoma Opening	405452	405447*	405449
● Large Stoma Opening	405453	N/A	405450

ONE-PIECE Systems

TWO BRANDS make up
ConvaTec's ONE-Piece Ostomy Systems



ESTEEM™

The Esteem™ + brand is an all-in-one solution combining the skin barrier and pouch in a single unit. Designed to be easy and convenient, it's one simple step to attach or remove. Available in multiple adhesive options, growing convexity options, and with ConvaTec's latest pouch technologies.



ActiveLife®

Active Life® is a legacy one-piece system brand that also combines the skin barrier and pouch in a single unit.

ESTEEM™ + POUCHES

 Recommended Product



 ESTEEM™+ POUCH

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Durahesive® Plus (*Extended Wear*)
- Acrylic Tape Collar
- Pouch Length = 12"
- Quantity = 10

HCPCS A5056 | Filter (Filter covers included)

Cutting Range	Transparent	Opaque
3/4 - 2 1/2" (19-64mm)	416976	416975



 ESTEEM™+ POUCH

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 12"
- Quantity = 10

HCPCS A5061 | No Filter

Cutting Range	Transparent	Opaque w/ Easy-View Window
13/16 - 2 3/4" (20-70mm)	416721	421853

HCPCS A4424 | Filter (Filter covers included)

Cutting Range	Transparent	Opaque w/ Easy-View Window	Opaque
13/16 - 2 3/4" (20-70mm)	416718	421858	416719

 Recommended Product

 ESTEEM™+ POUCH

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (Standard Wear)*
- Without Tape Collar
- Pouch Length = 14"
- Quantity = 10

HPCPS A5061 | No Filter
Cutting Range

Up to 4" (100mm)

Transparent

416908

*Double-layer adhesive (Stomahesive body-side, Durahesive top-side) which combines the skin friendly aspects of Stomahesive with the turtlenecking protection of Durahesive.


 ESTEEM™+ POUCH

Cut-to-Fit | Drainable | Clipped Closure

- Modified Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 12"
- 2 Tail Closures per box
- Quantity = 10

HPCPS A4424 | Filter (Filter covers included)*
Cutting Range

13/16 - 2 3/4" (20-70mm)

Transparent

421848


 ESTEEM™+ POUCH

Pre-Cut | Drainable | InvisiClose™ Tail Closure

- Modified Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 12"
- Quantity = 10

HPCPS A5061 | No Filter

Stoma Size	Opaque w/ Easy-View Window
1" (25mm)	421854
1 3/16" (30mm)	421855
1 3/8" (35mm)	421856
1 9/16" (40mm)	421857

HPCPS A4424 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque w/ Easy-View Window	Opaque
1" (25mm)	416724	421859	416725
1 3/16" (30mm)	416728	421860	416729
1 3/8" (35mm)	416732	421861	416733
1 9/16" (40mm)	416736	421862	416737


 ESTEEM™+ POUCH

Moldable | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 12"
- Quantity = 10

HPCPS A5061 | No Filter

Stoma Size	Transparent
13/16 - 1 3/16" (20-30mm)	413518
1 3/16 - 1 9/16" (30-40mm)	413519
1 9/16 - 1 15/16" (40-50mm)	413520

HPCPS A4424 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque w/ Easy-View Window
13/16 - 1 3/16" (20-30mm)	413515	413521
1 3/16 - 1 9/16" (30-40mm)	413516	413522
1 9/16 - 1 15/16" (40-50mm)	413517	413523



ESTEEM™+ POUCH

Pre-Cut | Drainable | Convex | InvisiClose™ Tail Closure

- Durahesive® (Extended Wear)
- Acrylic Tape Collar
- Convex Wafer Depth = 5mm
- Pouch Length = 12"
- Quantity = 10

HCPCS A5057 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque
3/4" (19mm)	416738	416739
7/8" (22mm)	416740	416741
1" (25mm)	416742	416743
1 1/8" (28mm)	416744	416745
1 1/4" (32mm)	416746	416747
1 3/8" (35mm)	416748	416749
1 1/2" (38mm)	416750	416751



ESTEEM™+ FLEX CONVEX POUCH

Cut-to-Fit | Drainable | Convex

- Esteem™ + Flex Convex Adhesive (Standard Wear)
- Without Tape Collar
- Convex Wafer Depth = 7mm
- Safe Seal Clipless Closure (1 Tail Closure per box)
- Pouch Length = 8 2/3"
- Includes Belt Tabs
- Quantity = 10

HCPCS A4389 | Filter (Filter covers included)


Cutting Range	Plateau Size	Opaque w/ Easy-View Window
3/4 - 1" (20-25mm)	V3	421622
3/4 - 1 3/8" (20-35mm)	V2	421619
3/4 - 1 2/3" (20-43mm)	V1	421615
3/4 - 2 6/25" (20-57mm)	V0	422159

Pre-Cut | Drainable | Convex

HCPCS A4389 | Filter (Filter covers included)

Stoma Size	Plateau Size	Opaque w/ Easy-View Window
1" (25mm)	V2	421620
1 1/6" (30mm)	V2	421621
1" (25mm)	V1	421616
1 1/6" (30mm)	V1	421617
1 3/8" (35mm)	V1	421618

Plateau Width Guidance for Esteem™+ Flex and Soft Convex

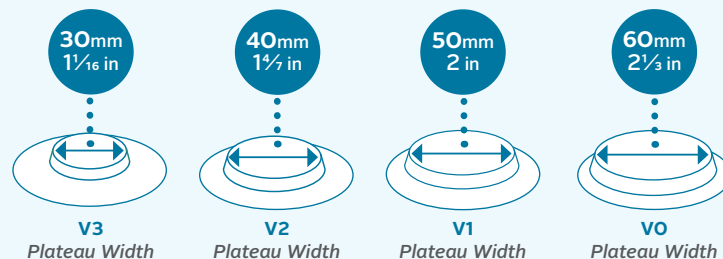

ESTEEM™+ SOFT CONVEX POUCH

Cut-to-Fit | Drainable | Convex

- Esteem™ + Soft Convex Adhesive (Standard Wear)
- Without Tape Collar
- Convex Wafer Depth = 4mm
- Safe Seal Clipless Closure (1 Tail Closure per box)
- Pouch Length = 8 2/3"
- Includes Belt Tabs
- Quantity = 10

HCPCS A4389 | Filter (Filter covers included)

Cutting Range	Plateau Size	Opaque w/ Easy-View Window
3/8 - 1 1/8" (10-28mm)	V3	422365
3/5 - 1 3/5" (15-40mm)	V2	422364
3/4 - 1 7/8" (20-47mm)	V1	422363

Pre-Cut | Drainable | Convex

HCPCS A4389 | Filter (Filter covers included)

Stoma Size	Plateau Size	Opaque w/ Easy-View Window
1" (25mm)	V2	422545
1 1/6" (30mm)	V2	422546
1" (25mm)	V1	422539
1 1/6" (30mm)	V1	422540
1 3/8" (35mm)	V1	422541



ESTEEM™+ POUCH

Cut-to-Fit | Closed-End

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 8"
- Quantity = 30

HCPCS A4416 | Filter (Filter covers included)

Cutting Range	Transparent	Opaque w/ Easy-View Window	Opaque
13/16 - 2 3/4" (20-70mm)	416700	421821	416701

ESTEEM™+ POUCH

Cut-to-Fit | Closed-End

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 6"
- Quantity = 30

HCPCS A4416 | Filter (Filter covers included)

Cutting Range	Opaque w/ Easy-View Window
13/16 - 2 3/4" (20-70mm)	421816




ESTEEM™+ POUCH

Pre-Cut | Closed-End

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 8"
- Quantity = 30

HCPCS A4416 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque w/ Easy-View Window	Opaque
1" (25mm)	416703	421687	416704
1 3/16" (30mm)	416709	421688	416707
1 3/8" (35mm)	416706	421689	416710
1 9/16" (40mm)	416712	421822	416713
2" (50mm)	N/A	421823	416715

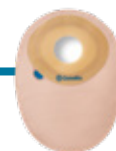

ESTEEM™+ POUCH

Moldable | Closed-End

- Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 8"
- Quantity = 30

HCPCS A4416 | Filter (Filter covers included)

Stoma Size	Opaque w/ Easy-View Window
13/16 - 1 3/16" (20-30mm)	413509
1 3/16 - 1 9/16" (30-40mm)	413510



 Recommended Product

 ESTEEM™+ POUCH

Urostomy | Cut-to-Fit | Accuseal™ Tap with Valve

- Durahesive® (Extended Wear)
- Acrylic Tape Collar
- Anti-Reflux Valve
- Pouch Length = 10"
- Includes Belt Tabs
- 2 Accuseal Night Drainage Adapters per box
- Quantity = 10

HCPCS A4428
Cutting Range
Transparent

1/2 - 1 3/4" (13-45mm)

421631


 ESTEEM™+ FLEX CONVEX POUCH

Urostomy | Cut-to-Fit | Convex

- Esteem™ + Flex Convex Adhesive (Standard Wear)
- Without Tape Collar
- Convex Wafer Depth = 7mm
- Soft Tap Closure
- Pouch Length = 7 2/3"
- Includes Belt Tabs
- 2 Esteem™ + Flex Convex Night Drainage Adapters per box
- Quantity = 10

HCPCS A4392
Cutting Range
Plateau Size
Opaque w/ Easy-View Window

3/8 - 1" (10-25mm)

V3

421628

3/8 - 1 3/8" (10-35mm)

V2

421626

3/8 - 1 2/3" (10-43mm)

V1

421623

Urostomy | Pre-Cut | Convex

HCPCS A4392
Stoma Size
Plateau Size
Opaque w/ Easy-View Window

9/16" (15mm)

V3

421629

9/16" (15mm)

V2

421627

9/16" (15mm)

V1

421624

4/5" (20mm)

V1

421625


 ESTEEM™+ SOFT CONVEX POUCH

Urostomy | Cut-to-Fit | Convex

- Esteem™ + Soft Convex Adhesive (Standard Wear)
- Without Tape Collar
- Convex Wafer Depth = 4mm
- Soft Tap Closure
- Pouch Length = 7 2/3"
- Includes Belt Tabs
- 2 Esteem™ + Soft Convex Night Drainage Adapters per box
- Quantity = 10

HCPCS A4392
Cutting Range
Plateau Size
Opaque w/ Easy-View Window

3/8 - 1 1/8" (10-28mm)

V3

422562

3/5 - 1 3/5" (15-40mm)

V2

422557

3/4 - 1 7/8" (20-47mm)

V1

422551

Urostomy | Pre-Cut | Convex

HCPCS A4392
Stoma Size
Plateau Size
Opaque w/ Easy-View Window

4/5" (20mm)

V3

422564

1" (25mm)

V3

422565

4/5" (20mm)

V2

422558

1" (25mm)

V2

422559

1 1/6" (30mm)

V2

422560

1" (25mm)

V1

422552

1 1/6" (30mm)

V1

422553

ACTIVELIFE® POUCHES

ACTIVELIFE® POUCH

Cut-to-Fit | Drainable | Clipped Closure

- Stomahesive® (Standard Wear)
- Acrylic Tape Collar
- Pouch Length = 12"
- 1 Tail Closure per box
- Quantity = See table below



HCPCS A5061 | No Filter

Cutting Range	Transparent Qty = 1	Transparent Qty = 10	Transparent Qty = 20
3/4 - 2 1/2" (19-64mm)	025433	022771	125330

ACTIVELIFE® POUCH

Cut-to-Fit | Drainable | Clipped Closure

- Durahesive® Plus (Extended Wear)
- Acrylic Tape Collar
- Pouch Length = 12"
- 1 Tail Closure per box
- Quantity = 10



HCPCS A4388 | No Filter

Cutting Range	Transparent	Opaque
3/4 - 2 1/2" (19-64mm)	400598	400599

ACTIVELIFE® POUCH

Pre-Cut | Drainable | Clipped Closure

- Stomahesive® (Standard Wear)
- Acrylic Tape Collar
- Pouch Length = 12" (*pouch length = 10")
- 1 Tail Closure per box
- Quantity = See table below



HCPCS A5061 | No Filter

Stoma Size	Transparent Qty = 10	Opaque Qty = 10	Opaque Qty = 20	Opaque* Qty = 10	Opaque* Qty = 20
3/4" (19mm)	022764	022757	125338	022750*	125331*
1" (25mm)	022765	022758	125339	022751*	N/A
1 1/4" (32mm)	022766	022759	125340	022752*	125333*
1 1/2" (38mm)	022767	022760	125341	022753*	125334*
1 3/4" (45mm)	022768	022761	125342	022754*	125335*
2" (50mm)	022769	022762	125343	022755*	125336*
2 1/2" (64mm)	022770	022763	125344	022756*	125337*

ACTIVELIFE® POUCH

Pre-Cut | Drainable | Convex | Clipped Closure

- Durahesive® (Extended Wear)
- Acrylic Tape Collar
- Pouch Length = 12"
- 1 Tail Closure per box
- Quantity = See table below



HCPCS A4390 | No Filter

Stoma Size	Transparent Qty = 5
3/4" (19mm)	175777
7/8" (22mm)	175778
1" (25mm)	175779
1 1/8" (28mm)	175780
1 1/4" (32mm)	175781
1 3/8" (35mm)	175782
1 1/2" (38mm)	175783
1 3/4" (45mm)	175784
2" (50mm)	175785

ACTIVELIFE® POUCH

**Cut-to-Fit | Stoma Cap**

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Quantity = 30

HCPCS A5055 | Filter (Filter covers included)

Cutting Range	Opaque
3/4 - 2" (19 - 50mm)	175611

ACTIVELIFE® POUCH

**Urostomy | Pre-Cut | Accuseal™ Tap with Valve**

- Durahesive® (Extended Wear)
- Acrylic Tape Collar
- Anti-Reflux Valve
- Pouch Length = 8"
- Includes Belt Tabs
- 2 Accuseal Night Drainage Adapters per box
- Quantity = 10

HCPCS A4428

Stoma Size	Transparent
3/4" (19mm)	650828
7/8" (22mm)	650829
1" (25mm)	650830
1 1/8" (28mm)	650831
1 1/4" (32mm)	650832
1 1/2" (38mm)	650833

ACTIVELIFE® POUCH

**Urostomy | Cut-to-Fit | Fold Over Tap**

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 9"
- Includes Belt Tabs
- 1 Universal Night Drainage Adapter per box
- Quantity = 10

HCPCS A5071

Cutting Range	Transparent
3/4 - 1 3/4" (19-45mm)	064927

ACTIVELIFE® POUCH

**Urostomy | Pre-Cut | Convex | Fold Over Tap**

- Durahesive® (Extended Wear)
- Acrylic Tape Collar
- Pouch Length = 8"
- Includes Belt Tabs
- 1 Universal Night Drainage Adapter per box
- Quantity = See table below

HCPCS A4393







Stoma Size	Transparent Qty = 5	Transparent Qty = 10
1/2" (13mm)	175790	125361
5/8" (16mm)	175791	125362
3/4" (19mm)	175792	125363
7/8" (22mm)	175793	125364
1" (25mm)	175794	125365
1 1/8" (28mm)	175795	125366
1 1/4" (32mm)	175796	125367
1 3/8" (35mm)	175797	125368
1 1/2" (38mm)	175798	125369

ACUTE KITS

ConvaTec's range of Acute Kits includes Post-Operative/Surgical kits in Sterile and Non-Sterile combinations, and Unit Dose kits.

Information on the components included in the Acute Kits can be found below the product information tables.

Color Key & Sizing Guide

COLOR MATCH						
FLANGE SIZE	1 1/4" (32mm)	1 1/2" (38mm)	1 3/4" (45mm)	2 1/4" (57mm)	2 3/4" (70mm)	4" (100mm)





POST-OPERATIVE/ SURGICAL KIT TWO-PIECE SYSTEM

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Double Layer Skin Barrier† (*Standard Wear*)
- Without Tape Collar
- Pouch Color: Transparent
- Pouch Length = 12" (**pouch length is 14"*)
- 1 Low-Pressure Adaptor
- Quantity = 5 Kits

HCPCS N/A

	Flange Size	Non-Sterile Kit	Sterile Kit	Ostomy Loop Rod
	1 3/4" (45mm)	416919	416920	None
	2 1/4" (57mm)	416923	416924	None
	2 3/4" (70mm)	416927	416928	65mm
	4" (100mm)	416929	416930	90mm

*Double-layer adhesive (*Stomahesive body-side, Durahesive top-side*) which combines the skin friendly aspects of Stomahesive with the turtlenecking protection of Durahesive.

Above kits include the following components:

Kit	Comparable to Skin Barrier in Kit* Stomahesive®	Comparable to Skin Barrier in Kit* Durahesive®	Pouch In Kit	Low Pressure Adaptor	Ostomy Loop Rod
416919/ 416920	401575	413155	416418	401993	None
416923/ 416924	401576	413156	416421	401994	None
416927/ 416928	401577	413157	416424	401995	022355
416929/ 416930	401906	401905	416472	401996	022356

*The double-layer skin barrier in the kit is not available in the retail market. Comparable skin barriers are in the chart above.

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Moldable | Drainable | InvisiClose™ Tail Closure

- Durahesive® (*Extended Wear*)
- Hydrocolloid Tape Collar (*White*)
- Pouch Color: Transparent
- Pouch Length = 12"
- 1 Low-Pressure Adaptor
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416917
●	2 1/4" (57mm)	416921
●	2 3/4" (70mm)	416925

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit	Low Pressure Adaptor
416917	411802	416418	401993
416921	411804	416421	401994
416925	411806	416424	401995

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Moldable | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (*Standard Wear*)
- Hydrocolloid Tape Collar (*White*)
- Pouch Color: Transparent
- Pouch Length = 12"
- 1 Low-Pressure Adaptor
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416918
●	2 1/4" (57mm)	416922
●	2 3/4" (70mm)	416926

Continued 

Natura™ Skin Barrier & Natura™+ Pouch Kit
Moldable | Drainable | InvisiClose™ Tail Closure

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit	Low Pressure Adaptor
416918	411803	416418	401993
416922	411805	416421	401994
416926	411807	416424	401995

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Moldable | Accordion Flange | Drainable | InvisiClose™ Tail Closure

- Durahesive® (*Extended Wear*)
- Hydrocolloid Tape Collar (*Tan*)
- Pouch Color: Transparent
- Pouch Length = 12"
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	2 1/4" (57mm)	421045
●	2 1/4" (57mm)	421046
●	2 3/4" (70mm)	421047

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit
421045	421039	416421
421046	421040	416421
421047	421041	416424

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Moldable | Accordion Flange | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (Standard Wear)
- Hydrocolloid Tape Collar (Tan)
- Pouch Color: Transparent
- Pouch Length = 12"
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	2 1/4" (57mm)	421042
●	2 1/4" (57mm)	421043
●	2 3/4" (70mm)	421044

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit
421042	421033	416421
421043	421034	416421
421044	421035	416424

NATURA™ SKIN BARRIER & SUR-FIT® NATURA™ + POUCH KIT

Moldable | Urostomy | Accuseal™ Tap with Valve

- Durahesive® (Extended Wear)
- Hydrocolloid Tape Collar (White)
- Anti-Reflux Valve
- Pouch Color: Transparent
- Pouch Length = 10"
- 1 Low Pressure Adaptor
- 1 Accuseal Night Drainage Adapter per box
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416932
●	2 1/4" (57mm)	416940
●	2 3/4" (70mm)	416948

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit	Low Pressure Adaptor
416932	411802	401544	401993
416940	411804	401545	401994
416948	411806	401546	401995

NATURA™ SKIN BARRIER & SUR-FIT® NATURA™ + POUCH KIT

Moldable | Urostomy | Accuseal™ Tap with Valve

- Stomahesive® (Standard Wear)
- Hydrocolloid Tape Collar (White)
- Anti-Reflux Valve
- Pouch Color: Transparent
- Pouch Length = 10"
- 1 Low Pressure Adaptor
- 1 Accuseal Night Drainage Adapter per box
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416934
●	2 1/4" (57mm)	416942
●	2 3/4" (70mm)	416950

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit	Low Pressure Adaptor
416934	411803	401544	401993
416942	411805	401545	401994
416950	411807	401546	401995

NATURA™ SKIN BARRIER & SUR-FIT® NATURA™ + POUCH KIT

Cut-to-Fit | Urostomy | Accuseal™ Tap with Valve

- Double Layer Skin Barrier[†] (Standard Wear)
- Acrylic Tape Collar (White)
- Anti-Reflux Valve
- Pouch Color: Transparent
- Pouch Length = 10"
- 1 Low Pressure Adaptor
- 1 Accuseal Night Drainage Adapter per box
- Quantity = 5 Kits

HCPCS N/A | Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416938
●	2 1/4" (57mm)	416946
●	2 3/4" (70mm)	416954

[†]Double-layer adhesive (Stomahesive body-side, Durahesive top-side) which combines the skin friendly aspects of Stomahesive with the turtlenecking protection of Durahesive.

Above kits include the following components:

Kit	Comparable to Skin Barrier in Kit* Stomahesive®	Comparable to Skin Barrier in Kit* Durahesive®	Pouch In Kit	Low Pressure Adaptor
416938	401575	413155	401544	401993
416946	401576	413156	401545	401994
416954	401577	413157	401546	401995

^{*}The double-layer skin barrier in the kit is not available in the retail market. Comparable skin barriers are in the chart above.

CONVATEC LOOP OSTOMY ROD

Sterile

- Quantity = 10

HCPCS A4421

	Flange Size	Size	Product Code
●	2 3/4" (70mm)	65mm	022355
●	4" (100mm)	90mm	022356



POST-OPERATIVE/ SURGICAL KIT ONE-PIECE SYSTEM

ESTEEM™ + POUCH KIT

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Double Layer Skin Barrier† (*Standard Wear*)
- Acrylic Tape Collar (*White*)
- Color: Transparent
- Pouch Length = 14"
- Quantity = 5 Kits

HCPCS N/A

Cutting Range	Non-Sterile Kit	Sterile Kit
Can be enlarged to accommodate stoma sizes up to 4" (100mm)	416915	416916

*Double-layer adhesive (*Stomahesive* body-side, *Durahesive* top-side) which combines the skin friendly aspects of *Stomahesive* with the turtlenecking protection of *Durahesive*.

UNIT-DOSE KITS

SUR-FIT® NATURA™ SKIN BARRIER & SUR-FIT® NATURA™ + POUCH KIT

Cut-to-Fit | Urostomy | Accuseal™ Tap with Valve

- Durahesive® (*Extended Wear*)
- Acrylic Tape Collar (*White*)
- Pouch Color: Transparent
- Pouch Length = 10"
- 1 Accuseal Night Drainage Adapter per box
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	401923
●	2 1/4" (57mm)	401924

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit
401923	413161	401544
401924	413162	401545

PEDIATRIC Products

LITTLE ONES is a pediatric Two-piece and One-piece system designed for babies smaller than 10 lbs, and children over 10 lbs.

The pouches have soft-welded edges and smaller pouch sizes designed to provide extra comfort for little ones.



Color Key & Sizing Guide

! When ordering a Pediatric Two-Piece System, **make sure you select the same flange/coupling size for the skin barrier and pouch.** Use the colored dots in the table to help.


COLOR MATCH	●	●	●	●
FLANGE/COUPLING SIZE	0 - 9/10" (0-23mm)	1/5 - 1 1/4" (4.8-31mm)	1 1/4" (32mm)*	1 3/4" (45mm)*

*Indicates Flange size measurement

LITTLE ONES[®] EXTRA SMALL TWO-PIECE OSTOMY SYSTEM

For babies smaller than 10lbs/4.5kg

LITTLE ONES[®] EXTRA SMALL SKIN BARRIER

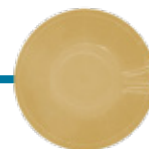

Cut-to-Fit

- Stomahesive[®] (Standard Wear)
- Adhesive Coupling Technology™
- No Starter Hole
- Quantity = 5

HCPCS A4414

	Cutting Range	Without Tape Collar
●	0 - 9/10" (0-23mm)	411632

LITTLE ONES[®] EXTRA SMALL POUCH

Drainable | InvisiClose™ Tail Closure

- Adhesive Coupling Technology™
- Pouch Length = 6"
- Quantity = 10



HCPCS A5063 | No Filter

	Stoma Size	Transparent
●	0 - 9/10" (0-23mm)	411631

LITTLE ONES® EXTRA SMALL ONE-PIECE OSTOMY SYSTEM

For babies smaller than 10lbs/4.5kg

LITTLE ONES® EXTRA SMALL POUCH

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (Standard Wear)
- Without Tape Collar
- No Starter Hole
- Pouch Length = 6"
- Quantity = 10



HCPCS A5061 | No Filter

Stoma Size	Transparent
0 - 9/10" (0-23mm)	411633

LITTLE ONES® TWO-PIECE OSTOMY SYSTEM

Suggested for children over 10lbs/4.5kg

LITTLE ONES® SKIN BARRIER

Cut-to-Fit

- Stomahesive® (Standard Wear)
- Adhesive Coupling Technology™
- Quantity = 5



HCPCS A4414

Cutting Range	Hydrocolloid Tape Collar WHITE	No Tape Collar
1/5 - 1 1/4" (4.8-31mm)	411642	411641

LITTLE ONES® POUCH

Drainable | InvisiClose™ Tail Closure

- Adhesive Coupling Technology™
- Pouch Length = 8"
- Quantity = 10



HCPCS A4425 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque
1/5 - 1 1/4" (4.8-31 mm)	411637	411638

LITTLE ONES® POUCH

Closed-End

- Adhesive Coupling Technology™
- Pouch Length = 6"
- Quantity = 10



HCPCS A4419 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque
1/5 - 1 1/4" (4.8-31 mm)	411640	411639

LITTLE ONES® SKIN BARRIER

Cut-to-Fit

- Stomahesive® (Standard Wear)
- Flange System
- Quantity = 5



HCPCS A4414

Flange Size	Cutting Range	Acrylic Tape Collar WHITE
1 1/4" (32mm)	1/2 - 3/4" (13-19mm)	401925
1 3/4" (45mm)	1 1/8 - 1 1/4" (28-32mm)	401926

LITTLE ONES® POUCH

Drainable | Clipped Closure

- Stomahesive® (Standard Wear)
- Flange System
- Pouch Length = 6"
- 10 Soft, Moldable Tail Closures per box (Designed for Pediatric Use)
- Quantity = 10




HCPCS A5063 | No Filter

	Flange Size	Transparent
●	1 1/4" (32mm)	401927
●	1 3/4" (45mm)	401928

LITTLE ONES® POUCH

Closed-End

- Flange System
- Pouch Length = 5"
- Quantity = 20

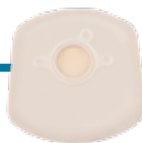

HCPCS A5054 | No Filter

	Flange Size	Opaque
●	1 1/4" (32mm)	401930
●	1 3/4" (45mm)	401931

LITTLE ONES® POUCH

Urostomy | Fold Over Tap

- Flange System
- Pouch Length = 6"
- 2 Universal Night Drainage Adapters per box
- Quantity = 10



HCPCS A5073 | No Filter

	Flange Size	Transparent
●	1 1/4" (32mm)	401929

LITTLE ONES®
**ONE-PIECE
OSTOMY SYSTEM**

Suggested for children over 10lbs/4.5kg

LITTLE ONES® POUCH

Cut-to-Fit | Drainable | Clipped Closure

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 6"
- 15 Soft, Moldable Tail Closures per box (Designed for Pediatric Use)
- Quantity = 15



HCPCS A5061 | No Filter

Cutting Range	Transparent
5/16 - 2" (8-50mm)	020922

LITTLE ONES® POUCH

Urostomy | Cut-to-Fit | Fold Over Tap

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 5"
- 2 Universal Night Drainage Adapters per box
- Quantity = 15



HCPCS A5071

Cutting Range	Transparent
5/16 - 1" (8-25mm)	020917

ACCESSORIES



SKIN PREPARATION

Solutions that help promote healthy skin by reducing irritation and damage.

Choosing the right accessories

Whether you are looking for better skin protection or added security to your ostomy system, with the right combination of accessories, you can feel even more comfortable and at ease.



**Skin
Preparation**



**Pouching
Support**



**Added
Security**



**Pouch
Removal**



SENSI-CARE™

Sting-Free Barrier Wipe

- Designed to create a protective layer that shields skin from damage
- Helps reduce skin stripping and irritation from body waste
- May help increase the contact or adhesion to skin
- Silicone based
- Alcohol free
- CHG Compatible¹⁷

HCPCS A5120

Liquid per Wipe	Product Code	Box of	Boxes per Case*
3 mL	413501	30	20

*Sold by the case

SENSI-CARE™

Sting-Free Barrier Spray

- Designed to create a protective layer that shields skin from damage
- Helps reduce skin stripping and irritation from body waste
- May help increase the contact or adhesion to skin
- Silicone based
- Alcohol free
- Non-Aerosol
- CHG Compatible¹⁷



HCPCS A4369

Spray Bottle Size	Product Code	Box of	Boxes per Case*
28 mL	420797	1	12
50 mL	413502	1	12

*Sold by the case

SENSI-CARE™

Sting-Free Barrier Foam Applicator

- Designed to create a protective layer that shields skin from damage
- Helps reduce skin stripping and irritation from body waste
- May help increase the contact or adhesion to skin
- Silicone based
- Alcohol free
- CHG Compatible¹⁷



HCPCS A5120 | Sterile

Applicator Size	Product Code	Individually Packaged Applicators per Box	Boxes per Case*
1 mL	420793	5	20
1 mL	420794	25	4
3 mL	420795	5	20
3 mL	420796	25	4

*Sold by the case

ALLKARE®

Protective Barrier Wipe

- Designed to create a protective layer that shields skin from damage
- Helps reduce skin stripping and irritation from body waste
- May help increase the contact or adhesion to skin
- Contains Alcohol



HCPCS A5120

Product Code	Box of
037439	50
037444	100

STOMAHESIVE®

Protective Powder

- Helps protect irritated skin by absorbing excess moisture
- For use on weeping skin around a stoma
- Lightly dust affected area with powder, wiping off excess. May be followed up with a skin barrier wipe or spray



HCPCS A4371

Size	Product Code	Box of
1 oz.	025510	1 bottle



POUCHING SUPPORT

Products that help create the best fit around the stoma, prevent leakage and help level out the application surface area for better pouch performance.

EAKIN COHESIVE®

Seals

- Moldable barrier that helps seal the area around hard-to-fit stomas from any leakage under the skin barrier
- Acts as a filler on uneven skin surfaces to help create a flatter area for the skin barrier to adhere to
- Designed to absorb moisture and prevent it coming in contact with, and causing harm to the skin—forms a gel that protects the skin



HCPCS A4385

Size	Product Code	Depth	Box of
Small 2"	839002	4.2 mm	20
Slims 2"	839005	3 mm	10
Large 4"	839001	2.1 mm	10
StomaWrap 3 1/3"	839006	3 mm	10

STOMAHESIVE®

Paste

- Helps to protect exposed skin between the base of the stoma and the opening of the skin barrier. Can also be used to fill crevices, skin folds and scars to help create a flat pouching surface for optimal appliance attachment
- Helps to prevent leakage under the skin barrier that may result in skin breakdown
- Pectin-based
- Contains alcohol



HCPCS A4406

Size	Product Code	Box of
2 oz.	183910	1 tube

EAKIN COHESIVE®

Paste

- Use to fill crevices, skin folds and scars to help create an even pouching surface
- Leakage prevention by helping form a waterproof barrier around the stoma
- Alcohol free



HCPCS A4405

Size	Product Code	Box of
2.1 oz.	839010	1 tube

STOMAHESIVE®

Skin Barrier

- Helps fill uneven skin surfaces to create a flat pouching surface for optimal appliance attachment
- Designed to help protect the skin from contact with bodily fluids
- Large cutting surface to accommodate large and swollen stomas
- Sterile Stomahesive® skin barrier can be used post-operatively



HCPCS A4362

Size	Product Code	Box of
4 x 4"	021712	5
4 x 4 (Sterile)	021901	5

HCPCS A5122

Size	Product Code	Box of
8 x 8"	021715	3

EAKIN COHESIVE®

Skin Barrier

- Designed to help protect the skin from contact with body fluids
- Acts as a filler on uneven skin surfaces to help create a flat pouching surface for optimal appliance attachment
- Moldable and moisture absorbing



HCPCS A4385

Size	Product Code	Box of
4 x 4"	839004	5
4 x 8"	839003	5

STOMAHESIVE®

Strips

- Helps to create a secure seal between the stoma and the appliance, to protect the skin against stoma output
- Can be cut, torn and shaped into a variety of shapes to help prevent leakage around the stoma by acting as a filler
- Designed to absorb moisture from the stoma, helping to create a snug fit around the stoma.



HCPCS A4362

Size	Product Code	Box of
3 cards of 5 strips each	025542	15 strips

SUR-FIT® NATURA™

Disposable Convex Inserts

- Convex inserts that can be added to a flat skin barrier to create a shallow convexity—insert into the inner diameter of the plastic flange on the pouch facing side of a SUR-FIT® Natura™ Durahesive® and Stomahesive® skin barrier
- Select an insert that is at least ¼" larger than the maximum diameter of the stoma
- For SUR-FIT® Natura™ skin barriers only. Not for use with ConvaTec Moldable Technology™ skin barriers, or any other product
- Quantity = 5



HCPCS A5093

	For Use With Flange Size	Stoma Size no Wider Than	Convex Insert Opening Size	Product Code
●	1 1/2" (38mm)	1/2"	3/4" (19mm)	404006
●	1 1/2" (38mm)	3/4"	1" (25mm)	404008
●	1 3/4" (45mm)	7/8"	1 1/8" (29mm)	404009
●	1 3/4" (45mm)	1"	1 1/4" (32mm)	404010
●	1 3/4" (45mm)	1 1/8"	1 3/8" (35mm)	404011
●	2 1/4" (57mm)	1 1/4"	1 1/2" (38mm)	404012
●	2 1/4" (57mm)	1 3/8"	1 5/8" (41mm)	404013



ADDED SECURITY

An extra layer of protection to keep your ostomy system in place.

EASE™



Strips

- Helps keep your skin barrier in place
- Can be applied to cover the edges of the skin barrier
- Designed to be thin and flexible
- Water-resistant

HCPCS A4362 (4 strips = 1 billable unit of A4362)

Size	Product Code	Quantity
1 1/5 x 3 1/2" (3 x 9cm)	422163	80 strips

CONVATEC®



Ostomy Appliance Belt

- Used to help keep your ostomy system in place

HCPCS A4367

Size	Product Code	Box of
Adjustable up to 42"	175507	1

DIAMONDS™



Gelling and Odor Control Sachets

- Gels liquid contents
- Helps reduce or eliminate excess gas that can lead to pouch ballooning
- Contains activated charcoal which helps to eliminate odors

HCPCS A4422

Product Code	Box of
TR105	100 Sachets

DUOLOCK™



Curved Tail Closure

- Helps secure the end of ConvaTec® Drainable pouch

HCPCS A4363

Product Code	Box of
175652	10

CONVATEC®



Tail Closure

- Helps secure the end of ConvaTec® Drainable pouch

HCPCS A4363

Product Code	Box of
175635	10

CONVATEC®



Night Drainage Container Set*

- Plastic vessel to use for overnight urine collection
- Container capacity of 2000 cc/ml
- Reusable, can be washed out between uses

HCPCS A5102

Product Code	Quantity
027060	1 set

*Contains: 1 ConvaTec® Night Drainage Container with tubing, 1 ConvaTec® Night Drainage Container Cover and 1 Universal Adapter

CONVATEC®



Night Drainage Container Tubing

- Replacement tubing for the ConvaTec® Night Drainage Container

HCPCS A4331

Product Code	Box of
027062	1



POUCH REMOVAL

Solutions that help promote healthy skin by protecting it from irritation and damage caused by adhesives.

SENSI-CARE™



Sting-Free Adhesive Remover Wipe

- Designed to remove adhesive residue from the skin
- Usage with each appliance change helps prevent residue build-up
- Silicone based
- Alcohol free
- CHG Compatible¹⁷

HCPCS A4456

Liquid per Wipe	Product Code	Box of	Boxes per Case*
3 mL	413500	30	20

*Sold by the case

SENSI-CARE™



Sting-Free Adhesive Releaser Spray

- Tracks between adhesive & skin, gently releasing the adhesive, tape or appliance
- Helps minimize painful skin stripping on removal¹⁸
- Silicone based
- Alcohol free
- Non-Aerosol
- CHG Compatible¹⁷

HCPCS A4455

Spray Bottle Size	Product Code	Box of	Boxes per Case*
50 mL	413499	1	12
150 mL	420798 (Sterile)	1	12

*Sold by the case

ALLKARE®



Adhesive Remover Wipe

- Helps gently remove adhesives, including skin barriers, tapes and hydrocolloid dressings
- Contains alcohol

HCPCS A4456

Product Code	Box of
037436	50
037443	100

IRRIGATION

VISI-FLOW®

Irrigator with Stoma Cone

HCPCS A4398 & A4399

Product Code	Box of
401989	1



VISI-FLOW®

Irrigation Sleeve

- Compatible with Natura™ and SUR-FIT® Natura™ skin barriers
- Color: Transparent
- 1 Tail Closure per box
- Quantity = 5

HCPCS A4397

	Flange Size	Product Code
●	1 1/2" (38mm)	401911
●	1 3/4" (45mm)	401912
●	2 1/4" (57mm)	401913
●	2 3/4" (70mm)	401914



VISI-FLOW®

Stoma Cone

HCPCS A4399

Product Code	Box of
022736	1



VISI-FLOW®

Irrigation Sleeve Tail Closure

HCPCS A4363

Product Code	Box of
175650	10



VISI-FLOW®

Irrigation Adaptor Faceplate

- Compatible with Natura™ and SUR-FIT® Natura™ skin barriers
- Quantity = 1

HCPCS A4361

	Flange Size	Product Code
●	1 3/4" (45mm)	401918
●	2 3/4" (70mm)	401919



eakin® WOUND & FISTULA Pouches

Eakin wound pouches can help contain fistula and handle exuding wounds. Pouch options include the ability to connect tubing for remote drainage, or a fold & tuck option.



EAKIN®

Wound Pouch (1.2 x 1.8" Cutting Area)

- Quantity = 10

HCPCS A6154

Cutting Range

Fold and Tuck Closure

1.2 x 1.8" (3 x 4.5cm)

839250



EAKIN®

Wound Pouch (3.14 x 2.4" Cutting Area)

- Quantity = 10

HCPCS A6154

Cutting Range

Remote Drainage Attachment and Tap Closure

3.14 x 2.4" (8 x 6.1cm)

839267



EAKIN®

Wound Pouch (3.0 x 4.3" Cutting Area)

- Quantity = 10

HCPCS A6154

Cutting Range

Remote Drainage Attachment and Tap Closure

Fold and Tuck Closure

3.0 x 4.3" (7.6 x 10.9cm)

839261

839251



EAKIN®

Wound Pouch (4.3 x 6.9" Cutting Area)

- Quantity = 10

HCPCS A6154

Cutting Range

Remote Drainage Attachment and Tap Closure

Fold and Tuck Closure

4.3 x 6.9" (10.9 x 17.5cm)

839262

839252



EAKIN®

Horizontal Cutting Wound Pouch

(9.7 x 6.3" Cutting Area)

- Quantity = 5



HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure with Access Window	Fold and Tuck Closure
9.7 x 6.3" (24.6 x 16cm)	839264	839253

EAKIN®

Vertical Cutting Wound Pouch

(6.3 x 9.7" Cutting Area)

- Quantity = 5



HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure	Remote Drainage Attachment and Tap Closure with Access Window
6.3 x 9.7" (16 x 24.6cm)	839265	839268

EAKIN®

Vertical Cutting Wound Pouch

(5.1 x 11.4" Cutting Area)

- Quantity = 5



HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure
5.1 x 11.4" (13 x 29cm)	839266

EAKIN®

Eakin Access Window

- Add to Eakin wound and fistula pouches to allow instant access to the wound (*may be used with 839252, 839253, 839262, 839263, 839265, 839266*)
- Quantity = 5



HCPCS N/A

Size	Remote Drainage Attachment and Tap Closure
N/A	839280


DISCONTINUED PRODUCTS

SPECIFIC PRODUCT NUMBERS HAVE CHANGED[†]

We've made changes to our Product Portfolio

- Some legacy pouches have been upgraded to our newest technology
- We have provided a cross reference (*see the following tables*) to map the discontinued products with its best replacement





Use this guide to identify the product numbers being upgraded to + pouches or being discontinued.

A  sign indicates these products feature + pouch technology. Other modifications noted in the tables include: ConvaTec filter technology, package quantities, pouch lengths, color, adhesive, larger cutting areas and pre-cut holes.

Natura™ + Color Key & Sizing Guide

COLOR MATCH						
FLANGE SIZE	1 1/4" (32mm)	1 1/2" (38mm)	1 3/4" (45mm)	2 1/4" (57mm)	2 3/4" (70mm)	4" (100mm)

Esteem Synergy®+ Color Key & Sizing Guide

COLOR MATCH				
COUPLING SIZE	Small Stoma Opening	Medium Stoma Opening	Large Stoma Opening	Extra-Large Stoma Opening

[†] Only product numbers listed are being converted to + Pouches. If an item is not listed, it is not changing.

Previous SKU #	Previous QTY	New SKU #	New QTY	Flange or Stoma Size	Differences	New SKU Page #
125352	10	175777	5	19mm	Pack size change	56
125353	10	175778	5	22mm	Pack size change	56
125354	10	175779	5	25mm	Pack size change	56
125355	10	175780	5	28mm	Pack size change	56
125356	10	175781	5	32mm	Pack size change	56
125357	10	175782	5	35mm	Pack size change	56
125358	10	175783	5	35mm	Pack size change	56
125359	10	175784	5	45mm	Pack size change	56
175761	15	416700	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	50
175762	15	416703	30	25mm	⊕ Modified Stomahesive®	52
175763	15	416706	30	35mm	⊕ 3mm Larger Pre-cut Hole Modified Stomahesive®	52
175764	15	416712	30	40mm	⊕ 2mm Larger Pre-cut Hole Modified Stomahesive®	52
175765	15	421823	30	50mm	⊕ 5mm Larger Pre-cut Hole Opaque Modified Stomahesive®	52
175766	15	421823	30	50mm	⊕ Opaque Modified Stomahesive®	52
175768	15	421821	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	50
175769	15	421687	30	25mm	⊕ Modified Stomahesive®	52
175770	15	421689	30	30mm	⊕ 3mm Larger Pre-cut Hole Modified Stomahesive®	52
175771	15	421822	30	40mm	⊕ 2mm Larger Pre-cut Hole Modified Stomahesive®	52
175772	15	421823	30	50mm	⊕ 5mm Larger Pre-cut Hole Modified Stomahesive®	52
175773	15	421823	30	50mm	⊕ Modified Stomahesive®	52
401521	30	421679	30	● 38mm flange	⊕	37
401525	30	421678	30	● 38mm flange	⊕	37
401529	30	421677	30	● 32mm flange	⊕ Length increase from 5" to 8"	37
401530	20	421679	30	● 38mm flange	⊕ Length increase from 5" to 8"	37
401531	20	421796	30	● 45mm flange	⊕ Length increase from 5" to 6" Filter HCPCS now A4419	38
401532	20	421797	30	● 57mm flange	⊕ Length increase from 5" to 6" Filter HCPCS now A4419	38
401916	1	See Notes		● 45mm flange	Irrigation components available separately	
401917	1	See Notes		● 70mm flange	Irrigation components available separately	
401920	5	416919	5	● 45mm flange	Kit upgraded to our improved InvisiClose closure system	60
401921	5	416927	5	● 70mm flange	Kit upgraded to our improved InvisiClose closure system	60
401922	5	416929	5	● 100mm flange	Kit upgraded to our improved InvisiClose closure system	60
401993	10	None		● 45mm flange	Refer to our line of Accordion products for low-pressure support	
401994	10	None		● 57mm flange	Refer to our line of Accordion products for low-pressure support	
401995	10	None		● 70mm flange	Refer to our line of Accordion products for low-pressure support	
401996	10	None		● 100mm flange	Discontinued	
402239	5	416923	5	● 57mm flange	Kit upgraded to our improved InvisiClose closure system	60
404007	5	413180	10	● 38mm flange	Look at our extended convex line	30
405406	10	See Notes		● Extra Large	Discontinued; Try 100mm Natura™ pouch and wafer	
405483	10	See Notes		● Extra Large	Discontinued; Try 100mm Natura™ pouch and wafer	
409275	5	416929	5	● 100mm flange	Kit upgraded to our improved InvisiClose closure system	60
411309	20	416415	10	● 38mm flange	⊕	34
411312	20	416423	10	● 70mm flange	⊕	34
413141	60	421821	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	50

Previous SKU #	Previous QTY	New SKU #	New QTY	Flange or Stoma Size	Differences	New SKU Page #
413142	60	421687	30	25mm	⊕ Modified Stomahesive®	52
413143	60	421689	30	35mm	⊕ 3mm Larger Pre-cut Hole Modified Stomahesive®	52
413144	60	421822	30	40mm	⊕ 2mm Larger Pre-cut Hole Modified Stomahesive®	52
413145	60	421823	30	50mm	⊕ 5mm Larger Pre-cut Hole Modified Stomahesive®	52
413146	60	421823	30	50mm	⊕ Modified Stomahesive®	52
413147	60	416700	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	50
413148	60	416703	30	25mm	⊕ Modified Stomahesive®	52
413149	60	416706	30	35mm	⊕ 3mm Larger Pre-cut Hole Modified Stomahesive®	52
413150	60	416712	30	40mm	⊕ 2mm Larger Pre-cut Hole Modified Stomahesive®	52
413151	60	421823	30	50mm	⊕ 5mm Larger Pre-cut Hole Opaque Modified Stomahesive®	52
413152	60	421823	30	50mm	⊕ Opaque Modified Stomahesive®	52
413169	60	421679	30	● 38mm flange	⊕	37
413172	60	421893	30	● 70mm flange	⊕	37
413173	60	421678	30	● 38mm flange	⊕	37
413176	60	421798	30	● 70mm flange	⊕	37
413335	20	421739	10	● 57mm flange	⊕	34
413336	20	421740	10	● 70mm flange	⊕	34
413435	10	421074	10	● 32mm flange	Length decreases from 10" to 9"	40
413436	10	None		● 38mm flange	Discontinued	
413511	30	None		40-50mm	Discontinued	
416400	30	421676	30	● 32mm flange	⊕ Enhanced with see-through window	37
416401	30	421677	30	● 32mm flange	⊕ Enhanced with see-through window	37
416402	30	421677	30	● 32mm flange	⊕ Length increase 6" to 8" Now with a filter	37
416404	30	421679	30	● 38mm flange	⊕ Enhanced with see-through window	37
416405	30	421679	30	● 38mm flange	⊕ Length increase 6" to 8" Now with a filter	37
416407	30	421681	30	● 45mm flange	⊕ Enhanced with see-through window	37
416408	30	421796	30	● 45mm flange	⊕ Now with a filter HCPCS now A4419	38
416410	30	421683	30	● 57mm flange	⊕ Enhanced with see-through window	37
416411	30	421797	30	● 57mm flange	⊕ Now with a filter HCPCS now A4419	38
416413	30	421893	30	● 70mm flange	⊕ Enhanced with see-through window	37
416702	30	421816	30	20-70mm	⊕	50
416780	30	416781	30	● Medium	Opaque	44
416797	10	416795	10	● Small	Now with a filter HCPCS now A4425	43
416900	10	416798	10	● Medium	Now with a filter HCPCS now A4425	43
421075	10	None		● 38mm flange	Discontinued	
421076	10	413437	10	● 45mm flange	Length increases from 9" to 10"	40
421612	30	None		20-25mm	Discontinued	
421613	30	None		20-35mm	Discontinued	
421614	30	None		20-43mm	Discontinued	
650408	30	421821	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	50
650422	30	416700	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive® No tape collar	50
839260	10	839267	10	4.3x3"	Slightly larger version	88
839263	5	839265	5	9.7x6.3"	Vertical orientation	89

OSTOMY RESOURCES AND SUPPORT

for Your Patient

ConvaTec

ConvaTec is a global medical products and technologies company developing ostomy products and solutions for almost 40 years.

Toll Free: 1-800-422-8811
www.convatec.com/ostomy

Ostomysecrets®

A ConvaTec company, offers a full line of stylish and functional undergarments, swim wear and accessories designed to support your ostomy solution.

Toll Free: 1-877-613-6246
www.ostomysecrets.com

OTHER OSTOMY RESOURCES

The United Ostomy Associations of America (UOAA)

National organization providing support, information and advocacy.

Toll Free: 1-800-826-0826
www.ostomy.org

American Cancer Society

Nationwide community-based, voluntary healthcare organization.

Toll Free: 1-800-ACS-2345
www.cancer.org

me+™

Support, insights & product information for those living with an ostomy. Enrolling in me+ is fast, easy and free.

Toll Free: 1-800-422-8811
e-mail: CIC@convatec.com
www.mepluscare.com

Wound Ostomy Continence Nurses Society™ (WOCN®)

Professional nursing society dedicated to providing expert health care to individuals with wound, ostomy and continence concerns.

Toll Free: 1-888-224-9626
www.wocn.org

Crohn's & Colitis Foundation of America

Non-profit, volunteer-driven organization dedicated to finding a cure for Crohn's Disease and ulcerative colitis.

Toll Free: 1-800-932-2423
www.ccfa.org

HCPCS	HCPCS Descriptors	Medicare Usual Maximum Qty.
A4331	Extension Drainage Tube	1 per 6 months
A4361	Ostomy Faceplate	3 per 6 months
A4362	Skin Barrier, Solid, 4x4 or Equivalent	20 per month
A4363	Ostomy Clamp, Any Type, Replacement Only	As needed
A4367	Ostomy Belt	1 per month
A4369	Ostomy Skin Barrier Liquid (spray or brush)	2 oz per month
A4371	Ostomy Skin Barrier Powder	10 oz per 6 months
A4385	Skin Barrier, Solid 4x4 or Equivalent, Extended Wear	As needed
A4388	Pouch, Drainable, w/Extended Wear Barrier Attached	As needed
A4389	Ostomy Pouch, Drainable, w/ Barrier Attached, w/ Built-In Convexity, w/ Filter (1-piece)	20 per month
A4390	Pouch, Drainable, w/Extended Wear Barrier, Attached, w/convexity	As needed
A4392	Ostomy Pouch, Urinary, w/ Barrier Attached, w/ Built-In Convexity (1-piece)	20 per month
A4393	Ostomy Pouch, Urinary, w/ Extended Wear Barrier Attached, w/ Built-In Convexity (1-piece)	20 per month
A4397	Irrigation Sleeve	4 per month
A4398	Ostomy Irrigation Bag	2 per 6 months
A4399	Ostomy Irrigation Cone/Catheter, and Brush	2 per 6 months
A4402	Lubricant	4 oz per month
A4405	Ostomy Skin Barrier, Non-Pectin Based, Paste	4 oz per month
A4406	Ostomy Skin Barrier, Pectin-Based, Paste	4 oz per month
A4407	Skin Barrier, w/Flange, Extended Wear w/ Convexity	As needed
A4409	Skin Barrier, w/Flange, Extended Wear <= 4x4	As needed
A4410	Skin Barrier, w/Flange, Extended Wear > 4x4	As needed
A4413	Pouch, Drainable, High Output, for Barrier w/Flange	As needed
A4414	Skin Barrier, w/Flange, <= 4x4	20 per month
A4415	Skin Barrier, w/Flange, > 4x4	20 per month
A4416	Pouch, closed, w/Barrier Attached, w/Filter	60 per month
A4417	Ostomy Pouch, Closed, w/ Barrier Attached, w/ Built-In Convexity, w/ Filter (1-piece)	60 per month
A4419	Pouch, closed for use on Barrier w/Flange, w/Filter	60 per month
A4421	Ostomy Supply, Miscellaneous	As needed
A4422	Ostomy Absorbent Material (Sheet/Pad/Crystal Packet)	No Fee Schedule
A4424	Ostomy Pouch, Drainable, w/ Barrier Attached, w/ Filter (1-piece)	20 per month
A4425	Pouch, Drainable, with filter for use w/Flange	20 per month
A4428	Pouch, Urinary, w/ Extended Wear Barrier w/Faucet-type Tap	As needed
A4432	Pouch, Urinary, for use on Barrier w/Flange, w/Faucet-type Tap	20 per month
A4455	Adhesive Remover or Solvent	2 oz per month
A4456	Adhesive Remover Wipes, Any Type	50 per month
A5054	Pouch Closed, for Barrier w/Flange	60 per month
A5055	Stoma Cap, each	31 per month
A5056	Pouch, Drainable, w/Extended Wear Barrier Attached, w/filter	40 per month
A5057	Pouch, Drainable, w/Extended Wear Barrier, w/convexity, w/filter	40 per month
A5061	Pouch, Drainable, w/Barrier Attached	20 per month
A5063	Pouch, Drainable, for use w/Flange	20 per month
A5071	Pouch, Urinary, w/Barrier Attached	20 per month
A5073	Pouch, Urinary, for Barrier w/Flange	20 per month
A5093	Ostomy Accessory, Convex Insert	10 per month
A5102	Bedside Drainage Bottle w/ or wo/ tubing, rigid or expandable	2 per 6 months
A5120	Skin Barrier, Wipes, each	150 per 6 months
A5122	Skin Barrier, Solid, 8 X 8 or Equivalent	20 per month
A6154	Wound Pouch	As needed

020917.....	74	125352.....	91	401525.....	91	404007.....	91
020922.....	74	125353.....	91	401526.....	38	404008.....	80
021712.....	79	125354.....	91	401527.....	38	404009.....	80
021715.....	79	125355.....	91	401528.....	38	404010.....	80
021901.....	79	125356.....	91	401529.....	91	404011.....	80
022355.....	66	125357.....	91	401530.....	91	404012.....	80
022356.....	66	125358.....	91	401531.....	91	404013.....	80
022736.....	85	125359.....	91	401532.....	91	404014.....	36
022750.....	56	125361.....	58	401533.....	40	404015.....	36
022751.....	56	125362.....	58	401534.....	40	404016.....	36
022752.....	56	125363.....	58	401535.....	40	404017.....	36
022753.....	56	125364.....	58	401536.....	40	404018.....	36
022754.....	56	125365.....	58	401537.....	40	404025.....	36
022755.....	56	125366.....	58	401538.....	40	404026.....	36
022756.....	56	125367.....	58	401539.....	40	404027.....	36
022757.....	56	125368.....	58	401540.....	40	404028.....	36
022758.....	56	125369.....	58	401541.....	40	404029.....	36
022759.....	56	175507.....	81	401542.....	39	404030.....	35
022760.....	56	175611.....	57	401543.....	39	404031.....	35
022761.....	56	175635.....	82	401544.....	39	404032.....	35
022762.....	56	175650.....	86	401545.....	39	404033.....	35
022763.....	56	175652.....	82	401546.....	39	404034.....	35
022764.....	56	175761.....	91	401547.....	39	404592.....	30
022765.....	56	175762.....	91	401548.....	39	404593.....	30
022766.....	56	175763.....	91	401549.....	39	404594.....	30
022767.....	56	175764.....	91	401550.....	39	405406.....	91
022768.....	56	175765.....	91	401551.....	39	405446.....	44
022769.....	56	175766.....	91	401552.....	39	405447.....	44
022770.....	56	175768.....	91	401553.....	39	405448.....	44
022771.....	55	175769.....	91	401554.....	39	405449.....	44
025433.....	55	175770.....	91	401555.....	39	405450.....	44
025510.....	77	175771.....	91	401557.....	37	405451.....	44
025542.....	80	175772.....	91	401558.....	37	405452.....	44
027060.....	82	175773.....	91	401559.....	37	405453.....	44
027062.....	82	175777.....	56	401573.....	31	405456.....	41
037436.....	84	175778.....	56	401574.....	31	405457.....	41
037439.....	77	175779.....	56	401575.....	31	405458.....	41
037443.....	84	175780.....	56	401576.....	31	405467.....	41
037444.....	77	175781.....	56	401577.....	31	405468.....	41
064927.....	58	175782.....	56	401905.....	27	405473.....	42
125257.....	31	175783.....	56	401906.....	31	405474.....	42
125258.....	31	175784.....	56	401909.....	39	405475.....	42
125259.....	31	175785.....	56	401911.....	85	405476.....	42
125260.....	31	175790.....	58	401912.....	85	405477.....	42
125261.....	31	175791.....	58	401913.....	85	405478.....	42
125262.....	31	175792.....	58	401914.....	85	405479.....	42
125263.....	31	175793.....	58	401916.....	91	405480.....	42
125264.....	31	175794.....	58	401917.....	91	405481.....	42
125265.....	31	175795.....	58	401918.....	86	405482.....	42
125266.....	31	175796.....	58	401919.....	86	405483.....	41
125267.....	32	175797.....	58	401920.....	91	409268.....	42
125268.....	32	175798.....	58	401921.....	91	409269.....	42
125269.....	32	183910.....	78	401922.....	91	409270.....	42
125270.....	32	400598.....	55	401923.....	68	409275.....	91
125271.....	32	400599.....	55	401924.....	68	411264.....	36
125272.....	32	401500.....	36	401925.....	72	411265.....	36
125273.....	32	401501.....	36	401926.....	72	411266.....	36
125274.....	32	401502.....	36	401927.....	73	411267.....	36
125275.....	32	401503.....	36	401928.....	73	411288.....	36
125276.....	32	401504.....	36	401929.....	73	411289.....	36
125277.....	32	401505.....	35	401930.....	73	411290.....	36
125330.....	55	401506.....	35	401931.....	73	411291.....	36
125331.....	56	401507.....	35	401932.....	36	411292.....	36
125333.....	56	401508.....	35	401933.....	36	411309.....	91
125334.....	56	401509.....	35	401934.....	36	411310.....	35
125335.....	56	401510.....	36	401935.....	36	411311.....	35
125336.....	56	401511.....	36	401936.....	36	411312.....	91
125337.....	56	401512.....	36	401989.....	85	411360.....	35
125338.....	56	401513.....	36	401993.....	91	411361.....	35
125339.....	56	401514.....	36	401994.....	91	411362.....	35
125340.....	56	401515.....	36	401995.....	91	411489.....	36
125341.....	56	401521.....	91	401996.....	91	411490.....	36
125342.....	56	401522.....	38	402239.....	91	411491.....	36
125343.....	56	401523.....	38	403948.....	41	411492.....	36
125344.....	56	401524.....	38	404006.....	80	411493.....	36

411631.....	70	413418.....	29	416733.....	48	421034.....	33
411632.....	70	413419.....	29	416736.....	48	421035.....	33
411633.....	71	413420.....	29	416737.....	48	421039.....	28
411637.....	72	413421.....	33	416738.....	49	421040.....	28
411638.....	72	413422.....	33	416739.....	49	421041.....	28
411639.....	72	413423.....	33	416740.....	49	421042.....	63
411640.....	72	413424.....	33	416741.....	49	421043.....	63
411641.....	71	413435.....	93	416742.....	49	421044.....	63
411642.....	71	413436.....	93	416743.....	49	421045.....	62
411657.....	41	413437.....	40	416744.....	49	421046.....	62
411658.....	41	413438.....	40	416745.....	49	421047.....	62
411659.....	41	413439.....	40	416746.....	49	421074.....	40
411800.....	29	413499.....	83	416747.....	49	421075.....	93
411801.....	33	413500.....	83	416748.....	49	421076.....	93
411802.....	29	413501.....	76	416749.....	49	421454.....	28
411803.....	33	413502.....	76	416750.....	49	421456.....	32
411804.....	29	413509.....	52	416751.....	49	421458.....	28
411805.....	33	413510.....	52	416774.....	43	421460.....	32
411806.....	29	413511.....	93	416775.....	43	421462.....	28
411807.....	33	413515.....	48	416776.....	44	421464.....	32
413141.....	91	413516.....	48	416777.....	93	421612.....	93
413142.....	93	413517.....	48	416778.....	43	421613.....	93
413143.....	93	413518.....	48	416779.....	43	421614.....	93
413144.....	93	413519.....	48	416780.....	93	421615.....	49
413145.....	93	413520.....	48	416781.....	44	421616.....	49
413146.....	93	413521.....	48	416782.....	43	421617.....	49
413147.....	93	413522.....	48	416783.....	43	421618.....	49
413148.....	93	413523.....	48	416795.....	43	421619.....	49
413149.....	93	416400.....	93	416796.....	43	421620.....	49
413150.....	93	416401.....	93	416797.....	93	421621.....	49
413151.....	93	416402.....	93	416798.....	43	421622.....	49
413152.....	93	416403.....	37	416799.....	43	421623.....	53
413153.....	27	416404.....	93	416900.....	93	421624.....	53
413154.....	27	416405.....	93	416901.....	43	421625.....	53
413155.....	27	416406.....	37	416902.....	43	421626.....	53
413156.....	27	416407.....	93	416903.....	43	421627.....	53
413157.....	27	416408.....	93	416908.....	47	421628.....	53
413159.....	27	416409.....	37	416915.....	67	421629.....	53
413160.....	27	416410.....	93	416916.....	67	421631.....	53
413161.....	27	416411.....	93	416917.....	61	421638.....	29
413162.....	27	416412.....	37	416918.....	61	421639.....	29
413163.....	27	416413.....	93	416919.....	60	421640.....	29
413164.....	27	416414.....	34	416920.....	60	421641.....	29
413165.....	27	416415.....	34	416921.....	61	421642.....	29
413166.....	27	416416.....	34	416922.....	29	421643.....	29
413167.....	27	416417.....	34	416923.....	60	421676.....	37
413168.....	27	416418.....	34	416924.....	60	421677.....	37
413169.....	93	416419.....	34	416925.....	61	421678.....	37
413170.....	38	416420.....	34	416926.....	61	421679.....	37
413171.....	38	416421.....	34	416927.....	60	421680.....	37
413172.....	93	416422.....	34	416928.....	60	421681.....	37
413173.....	93	416423.....	34	416929.....	60	421682.....	37
413174.....	38	416424.....	34	416930.....	60	421683.....	37
413175.....	38	416472.....	34	416932.....	64	421687.....	52
413176.....	93	416700.....	50	416934.....	65	421688.....	52
413177.....	30	416701.....	50	416938.....	66	421689.....	52
413178.....	30	416702.....	93	416940.....	64	421738.....	34
413179.....	30	416703.....	52	416942.....	65	421739.....	34
413180.....	30	416704.....	52	416946.....	66	421740.....	34
413181.....	30	416706.....	52	416948.....	64	421746.....	34
413182.....	30	416707.....	52	416950.....	65	421747.....	34
413183.....	30	416709.....	52	416954.....	66	421748.....	34
413184.....	30	416710.....	52	416975.....	46	421749.....	34
413185.....	30	416712.....	52	416976.....	46	421796.....	38
413186.....	30	416713.....	52	420695.....	37	421797.....	38
413187.....	30	416715.....	52	420696.....	37	421798.....	37
413188.....	30	416718.....	46	420697.....	37	421816.....	50
413312.....	35	416719.....	46	420793.....	77	421821.....	50
413313.....	35	416721.....	46	420794.....	77	421822.....	52
413314.....	35	416724.....	48	420795.....	77	421823.....	52
413334.....	35	416725.....	48	420796.....	77	421848.....	47
413335.....	93	416728.....	48	420797.....	76	421853.....	46
413336.....	93	416729.....	48	420798.....	83	421854.....	48
413417.....	29	416732.....	48	421033.....	33	421855.....	48

421856.....48
 421857.....48
 421858.....46
 421859.....48
 421860.....48
 421861.....48
 421862.....48
 421893.....37
 422159.....49
 422163.....81
 422363.....50
 422364.....50
 422365.....50
 422539.....50
 422540.....50
 422541.....50
 422545.....50
 422546.....50
 422551.....54
 422552.....54
 422553.....54
 422557.....54
 422558.....54
 422559.....54
 422560.....54
 422562.....54
 422564.....54
 422565.....54
 650408.....93
 650422.....93
 650828.....57
 650829.....57
 650830.....57
 650831.....57
 650832.....57
 650833.....57
 839001.....78
 839002.....78
 839003.....79
 839004.....79
 839005.....78
 839006.....78
 839010.....79
 839250.....88
 839251.....88
 839252.....88
 839253.....89
 839260.....93
 839261.....88
 839262.....88
 839263.....93
 839264.....89
 839265.....89
 839266.....89
 839267.....88
 839268.....89
 839280.....89
 TR105.....81

1. Weiss, A., et. Al., Statistical Brief #154, Readmissions to U.S. Hospitals by Procedure, 2010, Healthcare Cost and Utilization Project, Agency for Healthcare Research and Quality, April 2013. Data derived from readmission rates of colostomy, ileostomy, and other enterostomy procedures.
2. Hyde, L.Z., Al-Masrou, A.M., Kuritzkes, B.A., Suradkar, K., Valizadeh, N., Kiran, R.P. (2018). Readmissions after colorectal surgery: not all are equal. *International Journal of Colorectal Disease*, 33 (12), 1667-1674.
3. Taneja, C., Netsch, D., Rolstad, B.S., Inglese, G., Lamerato, L., Oster, G. (2017). Clinical and Economic Burden of Peristomal Skin Complications in Patients with Recent Ostomies. *Journal of Wound, Ostomy Continence Nursing*, 44 (4), 350-357.
4. Ostomy, Advocacy page. United Ostomy Association of America website. <https://www.ostomy.org/tag/survey/>. Accessed August 8, 2018.
5. Bulkeley, J.E., McMullen, C.K., Grant, M., Wendel, C., Hornbrook, M.C. Krouse, R.S. (2018). Ongoing Ostomy Self-Care Challenges of Long-Term Rectal Cancer Survivors. *Supportive Care in Cancer*, 26 (11), 3933-3939.
6. ConvaTec Internal in vitro lab testing. OCT o112.
7. Maria Teresa Szewczyk, MD, PhD; Grazyna Majewska, RN, ETN; Mary V. Cabral, MS, FNP-BC, CWOCN-AP; and Karin Hölzel-Piontek, RN; The Effects of Using a Moldable Skin Barrier on Peristomal Skin Condition in Persons with an Ostomy: Results of a Prospective, Observational, Multinational Study, *Ostomy Wound Management* 2014;60(12):16-26.
8. Bosio G., Pisani F., Lucibello, L. A Proposal for Classifying Peristomal Skin Disorders: Results of a Multicenter Observational Study. *OWM*. 2007; 53(9): 38-43.
9. Clinical Evaluation of Sur Fit Natura Cut to Fit Skin Wafers Gillian Seaton. Technical data on file, reference document WHR13940 TA238 2014.
10. Rolstad, BS. Principles and techniques in the use of convexity. *Ostomy Wound Manage*. 1996. 42(1); 24-32.
11. Pontieri, Lewis V. Basics of ostomy care. *Medsurg Nurs*. 2006; 15(4): 199-202.
12. Hanley, J. Convex stoma appliances: are we getting it right? *British Journal of Nursing (Stoma Supplement)*. 2013; 22(16).
13. Turnbull, GB. The convexity controversy. *Ostomy Wound Manage*. 2003; 49(1): 16-www.o-wm.com/article/1203. Accessed 29 Feb, 2016.
14. DRG data 2015, in T5 markets: Japan, UK, USA, Germany, China .
15. Boyd, K. Clinical Protocols for Stoma Care – Use of convex appliances. *Nursing Standard*. 2004; 18(20): 37-38.
16. Burch, 2011: Management of Stoma Complications. *NursingTimes.net*. 10 Nov 2011.
17. Data on file. Trio Protocol. PMS 002. December 2007. ConvaTec Inc.
18. In-Vitro: Trio Products 110099. Test Report AC-0170. July, 2012. Data on file. ConvaTec Inc.



The me+™ program is designed to give patients the support, insights and products they need to live life in a way that works for them.
It's free to sign up.


Learn more at mepluscare.com



Health Care Professionals

**TO ENROLL YOUR PATIENTS
in the me+ program**

call us at 877-585-0470 or
email concierge@convatec.com

 You can also scan this code to enroll your patient in our me+ program

Not a Health Care Professional?

CONTACT US TODAY to join the me+ program
call 1-800-422-8811 or email CIC@convatec.com



ConvaTec

Ostomy Care

For general product information visit www.convatec.com/ostomy

Disclaimer:

The reimbursement information provided by ConvaTec is intended to provide general information relevant to coding and reimbursement of ConvaTec's products only. Coverage and payment policies for the same insurer can vary from one region to another and may change from time to time because of ongoing changes in government and insurance industry rules and regulations. Therefore, please confirm HCPCS codes with your local DME-MAC, private insurer, or Medicaid agency before processing claims. ConvaTec does not guarantee coverage or payment of its products listed herein.