

AMFAB™

1446 South 35th Street
Galesburg, MI 49053


Phone 269.665.6703
800.432.6322
Fax 269.665.6717



Automatic Overbed Table/Model 1010H, 1010U

Long-Term Care Tables

Model 1010H (H-base)

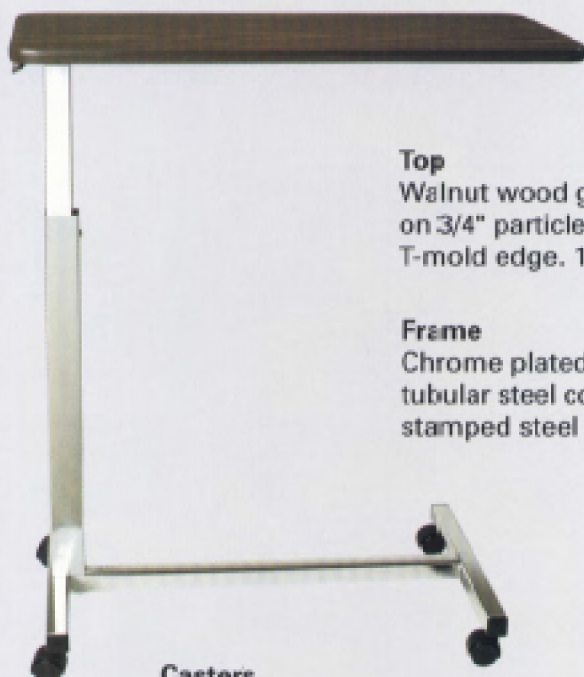


Casters
1 1/2" twin wheel

Top
Walnut wood grained melamine laminate on 3/4" particleboard core. Flush mounted T-mold edge. 15"W x 30"L.

Frame
Chrome plated 16 and 17 gage welded tubular steel column and base. 13 gage stamped steel support bracket.

Model 1010U (U-base)

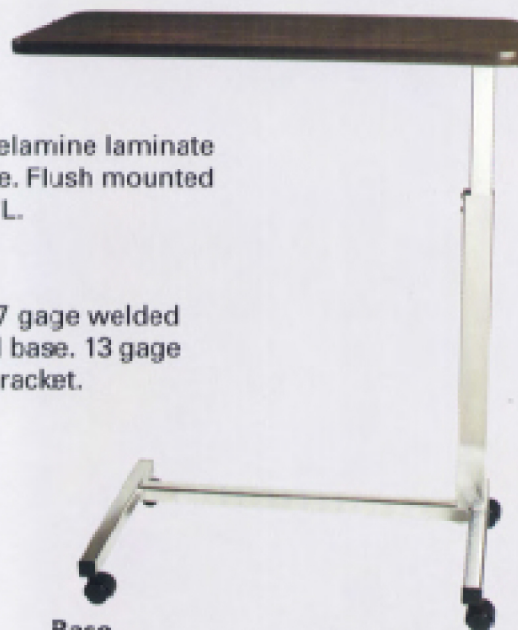

Base
U-base for easy use with chair.

Height
28" - 45" Infinitely adjustable.

AmFab's Automatic Overbed Tables Model 1010H & 1010U are economical, easy-to-use units for nursing home, hospital, and home convalescent needs.

Table top raises with the slightest upward pressure, lowers with a touch of the convenient trip lever, and locks automatically into position with AmFab's spring loaded mechanism.

The 1010U features a welded "U" base with the casters mounted at the extreme points on the legs for maximum stability. The 24" open area between the legs will accommodate most patient room seating to provide the seated resident with convenient access to the table.

A wide variety of high pressure laminate tops are available to complement your room decor and remain attractive throughout years of use. They are designed to resist damage by alcohol, boiling water, food acids, or ink spills.

Weight
27 lbs. Cube 1.1. UPS shippable.

Shipping
One table per carton, knocked down. Table easily assembled with four bolts.

- Options** (Available at additional cost)
- Alternate high pressure laminates*
 - Vanity & mirror (factory installed when ordered with table)
 - Painted frame*

* Please allow six weeks additional for delivery.

