



PHILIPS

DreamStation Go

Portable PAP therapy system

Innovation for people on the Go

Convenience
without compromise



The **DreamStation Go** is a compact and highly portable PAP system designed for patients who want to maintain their therapy at home and on the road.



Powerful, portable performance

At half the size of our System One and DreamStation platforms, DreamStation Go offers the same therapy and comfort features used by more than 5 million PAP users.



Mask compatibility

DreamStation Go gives your patients the freedom to use the mask they want.¹



Optional water-based heated humidification

The thinnest heated humidifier amongst leading PAP brands² and uses tap, bottled or distilled water.³



Optional battery

DreamStation Go's overnight battery can help your patients stay committed to their therapy when they are off the grid. Set at 10 cm of pressure, DreamStation Go displayed an average run time of 13 hours when operating on its overnight battery.




12 mm micro-flexible tubing

DreamStation Go features the smallest PAP tubing we've ever offered, providing enhanced flexibility, portability and connectivity. In a recent head-to-head trial, patients found the DreamStation 12 mm micro-flexible tubing to offer better flexibility, more comfort and greater ease of movement in comparison to standard tubing.⁴



Compliance reporting made easy

Like our other PAP devices, DreamStation Go is compatible with our patient data management system, which combines usage and compliance data from multiple Philips Respironics PAP devices into one simple report.

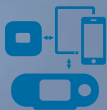

DreamMapper cross-device connectivity

Users of DreamMapper, our patient self-management tool, can sync data from multiple DreamMapper-enabled devices⁵ to one account. They also have the ability to track and monitor their combined progress across both home and portable PAP devices via one convenient dashboard on their smartphone, tablet or computer.

¹Masks with standard 22 mm outer diameter connection.

²When compared against ResMed AirSense, AirStart; Devilbiss IntelliPAP 1 and IntelliPAP 2; Fisher and Paykel – SleepStyle, Icon.

³In very hard water areas, it is recommended that you use distilled or bottled water.

⁴Head-to-head trial of DreamStation 15mm tubing vs. DreamStation 12 mm micro-flexible tubing, 2016.

⁵To see which therapy devices are compatible with DreamMapper, visit www.dreammapper.com/compatible.

Smart humidification, **all night long**

DreamStation Go's heated humidifier has been designed to use tap, bottled or distilled water¹ so travelers don't have the inconvenience of packing or purchasing distilled water when they arrive at their destination. It features water-saving technology designed to provide a full night of heated humidification² by actively monitoring ambient conditions and therapy settings and automatically adjusting moisture delivery when humidifier target time is enabled. Also, in the event your patient underfills or forgets to fill the humidifier, our heater plate shutdown feature will take over to prevent dry, uncomfortable air from being delivered to the patient for the remainder of the night.

Benefits

- Water-based heated humidifier is designed to use any potable, i.e. drinking water
- Water-saving technology designed to manage the device's water usage evenly throughout the night²
- Philips' first PAP heated humidifier to offer automatic heater plate shutdown
- The only heated humidifier with micro-flexible 12 mm tubing





¹In very hard water areas, it is recommended that you use distilled or bottled water.

²Adjusts humidification output when water-saver target time is enabled. Refer to the User Manual for more information.

Powerfully convenient features



Optional battery provides approximately 13 hours of therapy!¹



Intuitive color touchscreen offers quick, simple setup and easy navigation.



Built-in power supply reduces the number of parts managed by the patient.



Optional microSD card stores compliance and waveform data.²



DreamStation Go is the only travel PAP with a USB charging port,³ allowing for easy charging of devices.



¹PAP run-time test conditions: CPAP mode, pressure 10 cm H₂O, 12mm tubing, 37LPM of leak, 73.4 F room at an elevation of 1650 ft.

²Sold separately.

³When compared against HDM ZI, Somnetics Transcend and ResMed's AirMini devices, 2017.

⁴RTCA/DO-160G section 21, category M.

⁵In very hard water areas, it is recommended that you use distilled or bottled water.



Bluetooth is standard with every device.



Unit complies with FAA requirements for in-flight use.⁴



Device is compatible with DreamMapper, our patient self-management tool.

12 mm tubing

Micro-flexible tubing is our smallest and lightest ever.



Heated humidification

Water-based heated humidifier that uses tap, bottled or distilled water.⁵

Create a world of opportunities



Heated humidifier

Water-based heated humidifier that conveniently uses tap, bottled or distilled water.¹



Overnight battery

Compact and lightweight, the battery features a recharge time of just 4-5 hours.



12 mm micro-flexible tubing

This tubing offers more comfort and greater ease of movement than that of standard tubing.²



Small travel kit

Kit includes an easy-to-pack protective case for the DreamStation Go and a washable bag for the mask and tube.



Medium travel kit

Kit includes a protective case that accommodates the DreamStation Go with either the battery or heated humidifier as well as a washable bag for the mask and tube.



Air filters

Choice of a reusable pollen or disposable 30-day fine filter.

¹In very hard water areas, it is recommended that you use distilled or bottled water.

²Head-to-head trial of DreamStation 15mm tubing vs. DreamStation 12mm micro-flexible tubing, 2016.



DreamStation Go Auto specifications

Pressure range	4 to 20 cm H ₂ O
Pressure mode(s)	CPAP, APAP
Flex pressure relief	1 to 3
Flex modes	C-Flex, C-Flex+, A-Flex
Ramp time	0 to 45 min (5-min increments)
Starting ramp pressure	Patient adjustable, SmartRamp
Dimensions	15 cm L x 15 cm W x 5.8 cm H
Weight with power supply	844 g (1.86 lb)
Filters	Reusable pollen, disposable fine
Device setup	Touchscreen color display
Data storage capacity	On-board: > 1 year MicroSD card (optional): lifetime
Connectivity	Bluetooth
Altitude compensation	Automatic
Electrical requirements	100–240 VAC, 50/60 Hz
Warranty	Two years

DreamStation Go CPAP specifications

Pressure range	4 to 20 cm H ₂ O
Pressure mode(s)	CPAP
Flex pressure relief	1 to 3
Flex modes	C-Flex, C-Flex+
Ramp time	0 to 45 min (5-min increments)
Starting ramp pressure	Patient adjustable, SmartRamp
Dimensions	15 cm L x 15 cm W x 5.8 cm H
Weight with power supply	844 g (1.86 lb)
Filters	Reusable pollen, combination fine
Device setup	Touchscreen color display
Data storage capacity	On-board: > 1 year MicroSD card (optional): lifetime
Compliance meter	Breathing detection
Altitude compensation	Automatic
Electrical requirements	100–240 VAC, 50/60 Hz
Warranty	Two years

Ordering information

Flow generators	Part number
DreamStation Go CPAP	DSG400S11
DreamStation Go Auto CPAP	DSG500S11
Flow generators with heated humidifier	Part number
DreamStation Go CPAP with heated humidifier	DSG400H11
DreamStation Go Auto CPAP with heated humidifier	DSG500H11
Humidifier	Part number
DreamStation Go heated humidifier	DSGH
DreamStation Go humidifier replacement tank bottom	1140264
Power	Part number
DreamStation Go overnight battery	1133281
Travel kits	Part number
DreamStation Go small travel kit ¹	1133275
DreamStation Go medium travel kit ²	1133276
Consumables/Replacement	Part number
DreamStation micro-flexible 12 mm tubing	PR12
Power cord 6 ft white	1132551
Power cord 10 ft white	1132552
Beauty panel - front	1133741
Beauty panel - rear	1133742
Reusable pollen filter 2-pk	1133743
Disposable 30-day fine filter 6-pk	1133744

Heated humidifier specifications

Dimensions	14.7 cm L x 15.1 cm W x 5.9 cm H
Weight of humidifier	344.7 g (0.76 lb)
Warranty	Two years
Tubing compatibility	Micro-flexible 12 mm, 15 mm and 22 mm standard tubing
Tank water capacity	240 ml

*All other specifications are the same as DreamStation Go Auto CPAP

¹Small travel kit includes an easy-to-pack protective case for the DreamStation Go and a separate washable mask and tube bag.

²Medium travel kit includes a protective case that accommodates the DreamStation Go with either the overnight battery or heated humidifier. Also includes a separate washable mask and tube bag.

A close-up photograph of a white Philips DreamStation Go CPAP machine. A white, corrugated, flexible tube is connected to the top of the machine and arches over the top right. The machine has a clear plastic water chamber at the bottom with a 'MAX' fill line. The background is a warm, out-of-focus brown.

Learn more

For more information about DreamStation Go,
visit www.philips.com/dreamstationgo.



© 2019 Koninklijke Philips N.V. All rights reserved.
Specifications are subject to change without notice.

www.philips.com/respironics

Caution: U.S. federal law restricts these devices to
sale by or on the order of a physician.

RRDPGH EB 3/27/19 PN 1141233 MCI 4108863
1010 Murry Ridge Lane, Murrysville, PA 15668
800 345 6443 • 724 387 4000