

# Kangaroo OMNI™ Enteral Feeding Pump Overview

## Kangaroo OMNI™ Feeding Set Overview

The below is an example of a Kangaroo OMNI™ ENTelliSet™ Feeding Set with Flush Bag for illustration purposes only.

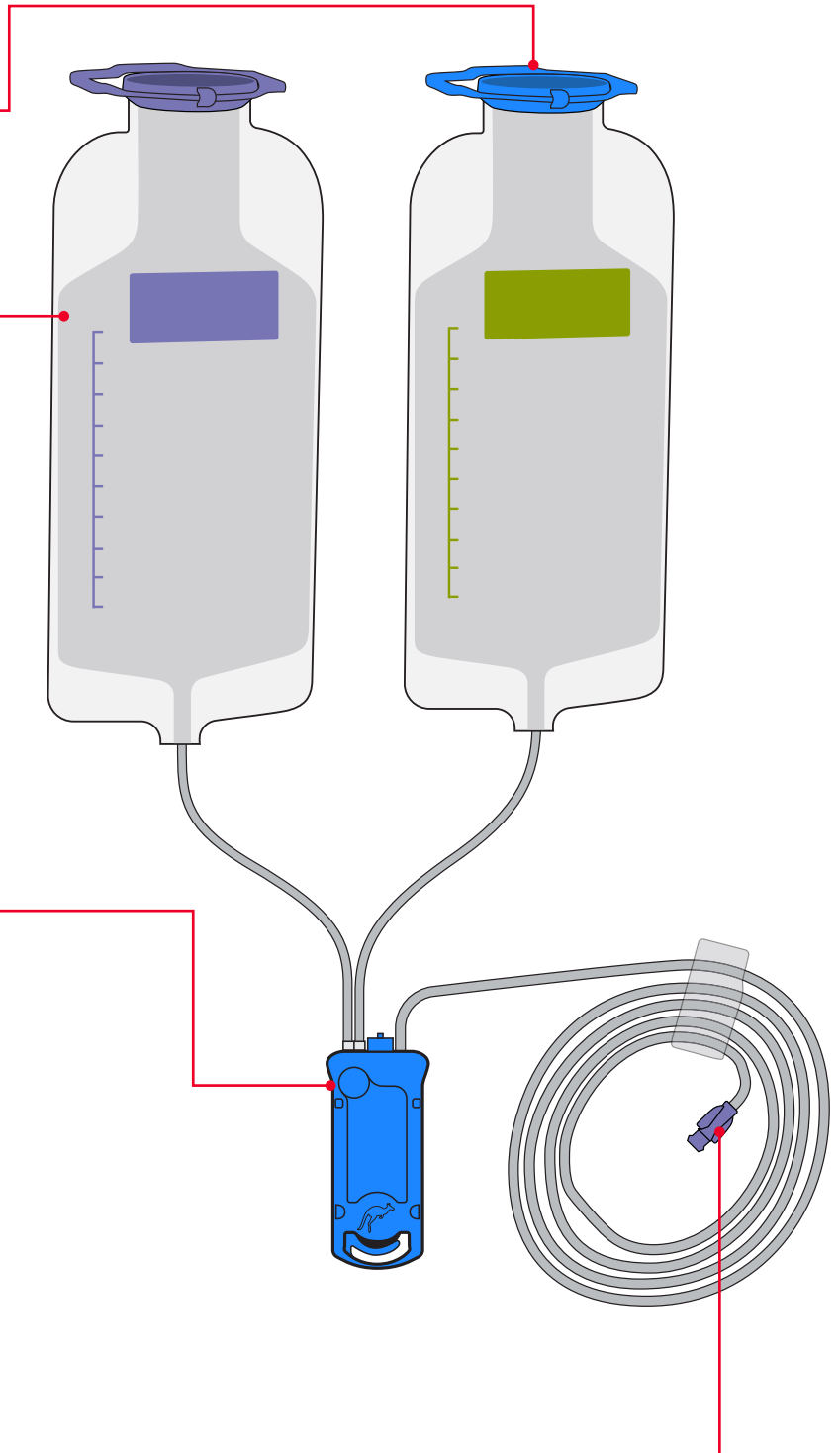
**Note:** product availability and type vary by region.

**Flush Bag:** (when applicable) intended for use in the administration of water to help manage patient hydration requirements and help keep the feeding tube patent to avoid tube clogs.

**Feeding Set Bag:** (also available in Spike Set option) intended for use in the administration of enteral fluids.

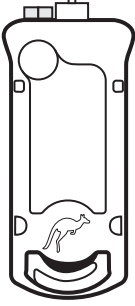
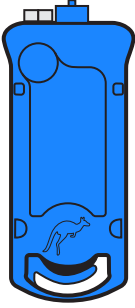

**Kangaroo OMNI™ Feeding Sets:** various feeding set options are compatible and available for use with the Kangaroo OMNI™ Enteral Feeding Pump Platform, see table on following page for the list of Kangaroo OMNI™ Feeding Sets.

**ENFit™ Feeding Set Connector:** connects the feeding set to the feeding tube. The ENFit™ enteral feeding connector features a unique design to reduce the risk of enteral feeding tube misconnections and improve patient safety using internationally recognized ISO 80369-3 guidelines.



# Kangaroo OMNI™ Enteral Feeding Pump Overview

## Kangaroo OMNI™ Feeding Set Overview

CASSETTE	COLOR	PRODUCT NAME	DESCRIPTION
	Clear	Kangaroo OMNI™ Feeding Set	Feeding set designed to deliver available formula (categorized as level 0-1 within the IDDSI framework)
	Blue	Kangaroo OMNI™ Feeding Set with Flush Bag	Feeding set designed to deliver available formula (categorized as level 0-1 within the IDDSI framework) and water for hydration/flushing
	Green	Kangaroo OMNI™ ENTelliSet™ Thick Formula Feeding Set Kangaroo OMNI™ ENTelliSet™ Thick Formula Feeding Set with Flush Bag	Feeding set designed to deliver homogenized, blenderized formulas (categorized as level 2, 3, or 4 drinks within the IDDSI framework) Feeding set designed to deliver homogenized, blenderized formulas (categorized as level 2, 3, or 4 drinks within the IDDSI framework) and water for hydration/flushing