



INVACARE® ETUDE HC HOMECARE BED





Etude HC Homecare Bed

The ETUDE HC Homecare bed, with its attractive aesthetics, flexible design, simple operation and easy handling is the ideal choice for community home care.

The ETUDE HC is the only advanced homecare bed to offer Shear Management Technology as a standard feature, reducing shear and friction providing superior comfort and safety. The SMT feature is provided through the auto-regression of the head section while repositioning in ETUDE HC through use of the 8 function hand control.

Features and Options

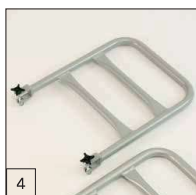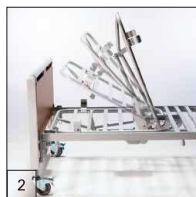

1. Built in Dual Height Brackets.
2. Shear Management Technology.
3. Transport and Storage Kit.
4. Support Handle.
5. Lifting Support Trapeze Kit.
6. 4\"/>

Modern and Flexible

ETUDE HC delivers a modern, streamlined and discreet design perfectly complimenting any interior. The ETUDE HC can be configured with attractive wood side panels and a choice of safety rails to accommodate specific transfer, safety and support requirements. ETUDE HC provides a deck height ranging from 7.8\"/>

Easy Handling & Storage

ETUDE HC recognizes that portability, storage and set-up are important features for users and caregivers alike. The simplistic design allows for simple assembly without tools and transport can be managed safely and comfortably as each section of ETUDE HC weighs less than 52 lbs. Smooth roll 3\"/>

Highest Safety Levels

The ETUDE HC conforms to the IEC bed standard IEC60601-2-52.

This equates to the highest degree of safety for both the client and the caregiver. The 8 function pendant enables the caregiver to lock out any of the beds functions for enhanced patient safety. ETUDE HC's selection of side rail supports have all met the demands for strength and entrapment guidelines required in the new safety standards for beds. Visit www.invacare.ca for more information.

Invacare Canada LP
www.invacare.ca



Width outside:
40"
Width inside:
36"



Length outside:
86.5"
Length inside:
82"



Height adjustable:
7.8" - 26.5"



Mattress support
dimensions:
36" - 80"



0 - 70°
Backrest angle



0 - 18°
Thigh angle



0 - 25°
Legrest lift



11°
Tilt option*



Product weight:
166 lbs.



Max. Bed weight capacity: 440 lbs.
Max. User weight capacity: 360 lbs.

Frame colour and wood work



Aluminium Grey



Beech

Warranty

5 years on frame, 3 years on actuators, 2 years on pendant and electronics.

Warranty against manufacturing defects and does not include normal wear, misuse or abuse.

Custom modifications to ETUDE HC void manufacturer Warranty.

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada (800)
668-5324

©2019
Invacare Corporation. All
rights reserved.
Trademarks are identified
by the symbols TM and [®].
All trademarks are owned
by or licensed to Invacare
Corporation unless
otherwise noted.
Specifications are subject
to change without notice.

Form No. 15-912C
Rev: 4/19

