

Live the Longer
Life You Love[®]



204I BEDCANE

A Bed Handle With Free Organizer Pouch



HANDY 4-POCKET ORGANIZER





EASILY SIT UP IN BED



EASILY STAND UP



ORGANIZER KEEPS ITEMS CLOSE BY

Features & Benefits

- **Ergonomic No-Slip Handle** - Allows for easy transfer in and out of bed.
- **4-Pocket Organizer Included** - Prevents entrapment and provides storage space for handy items.
- **Portable** - Folds for travel or storage.
- **Safety Strap** - Secures rail to any home or hospital bed.

Specifications

Rail Height: Adjustable 19" - 22"

Rail Width: 15"

Base Dimensions: 23" x 19"

Weight of product: 12 lbs

Weight Capacity: 350 lbs

Made for Bed Type: Attaches to any home or hospital bed with included Safety Strap

Assembly: Installs in minutes with 3 bolts and safety pin

Package Dimensions: 25" x 19" x 2"

Case of 4 Dimensions: 25" x 19" x 8"; 49 lbs

How Do We Compare

	Bedcane	Other Bed Rails
Organizer Pouch Included	✓	✗
No-Slip Grip	✓	✗
Sturdy Wood Base with Safety Strap	✓	✗
Folds for Travel	✓	✗
Provides Assistance In & Out of Bed	✓	✓
Durable Powder Coat Finish	✓	✓


Easy Assembly



Attach Handle to Bracket with Bolt and Safety Lock Pin



Attach Pocket Organizer and Safety Strap

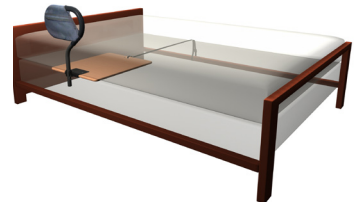

Secure to Bed Frame with Included Safety Strap