

**PHILIPS**  
RESPIRONICS

InnoSpire Mini

So your patient can  
**enjoy life**

The InnoSpire Mini compressor nebulizer system allows patients to effectively nebulize most commonly prescribed aerosol medications whenever they want and wherever they may be.

InnoSpire Mini is part of our InnoSpire family of compressor nebulizer systems that provides reliable drug delivery to patients.

An ideal traveling companion, InnoSpire Mini's compact size allows for easy transport while a breath-enhanced reusable nebulizer makes it easy to use.

**Proven nebulizer technology**

Coupled with the SideStream Plus reusable handset, the InnoSpire Mini provides highly efficient aerosol therapy in fast treatment times.

**Key features**

- Matched with SideStream Plus
- Compact and lightweight
- Battery powered (optional)
- Up to 90 min autonomy when battery powered
- Multiple charging options





**Reliable aerosol delivery**

SideStream Plus has been designed to draw in more air, via an active venturi system, resulting in fast drug delivery and short treatment times.

SideStream Plus' unique diamond jet design helps minimize wear from frequent use and aids consistent delivery of a high quality aerosol.


# Giving your patient **choice**

The InnoSpire Mini can be powered by the traditional wall outlet or by an optional rechargeable lithium ion battery, which can operate the compressor for up to 90 minutes between charges. The InnoSpire Mini can also be powered by a 12V DC car adapter while on the move; its portability enables patients to take their treatment virtually anytime and anywhere, at their convenience.

The versatile InnoSpire Mini is also compliant with selected FAA (US Federal Aviation Administration) environmental standards\*\*\* for airborne equipment.

### SideStream Plus

Equipped with a low resistance valve to deliver fast treatments with less waste.

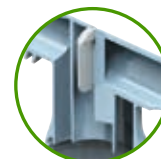

### Valve opens

Significantly less pressure is needed to open the inspiratory valve, making it ideal for weaker breathers.



### Filter

Our optional filter accessory also helps to reduce the unnecessary release of medication in ambient air, especially when nebulizing antibiotics.



### Ordering information

### Part no.

|   |         |
|---|---------|
| InnoSpire Mini without battery – EU plug** (with user guide†) | 1109291 |
| InnoSpire Mini with battery – EU plug** (with user guide†)    | 1109292 |
| InnoSpire Mini with battery – EU plug** (with user guide††)   | 1109438 |
| InnoSpire Mini with battery – EU plug** (with user guide††)   | 1109294 |
| InnoSpire Mini with battery – UK plug** (with user guide††)   | 1111453 |

### Accessories/replacement parts

### Part no.

|  |         |
|--|---------|
| InnoSpire Mini replacement air filter pack   | 1109488 |
| SideStream Plus with mouthpiece reusable kit | 1092001 |

### Performance data

|                               |          |
|-------------------------------|----------|
| MMAD                          | 4.77µm   |
| Fine particle fraction        | 51.44%   |
| Treatment time (2.5mL charge) | 6.88 min |

### Technical specifications

|                          |                                   |
|--------------------------|-----------------------------------|
| Compressor max. pressure | 207kPa (30psi)                    |
| Average flow rate        | 3.4LPM @ 12.7psi                  |
| Weight                   | Compressor: 420g<br>Battery: 172g |
| Sound level              | 50 ± 4 dBA                        |
| Operating temperature    | 5°C to 40°C                       |
| Operating humidity       | 10% RH to 95% RH                  |

\* Can be fully recharged in approximately 6 hours

\*\* Including rechargeable battery (except for PN 1109291), 12V DC car adapter, AC adapter, unit carrying case and (4) spare filters

\*\*\* DO-160G standards for airborne equipment and the international basic safety and essential performance standard IEC 60601-1 third edition

† Languages: EN, BG, CZ, EL, ES, HR, ID, IT, HU, PT, RO, SK, SL, SR, TR, and KO

†† Languages: EN, AR, DA, DE, EE, FI, FR, IL, KZ, LT, LV, NL, NO, PL, RU, SV, and TH

