


# Comprehensive Solutions For Advanced Wound Care





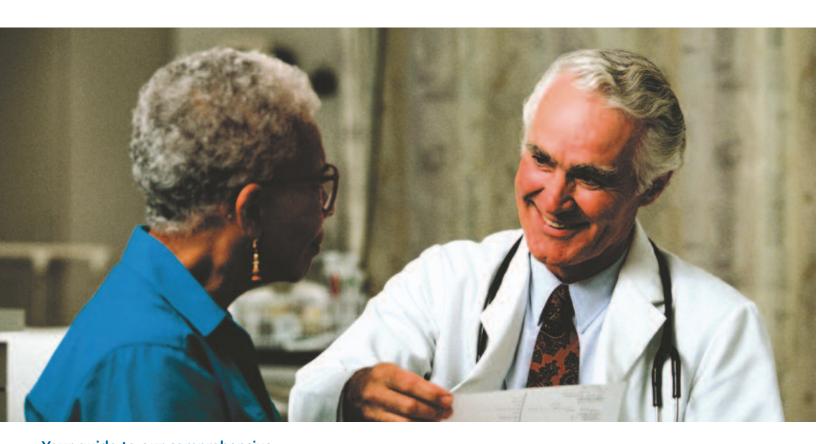

Your guide to our comprehensive solutions for advanced wound care, including a handy tearsheet of at-a-glance care plan guidelines

# Table of Contents

Innova	tive Services2
Hydro	Foam Back Dressing Film Back Dressing Sacral Dressing Thin Dressing
Island	and Barrier Island Dressings6
Unna I	Boots8
Transp	earent Film Dressings9
Foam I	Dressings
Hydro	gel Dressings
Calciur	m Alginate Dressings 14 Sheet Dressing Rope Dressing
Wound	d Cleansers
Derma	Gran®

**Care Plan Guidelines Tearsheet** 

1

# Benefit from innovative services from an industry leader

There's a reason more than 53,000 healthcare providers count on McKesson Medical-Surgical. We are the nation's leading distributor of medical-surgical supplies in the long-term care and physician office markets.

We are the industry frontrunner in developing quality products as well as innovative point-of-care and customer service programs. Every day, in every way, we empower customers like you to save time, save money and give your residents superior care that sets you apart.

- Quality products. Depend on the industry knowledge and product quality that have been part of the fabric of our company for more than 170 years.
- Consistent service. Expect us to work one-on-one with you so that we can understand which of our services will benefit you most. Also expect us to have the products you need when you need them.

  Satisfaction guaranteed.
- Customer savings programs.
  Save on the supplies you use most by having your purchases applied toward our savings programs.
- Inventory management programs. Have us design a custom program to help you maximize efficiencies, reduce costs and streamline inventories.

# Maximize savings with value-priced performance products

At McKesson Medical-Surgical, the company you've come to know for product innovation, we proudly introduce our line of advanced wound care products.

McKesson Brand offers you a comprehensive range of products competitively priced to deliver real value.

Please take a moment to look through our catalog. You will notice the "Latex Free" symbol by all latex-free products. If you can't find the product you need, keep in mind that the McKesson Brand product line is growing every day. Contact your McKesson Medical-Surgical Account Manager with questions, or call your local Customer Service Representative.





# Advanced Wound Care

For the customer service center nearest you, call **800.446.3008** 

www.mckesson.com



# Hydrocolloid Dressings

Hydrocolloid dressings are designed for the management of partial and full thickness wounds such as leg ulcers and pressure ulcers. Hydrocolloid dressings help maintain a moist wound-healing environment, protect the wound from external contamination and promote autolytic debridement.

Hydrocolloid dressings should not be used on patients with third-degree burns or known sensitivities to the components of the dressings.

Hydrocolloids can be used to help protect skin at risk for breakdown from damage due to friction.

Hydrocolloid dressings are designed to manage light to moderate wound exudate and provide a moist woundhealing environment for the wound bed. Conditions that lead to the formation of the wound (pressure, incontinence, friction and shearing forces) must be considered and addressed as part of the total wound care treatment plan.



#### Foam Back Dressing

Since patient comfort is key, hydrocolloids with foam backings are possibly the most popular type available, as the foam provides an added layer of comfort. The edge of our version is tapered to prevent bunching and rolling up, therefore maximizing wear time and increasing cost-effectiveness.

#### EXTENDED CARE

85142100	4"x 4"	5/Box
51662100	6"x 6"	5/Box
85182100	8"x 8"	5/Box

#### PRIMARY CARE

61-85144	4"x 4"	5/Box I 20 Boxes/Case
61-85166	6"x 6"	5/Box I 16 Boxes/Case
61-85188	8"x 8"	5/Box I 16 Boxes/Case



## Film Back Dressing

This dressing represents the latest in hydrocolloid technology. The film backing is designed to act as a barrier to external contaminants, and it is translucent, which enables clinicians to monitor the wound between dressing changes. The edge is also tapered, which reduces the risk of bunching and rolling up, maximizing wear time and increasing cost-effectiveness.

#### EXTENDED CARE

52222100	2"x 2"	20/Box	
85242100	4"x 4"	10/Box	
85262100	6"x 6"	5/Box	
85282100	8"x 8"	5/Box	

#### PRIMARY CARE

61-85222	2"x 2"	20/Box I 20 Boxes/Cas
61-85244	4"x 4"	10/Box I 20 Boxes/Cas
61-85266	6"x 6"	5/Box I 16 Boxes/Cas
61-85288	8"x 8"	5/Box I 16 Boxes/Cas



#### **Sacral Dressing**

Sacral wounds are some of the most difficult wounds to manage, as anatomical curves in combination with varying levels of exudate create great challenges. This sacral dressing is designed to meet these challenges, and is made with the same design components of our film back hydrocolloid dressing.

#### EXTENDED CARE

85192100	Sacrai	2/ROX

#### PRIMARY CARE

61-85192	Sacral	5/Roy I	16	Royac/Cas



## **Thin Dressing**

This very thin dressing enhances a clinician's ability to monitor the wound without having to remove the dressing. It is very conformable and can be used to help protect against damage due to friction.

#### XTENDED CARE

53222100	2"x 2"	20/Box
85342100	4"x 4"	10/Box
85362100	6"x 6"	5/Box

61-85322	ZXZ	20/Box I 20	Roxes/Case
61-85344	4"x 4"	10/Box I 20	Boxes/Case
61-85366	6"x 6"	5/Box I 16	Boxes/Case



# Island and Barrier Island Dressings

Island and barrier island dressings may be used as a primary or secondary dressing. The island dressing is a bordered gauze dressing consisting of an absorbent, non-adherent pad covered by a breathable adhesive border. The barrier island dressing is a composite dressing appropriate for wounds with light to moderate exudate when a semi-occlusive environment is desired to protect from contamination.



#### **Island Dressing**

When a gauze dressing with adhesive border is required, the island dressing offers the convenience of an all-in-one dressing. An absorbent non-adherent central pad is surrounded by a breathable adhesive border, then covered with a layer of soft, highly conformable gauze. The non-adherent pad in the center is designed to avoid disrupting the wound bed upon removal.

#### EXTENDED CARE

89022000	2"x 2"	25/Box
90442000	4"x 4"	25/Box
89062000	6"x 6"	25/Box
89042000	4"x 6"	25/Box
89482000	4"x 8"	25/Box
89682000	6"x 8"	25/Box
89412000	4"x 10"	25/Box
89422000	4"x 14"	25/Box

#### PRIMARY CARE

61-89022	2"x 2"	25/Box I 4 Boxes/Case
61-89044	4"x 4"	25/Box I 4 Boxes/Case
61-89066	6"x 6"	25/Box I 4 Boxes/Case
61-89046	4"x 6"	25/Box I 4 Boxes/Case
61-89048	4"x 8"	25/Box I 4 Boxes/Case
61-89068	6"x 8"	25/Box I 4 Boxes/Case
61-89041	4"x 10"	25/Box I 4 Boxes/Case
61-89042	4"x 14"	25/Box I 4 Boxes/Case



#### **Barrier Island Dressing**

Our barrier island dressing offers the convenience of an all-in-one dressing by combining an absorbent non-adherent central pad, with film backing. This film is designed to provide an additional barrier to external contaminants, making this dressing appropriate for patients suffering from incontinence and those who are able to bathe or shower.

## EXTENDED CARE

89142000	4"x 4"	25/Box
89162000	6"x 6"	25/Box
89182000	6"x 8"	25/Box

61-89144	4"x 4"	25/Box I 4 Boxes/Cas
61-89166	6"x 6"	25/Box I 4 Boxes/Cas
61-89168	6"x 8"	25/Box I 4 Boxes/Cas



# Unna Boots

An unna boot is a gauze bandage that has been evenly impregnated with a non-hardening paste of zinc oxide. It is traditionally used with physician supervision in the management of venous stasis leg ulcers, lymphatic edema, thrombophlebitis and skin grafts on burn wounds. It should be applied by a trained clinician. Available with or without calamine.

## **Unna Boots**

EXTENDED CARE			PRIMARY CARE		
		Each	65-30	3"x 10 yds.	12/0200
	3"x 10 yds.			,	
65042100	4"x 10 yds.	Each	65-40	4"x 10 yds.	
65302100*	3"x 10 yds.	Each	65-30C*	3"x 10 yds.	
65402100*	4"x 10 yds.	Each	65-40C*	4"x 10 yds.	12/Case

\* indicates with calamine



# Medi-Pak™ Performance

# Transparent Film Dressings

Transparent, adhesive film is moisture vapor permeable and conformable. These products may be used to provide a moist wound-healing environment for superficial wounds; also ideal for use as a secondary dressing or on IV sites to provide catheter fixation.

#### **Films**

#### EXTENDED CARE

31232100	2 3/8" x 2 3/4"	Tab Delivery	100/Bo
		,	
30232100	2 3/8"x 2 3/4"	Frame Delivery	100/Bo
30452100	4"x 4 3/4"	Frame Delivery	50/Box
30682100	6"x 8"	Frame Delivery	10/Box

51-83123	2 3/8"x 2 3/4"	Tab Delivery	100/Box I 10 Boxes/Case
51-83023	2 3/8"x 2 3/4"	Frame Delivery	100/Box I 4 Boxes/Case
51-83045	4"x 4 3/4"	Frame Delivery	50/Box I 4 Boxes/Case
51-83068	6"x 8"	Frame Delivery	10/Box   8 Boxes/Case







# Foam Dressings

Sterile foam dressings are made from polyurethane and are available with or without a protective film covering that helps shield the wound from outside contaminants. These dressings are designed to manage moderate to highly exudating, partial or full-thickness wounds. The polyurethane foam sheets are hydrophilic and resist breakdown, making dressing changes easier for both the healthcare professional and patient.



# Hydrocellular Foam Dressing with Film Backing

The non-adhesive film-backed foam sheet dressing is designed to manage moderate to heavy levels of exudate. Since there is no adhesive, a secondary dressing is required to secure it in place. The film backing allows for the dressing to be wiped down and cleaned easily when held in place by tape.

#### EXTENDED CARE

84542100	4"x 4"	10/Box
84562100	6"x 6"	10/Box

#### PRIMARY CARE

61-84544 4" x 4" 10/Box I 4 Boxes/Case 61-84566 6" x 6" 10/Box I 4 Boxes/Case



## **Hydrophilic Foam Dressing**

This non-adhesive, non-film-backed foam sheet dressing is more appropriate when a retention bandage is used as the secondary dressing holding the foam in place.

#### EXTENDED CARE

84322100	2"x 2"	10/Box
84342100	4"x 4"	10/Box
84312100	1"x 8"	10/Box

#### PRIMARY CARE

61-84322	2"x 2"	10/Box I 4 Boxes/Case
61-84344	4"x 4"	10/Box I 4 Boxes/Case
61-84318	1"x 8"	10/Box I 4 Boxes/Case



## Thin Adhesive Foam

The thin adhesive foam dressing, similar to the thin hydrocolloid dressing, is designed to provide protection to at-risk skin. Since it is so thin, it is designed to handle only low to moderate amounts of exudate, while providing a barrier to external contamination. The thin foam dressing is adhesive-based, so it does not require a secondary dressing to keep it in place.

#### EXTENDED CARE

84422100	2"X 3"	10/RO
84442100	4"x 4"	10/Bo

#### PRIMARY CARE

61-84423	2"x 3"	10/Box I 4 Boxes/Case
61-84444	4"× 4"	10/Rox L4 Roxes/Case



#### **Adhesive Foam**

The adhesive foam dressing with adhesive border has a thin foam backing designed to provide a barrier to external contamination with the addition of a thicker non-adhesive foam island in the center. This additional foam is designed to allow the dressing to handle greater amounts of exudate, while keeping the wound contact layer adhesive-free. Thus, disruption to the wound bed is kept at a minimum.

#### EXTENDED CARE

84242100	4"x 4"	10/Box
84262100	6"x 6"	10/Box

61-84244	4"x 4"	10/Box I 4 Boxes/Cas
61-84266	6"x 6"	10/Box I 4 Boxes/Cas



# Hydrogel Dressings

Hydrogels are useful for rehydrating sloughy or necrotic tissue and enhancing autolytic debridement. They are best used on low exuding or dehydrated wounds such as minor burns, abrasions/lacerations, donor sites and pressure ulcers.

Hydrogel sterile dressings are for the wound management of partial and full thickness pressure ulcers, stasis ulcers, first and second degree burns, cuts and abrasions.



#### Sheet

Hydrogel sheets are sterile and designed to act as a primary cover for superficial or shallow wounds. A secondary dressing is required to keep these dressings in place. A removable sheet allows this dressing to be either occlusive or not occlusive.

EXT	ΓEΝ	۱D	ED	CA	<b>ARE</b>

82022100	2"x 2"	10/Box
82032100	2"x 3"	10/Box
84022100	4"x 4"	10/Box
82062100	6"x 8"	10/Box



#### PRIMARY CARE

61-82022	2"x 2"	10/Box I 4 Boxes/Case
61-82023	2"x 3"	10/Box I 4 Boxes/Case
61-82044	4"x 4"	10/Box I 4 Boxes/Case
61-82068	6"x 8"	10/Box I 4 Boxes/Case



## **Impregnated Gauze**

The flexibility of the sterile hydrogel impregnated gauze dressings allows them to be used as either primary cover or filler dressings. A secondary dressing is required to secure the dressing in place.

#### EXTENDED CARE

53022100	2"x 2"	10/Box
53042100	4"x 4"	10/Box

#### PRIMARY CARE

61-53022	2"x 2"	10/Box I 4	Boxes/Case
61-53044	4"x 4"	10/Box I 4	Boxes/Case



#### **Amorphous**

The 1 oz. amorphous gel is packaged in an innovative Bellows (accordion-like) tube, designed to ease dispensing of the gel into difficult-to-reach areas. The 3 oz. size is available in a multipleuse standard tube container.

#### EXTENDED CARE

81112100	1	Fl.	Oz.	Each
81132100	3	Fl.	Ozs.	Each

61-81110	1	Fl.	Oz.	12/Cas
61-81130	3	Fl.	Ozs.	12/Cas



# Calcium Alginate Dressings

Calcium alginate dressings are primary dressings to be used on partial or full thickness wounds with moderate to heavy exudate. Highly absorbent, these dressings easily conform to the wound cavity or tunnel. Additionally, calcium alginates can be used to promote hemostasis when minor bleeding is present.



#### **Sheet Dressing**

The square and rectangular sheets can be used as either primary cover dressings or to fill full thickness wounds.

#### EXTENDED CARE

88022100	2"x 2"	10/Box
88042100	4"x 4"	10/Box
88482100	4"x 8"	10/Box

#### PRIMARY CARE

61-88022	2"x 2"	10/Box I 4 Boxes/Case
61-88044	4"x 4"	10/Box I 4 Boxes/Case
61-88048	4"x 8"	10/Box I 4 Boxes/Case



## **Rope Dressing**

Rope calcium alginate is used for packing deep or undermined wounds.

#### EXTENDED CARE

88122100	3/4"x 12"	10/Box
88102100	3/4"x 36"	Bottle

#### PRIMARY CARE

61-88112	3/4"x 12"	10/Box I 4 Boxes/Case
61-88100	3/4"x 36"	12 Bottles/Case

## Medi-Pak™ Performance

# Wound Cleansers

Our wound cleansers are specially formulated for the cleansing of pressure ulcers, stasis ulcers, diabetic ulcers, post-surgical wounds, first and second degree burns, cuts, abrasions and minor irritations of the skin. They effectively cleanse a wound before dressing application, and contain the surfactant Poloxamer 188, known to be gentle on tissue.



#### **Preservative-Free Cleanser**

This sterile (preservative-free) wound cleanser is for single use on chronic and acute wounds.

#### EXTENDED CARE

91012100 118 ml / 4 Fl. Ozs. Each

PRIMARY CAP

61-69101 118 ml / 4 Fl. Ozs. 12/Case



#### **Preserved Cleanser**


This non-sterile (preserved) wound cleanser is for multiple use on chronic and acute wounds. The nozzle can be set to either spray or stream to accommodate the desired pressure for dislodging materials from the wound bed.

#### EXTENDED CARE

92012100 473 ml / 16 Fl. Ozs. Each

PRIMARY CARE

61-69201 473 ml / 16 Fl. Ozs. 12/Case

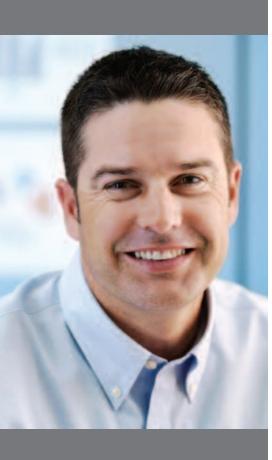

#### McKesson

# DermaGran®

The family of McKesson skin protectant products features zinc-nutrient balanced pH technology, which helps restore the acid mantle and optimizes the wound environment. This can mean the difference between successful healing and chronic wounds.

# **DermaGran® GP General Skin Protectant** (not pictured)

DermaGran® GP is recommended for use as a moisture barrier on external skin areas where repeated exposure to body excrements and exudates may cause skin breakdown. Contains allantoin and aloe vera gel. May be used as a skin barrier on friction points.



#### EXTENDED CARE

20051400 DermaGran® GP 4 Oz. Eac

PRIMARY CARE

61-S2054Z Dermagran® GP 4 Oz. 12/Case



#### **Hydrophilic Wound Dressing**

The DermaGran®-B Wound Dressing provides filler for the wound. It contains a patented zinc-nutrient formulation and along with its balanced pH technology, provides a moist, mildly acidic environment conducive to wound healing. The product contains ionic zinc, vitamin A, vitamin B6, calcium and magnesium. The DermaGran®-B Wound Dressing is used for skin ulcers (diabetic, venous stasis, stage II-IV pressure sores), surgical incisions, and superficial injuries such as partial thickness thermal burns, superficial lacerations, cuts, or abrasions.

#### **EXTENDED CARE**

31002100 DermaGran®-B 3 Oz. Each

PRIMARY CARE
61-SPD03 DermaGran®-B 3 Oz. 12/Case

#### **Ointment**

DermaGran® Ointments are pH balanced skin protectants.

These products are formulated with a patented zincnutrient formulation and along with their balanced pH technology provide a moist, mildly acidic environment conducive to wound healing. In addition, these products contain ionic zinc, vitamin A, calcium and magnesium. These ointments should be used for stage I and II pressure sores and preventative skin care.

#### EXTENDED CARE

 42001400
 DermaGran®
 4 Oz.
 Each

 42021400
 DermaGran®
 4 Oz. Jar
 Each

 PRIMARY CARE

 61-DT4
 DermaGran®
 4 Oz.
 12/Case

 61-DG4
 DermaGran®
 4 Oz. Jar
 12/Case

#### **Barrier Cream**

DermaGran® BC is a protective, long-lasting perineal barrier with a patented base formulation designed to be distinctly effective in shielding against urine, feces and other wetness due to incontinence. The non-greasy pH-balanced technology contains zinc and vitamin A for added protection against minor skin irritations.

#### EXTENDED CARE

35021400 DermaGran® BC 3.5 Fl. Oz. Each

PRIMARY CARE
61-ND35 DermaGran® BC 3.5 Fl. Oz. 12/Case

6 DermaGran is a registered trademark of Derma Sciences, Inc.



# **McKesson Medical-Surgical**

8741 Landmark Road Richmond, VA 23228

www.mckesson.com 800.446.3008