

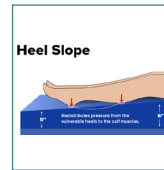


Gravity 8 Deluxe Long Term Care Pressure Redistribution Mattress

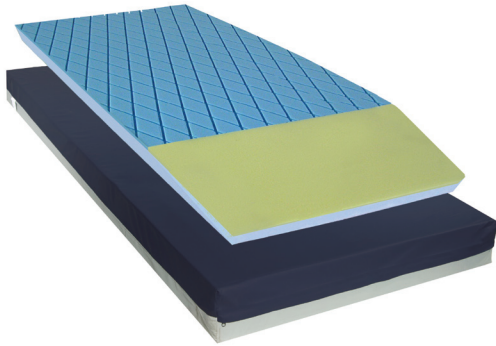

INSET B

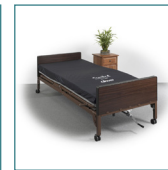

15886

- High quality therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- High end therapeutic foam mattress maximizes value and comfort with 3 full layers of latex free foam
- Full width top layer uses unique die-cutting to optimize zoned pressure redistribution
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High density foam perimeters provide solid edge support and facilitate safe patient transfers
- Available with or without raised side rails
- Heel section is sloped with channel cut foam center offloading pressure from vulnerable heels
- Removable multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633

PRODUCTS

Description

15886	Multi-Layered/Multi-zoned Foam Mattress, E0184, 1/ea
15870	Multi-Layered/Multi-zoned Foam Mattress, E0184, 1/ea
15884	Multi-Layered/Multi-zoned Foam Mattress, E0184, 1/ea
15888	Multi-Layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter and Cut-out, E0184, 1/ea
15885	Multi-Layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter and Cut-out, E0184, 1/ea

SPECIFICATIONS

Description

Cover Material	Stretch, Fluid-Resistant Nylon Top Waterproof Vinyl Bottom
Mattress Layers of Foam	3
Side Rails	Yes
Dedicated Heel Section	Yes
Fire Rating	CFR 16 Part 1633
Weight Capacity	400 lbs.
Warranty	10 Year Limited Foam 1 Year Cover

ACCESSORIES

Description

15870C	Replacement Cover for 15870, 1/ea
15884C	Replacement Cover for 15884, 1/ea
15885C	Replacement Cover for 15885, 1/ea
15886C	Replacement Cover for 15886, 1/ea
15888C	Replacement Cover for 15888, 1/ea
158709C	Replacement Cover