



## Multi-Ply 6500 Dynamic Elite Pressure Redistribution Foam Mattress

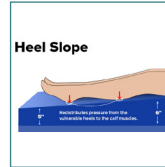

HEEL SLOPE

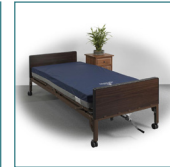

6500-DE-1-FB

- High quality therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- Drive's premium therapeutic foam mattress maximizes value and comfort with 4 full layers of latex free foam
- Full width gel infused, visco-elastic foam top layer dissipates heat providing a cool and dry environment to help preserve skin integrity
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High density foam perimeters provide solid edge support and facilitate safe patient transfers
- Available with or without raised side rails
- Heel section is sloped with patented LuraQuilt foam offloading pressure from vulnerable heels
- Removable premium, multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633

### PRODUCTS

| PRODUCTS        | Description   |
|-----------------|---|
| 6500-DE-2-FB    | Multi-layered/Multi-zoned Foam Mattress, 1/ea                                 |
| 6500-DE-2-RR-FB | Multi-Layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter Cut-out, 1/ea |

### SPECIFICATIONS

| SPECIFICATIONS          | Description   |
|-------------------------|---|
| Cover Material          | Fluid-Resistant Stretch Top   Waterproof Vinyl Bottom |
| Mattress Layers of Foam | 4   |
| Side Rails              | Yes   |
| Dedicated Heel Section  | Yes   |
| Fire Rating             | CFR 16 Part 1633 Federal Fire Code                    |
| Weight Capacity         | 400 lbs.  |
| Warranty                | 15 Year Limited Foam   1 Year Cover                   |