

ULTRA


Foam Dressings Featuring AFM® Technology



Every Milliken ULTRA and ULTRA Silver dressing utilizes Active Fluid Management® Technology created by Milliken Healthcare Products, LLC. This three-part design combines two layers of high-performance fabric using a micro-knit process, creating a moisture transfer mechanism that pulls away and traps excess exudate. The combination of these fabrics keeps the wound bed moist, the periwound skin healthy, and harmful exudate away from a healing wound.

The patented foam technology in ULTRA and ULTRA Silver increases the clinically relevant capacity of the dressing by pulling excess exudate through the dressing and across the entire foam substrate. This increased absorbing capacity allows for superior moisture management, longer wear times, and ultimately fewer dressing changes.

ULTRA Silver combines the moisture management benefits of AFM Technology's three-part design with a proprietary ceramic silver ion technology. This controlled release of ionic silver inhibits a broad spectrum of infection-causing microbes including antibiotic-resistant strains such as MRSA. The release of silver into the dressing is activated by the body's own sodium and begins inactivating pathogens within 24 hours of dressing application. Antimicrobial activity demonstrated up to three days.

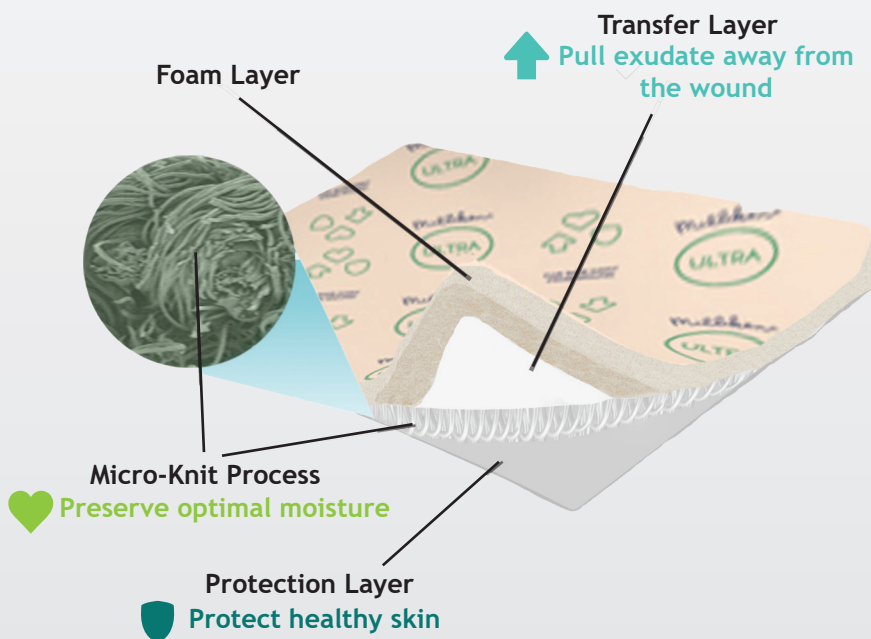

INDICATIONS FOR USE

- Partial thickness burns
- Incisions
- Skin grafts
- Donor sites
- Lacerations
- Abrasions
- Stage I-IV dermal ulcers (vascular, venous, pressure and diabetic)

BENEFITS

- Removes excess exudate through Active Fluid Management
- Minimizes risk of maceration
- Protects healthy, intact skin
- Provides significant cost savings
- Allows for fewer painful dressing changes





Milliken™



MILLIKEN HEALTHCARE PRODUCTS, LLC

ORDERING INFORMATION



ULTRA	Size	Code	HCPCS	Packaging
	4" x 5"	3000006840	A6210	10 ea/bx 10 bx/cs
	8" x 20"	3000112221	TBD	2 ea/bx 5 bx/cs



ULTRA Silver	Size	Code	HCPCS	Packaging
	4" x 5"	3000051018	A6210	10 ea/bx 10 bx/cs
	8" x 20"	3000112212	TBD	2 ea/bx 5 bx/cs

Three Easy Ways to Order:

1. Preferred Distributor
2. Preferred DME
3. Direct from Milliken



pull preserve protect

Pull exudate away from the wound.
Preserve optimal moisture.
Protect healthy skin.

medproducts@milliken.com | healthcare.milliken.com | 1-866-491-6556