

Live the Longer  
Life You Love®



## 5800 STABLE RAIL

*The Bed Rail With Legs*





EASILY SIT UP IN BED



EASILY STAND UP



ORGANIZER KEEPS ITEMS CLOSE BY

## Features & Benefits

- **Extendable Legs** – Provides added stability when standing up.
- **Ergonomic Cushion Handle** – Allows for easy transfer in and out of bed.
- **4-Pocket Organizer** – Provides storage space for handy items.
- **Anti-slip grips** – Secures rail in-between mattress and bed frame.
- **Height Adjustment** – Accomodates any home or hospital bed.
- **Easy Installation** – Installs in seconds with no tools required.

## Specifications

Height adjustment from floor to top of handle: 21" - 44"

Height adjustment from mattress base to top of handle: 17.5" - 22.5"

Rail width: 19"

Height range from floor to horizontal support structure under mattress: 13.5" - 21.5"

Depth of horizontal support structure under mattress: 25"

Weight of Product: 8 lbs

Weight capacity: 300 lbs

Made of: Steel frame with durable powder coat finish

Case of 4 Dimensions: 28" L x 20.25 W x 5.75" H; 32 lbs

## How Do We Compare

Tool-Free Assembly



Extendable Legs for Added Stability



4-Pocket Organizer Pouch Included



Provides Assistance In and Out of Bed



Durable Powder Coat Finish



## Easy Assembly (No Additional Tools Required)




Adjust base tubes to accommodate left or right side of bed



Attach base tubes to main handle



Insert base tubes between mattress and box spring

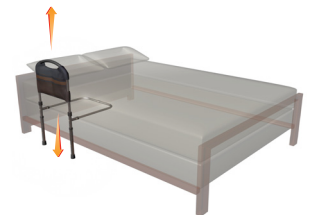

Adjust height of handle and legs as desired